


Cadenza CD-C18

250 x 2 mm
250 x 3 mm

Application


Column size comparison for elution of investigational medication

医薬品候補化合物の溶出に関するカラムサイズ比較

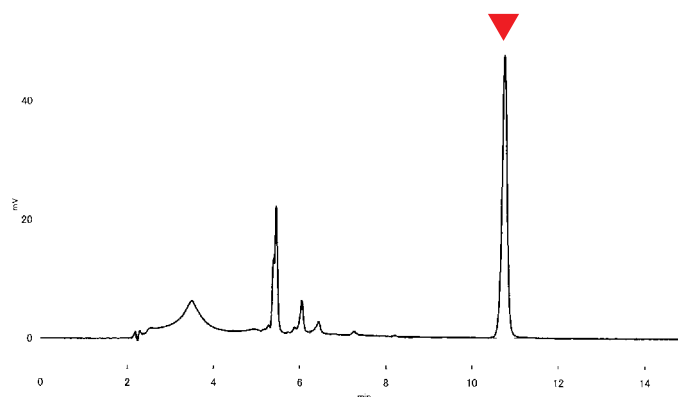


Cadenza CD-C18
250 x 2 mm

0.2 mL/min
500 µL inj.

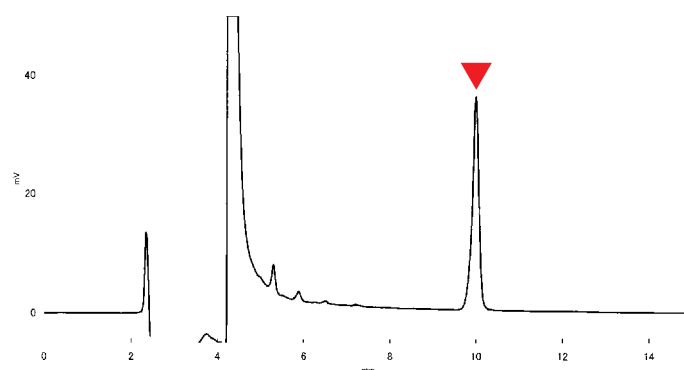

Cadenza CD-C18
250 x 3 mm

0.4 mL/min
500 µL inj.

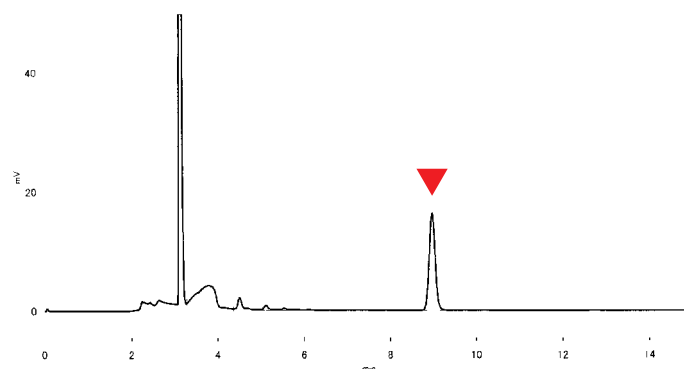

Brand-A
150 x 6 mm

1 mL/min
500 µL inj.

acetonitrile / water / phosphoric acid = 600 / 400 / 1, 40 °C, 280 nm

Courtesy of Hiroshi Okumura