

Cadenza CD-C18


100 x 2 mm

Application

LC-MS Application for Endocrine Disruptor

LC-MSアプリケーション (環境ホルモン, オクチルフェノール)

API4000 (TAKARA BIO INC.)

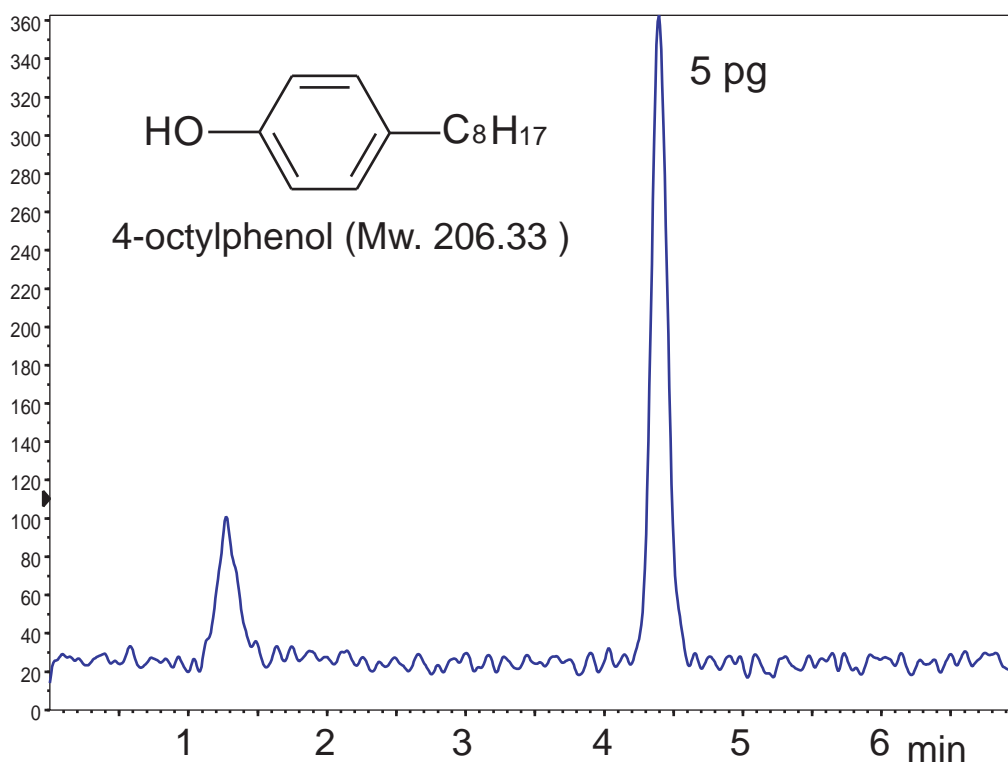

Cadenza CD-C18, 100 x 2 mm
 80% methanol, 0.2 mL/min, 40 °C, 5µL(1pg/µL)
 API4000: ESI, MRM Negative, Q1/Q3 : 205 / 106

Courtesy of J.Watanabe, TAKARA BIO INC.