
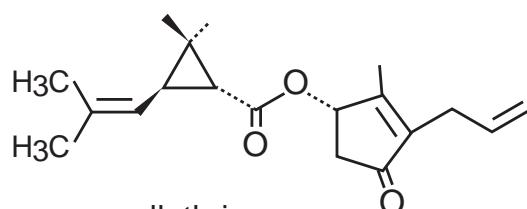
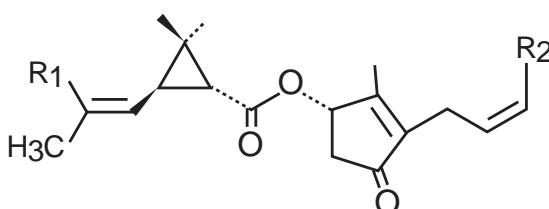


Unison UK-C18

100 x 4.6 mm

Application

Pyrethrins
天然ピレトリン類

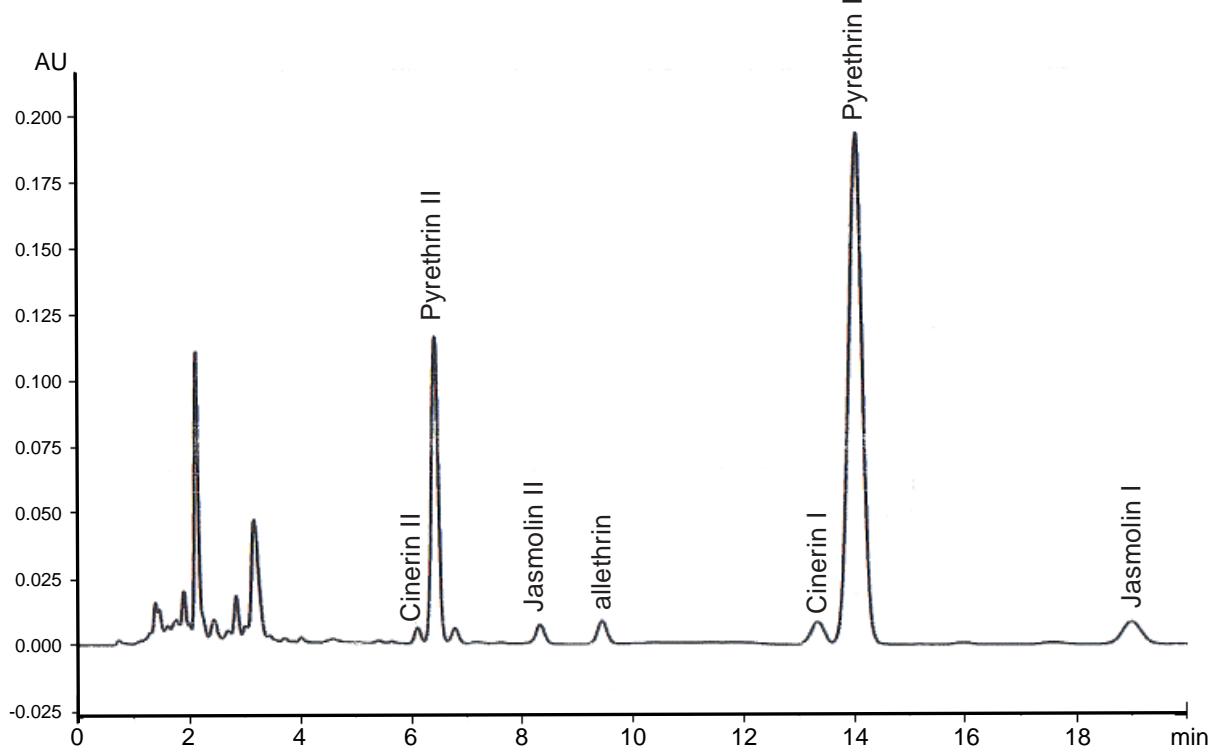

	R ₁	R ₂
Pyrethrin I	CH ₃	CH=CH ₂
Pyrethrin II	COOCH ₃	CH=CH ₂
Cinerin I	CH ₃	CH ₃
Cinerin II	COOCH ₃	CH ₃
Jasmolin I	CH ₃	CH ₂ CH ₃
Jasmolin II	COOCH ₃	CH ₂ CH ₃

Unison UK-C18, 100 x 4.6 mm
water / acetonitrile = 35 / 65
1.0 mL/min, 40 °C
UV at 220 nm, 1µL

Courtesy of Makoto IHARA and Prof. Koichiro KOMAI
Graduate School of Agriculture, Kinki University