








Unison UK-C18

75 x 3 mm

Application

Simultaneous Analysis Of Water-soluble And Fat-soluble Vitamins

水溶性ビタミンと脂溶性ビタミンの同時分析

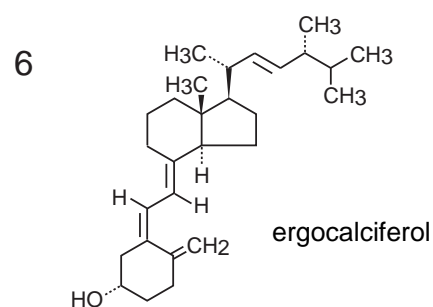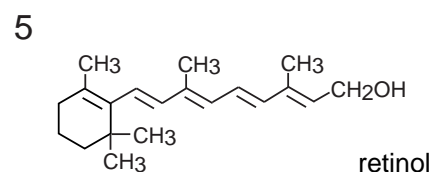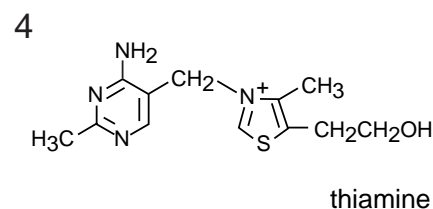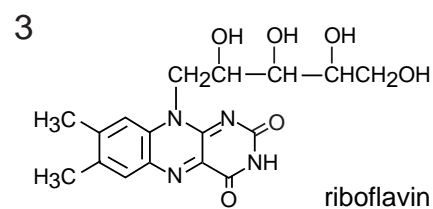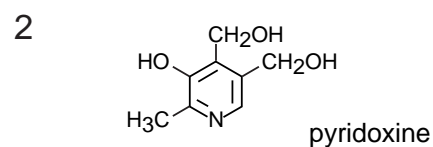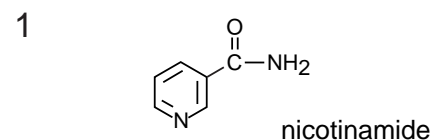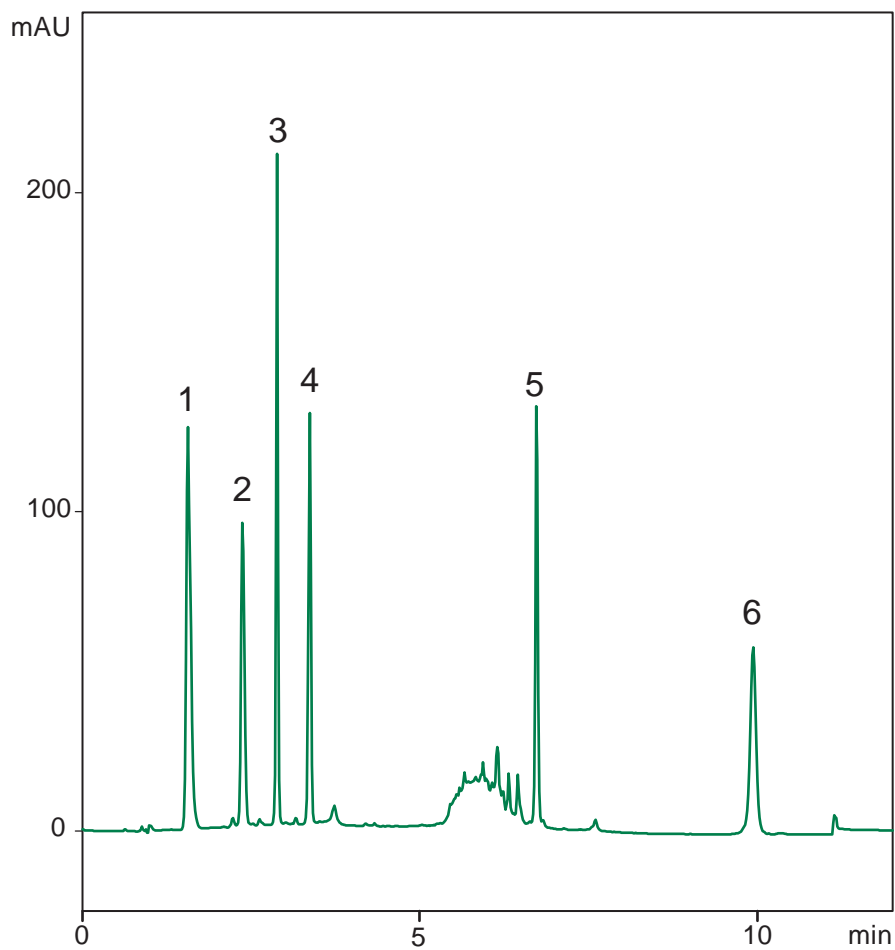

Unison UK-C18, 75 x 3 mm

A: water / heptafluorobutyric acid = 100 / 0.1

B: acetonitrile

10 - 40%B (0 - 4min)

40 - 100%B (4 - 4.5min)

100%B (4.5 - 10min)

0.5 mL/min, 37 °C, 260 nm, 0.24 - 2.4µg(2µL)