
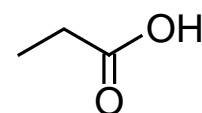


Unison UK-C18
 Unison UK-C8
 Unison UK-Phenyl


50 x 4.6 mm

Application


Short-Chain Fatty Acids
 短鎖脂肪酸

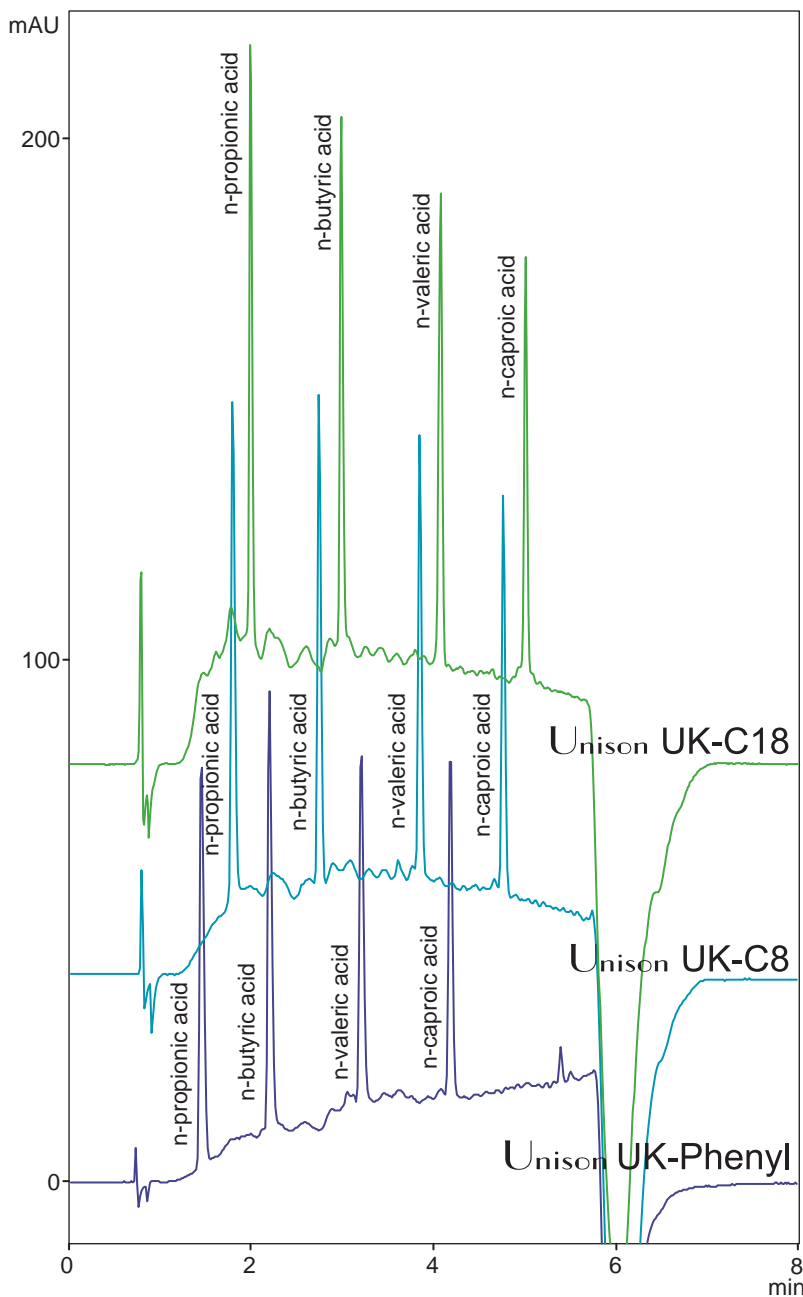


n-propionic acid

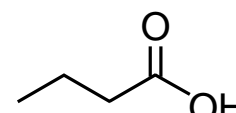

n-butyric acid

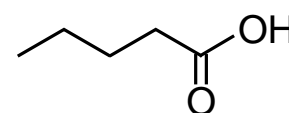

n-valeric acid

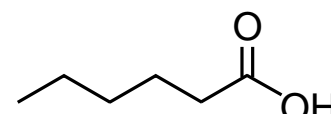

n-caproic acid

Unison, 50 x 4.6 mm
 A: water / TFA = 100 / 0.1
 B: acetonitrile / TFA = 100 / 0.1
 0.5 - 50%B (0-5min)
 1 mL/min, 37 deg.C, 210 nm