


Cadenza CL-C18  
 Unison UK-C18  
 Cadenza CD-C18  
 Unison UK-Silica

75 x 4.6 mm  
 250 x 4.6 mm


Application

### Antioxidants



酸化防止剤

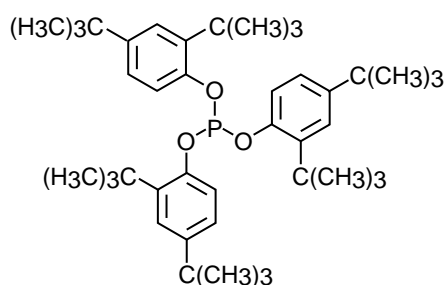

irgafos 168

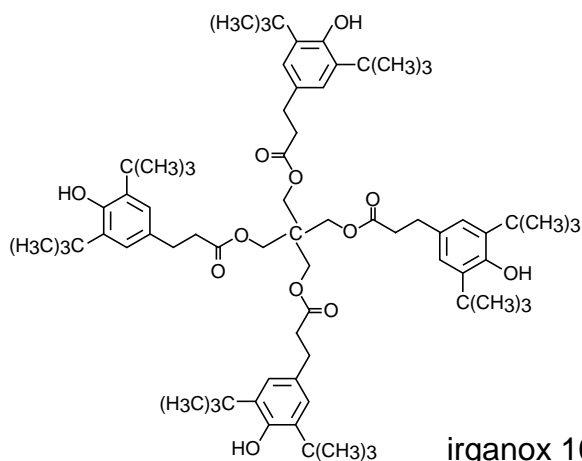

irganox 1010

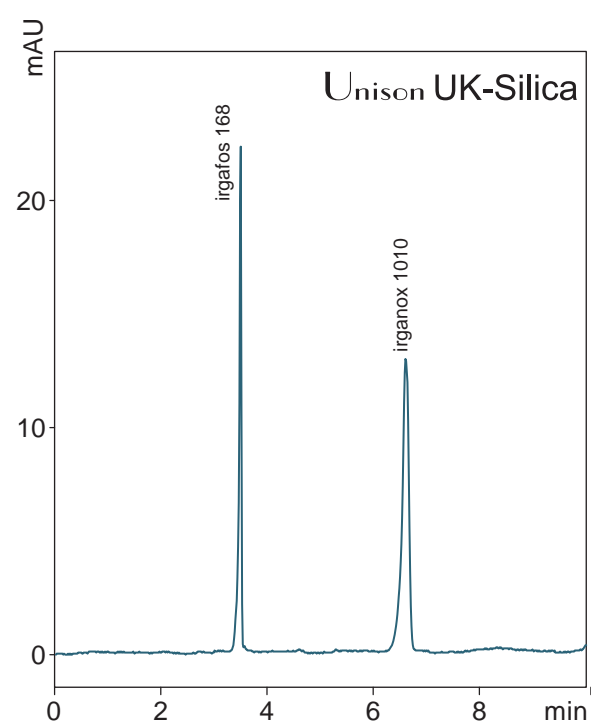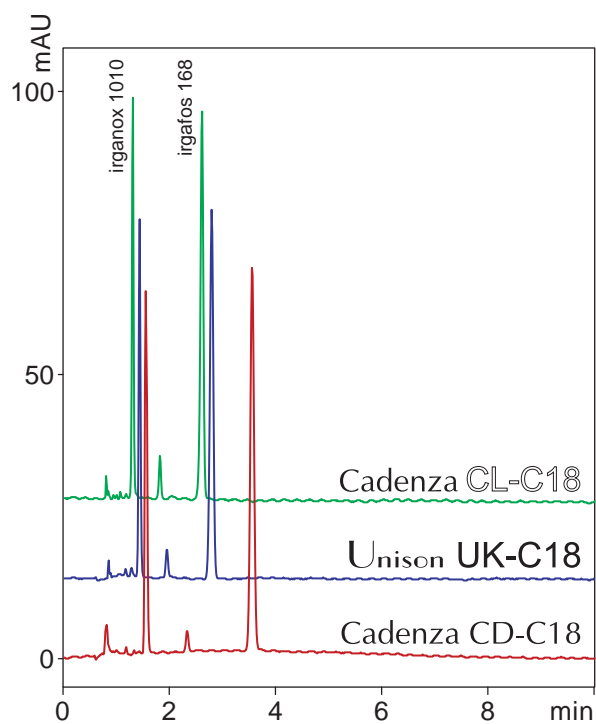

75 x 4.6 mm  
 acetonitrile / tetrahydrofuran  
 = 80 / 20  
 1 mL/min, 37 deg.C, 230nm

Unison UK-Silica  
 250 x 4.6 mm  
 hexane / ethyl acetate = 90 / 10  
 1 mL/min, 37 deg.C, 270nm