

Unison UK-C8




75 x 3 mm

Application

LC-MS Application for Coenzyme Q10

LC-MSアプリケーション (コエンザイム Q10)

API4000 (TAKARA BIO INC.)

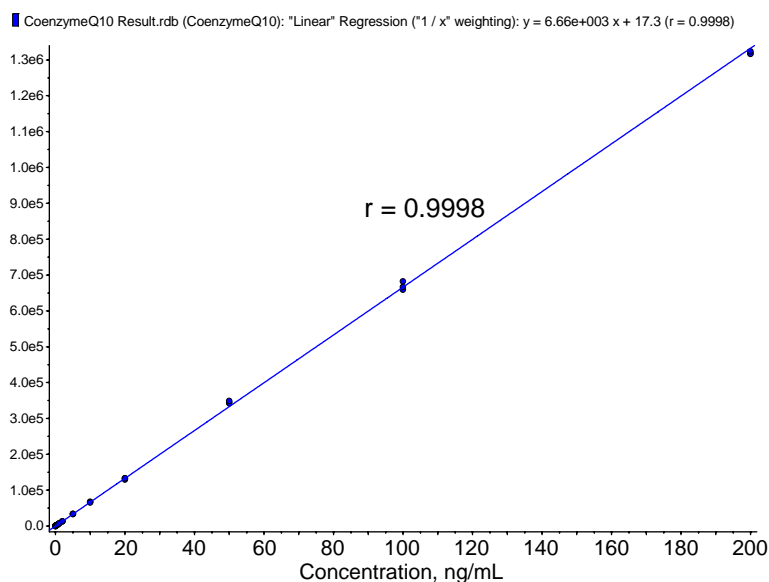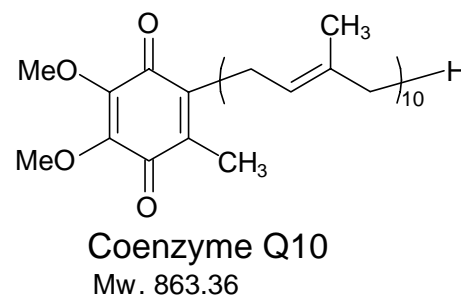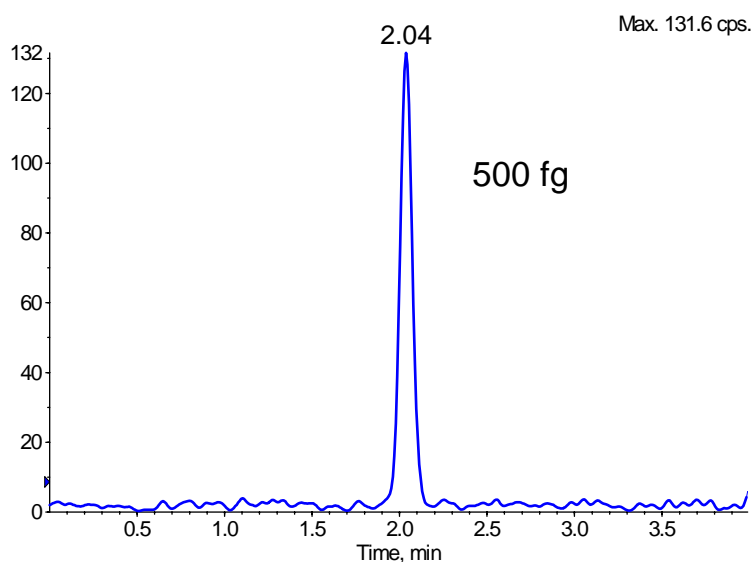

Unison UK-C8, 75 x 3 mm

0.1% formic acid / acetonitrile / isopropanol = 1 / 5 / 5

0.8 mL/min, 10 uL(50 fg/uL)

API4000: APCI, MRM Positive, Q1/Q3 : 863.5 / 197

Courtesy of J.Watanabe, TAKARA BIO INC.