



Unison UK-Phenyl
 Unison UK-C8
 Cadenza CD-C18
 Unison UK-C18
 Cadenza CL-C18

75 x 3 mm

Application

Veterinary Medicine (Antibiotic)

動物用医薬品(抗生物質)

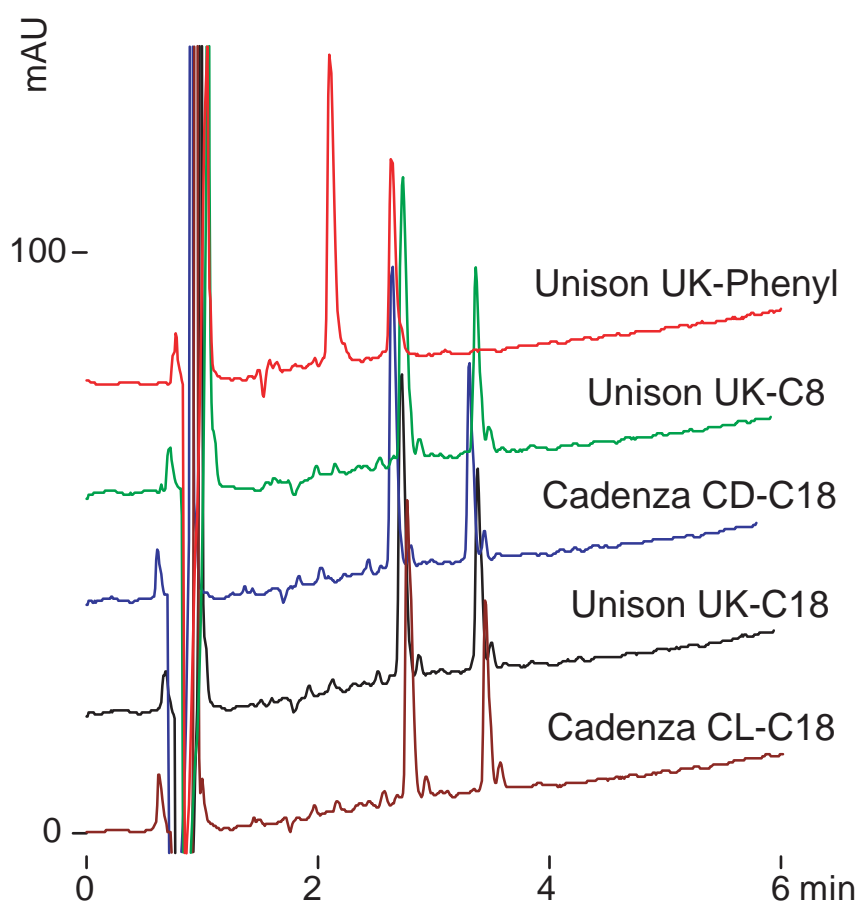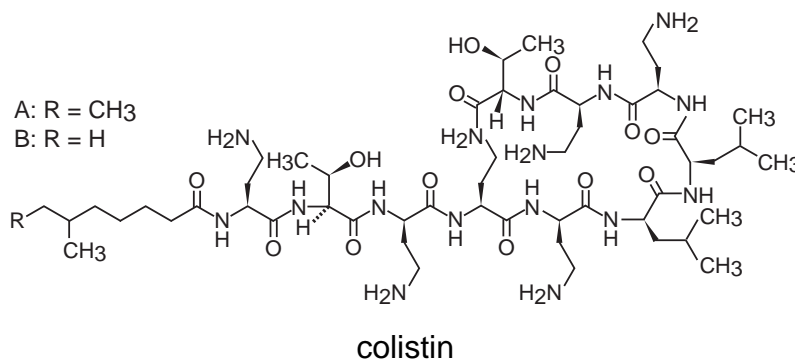

75 x 3 mm,
 A: water / trifluoroacetic acid = 100 / 0.1
 B: acetonitrile / trifluoroacetic acid = 100 / 0.1
 25-40 %B (0-5 min), 0.5 mL/min,
 37 deg.C, 220 nm, 1 ug/uL (20uL)