




Unison UK-Silica

250 x 4.6 mm

Application

Pyromellitic acids

ピロメリト酸

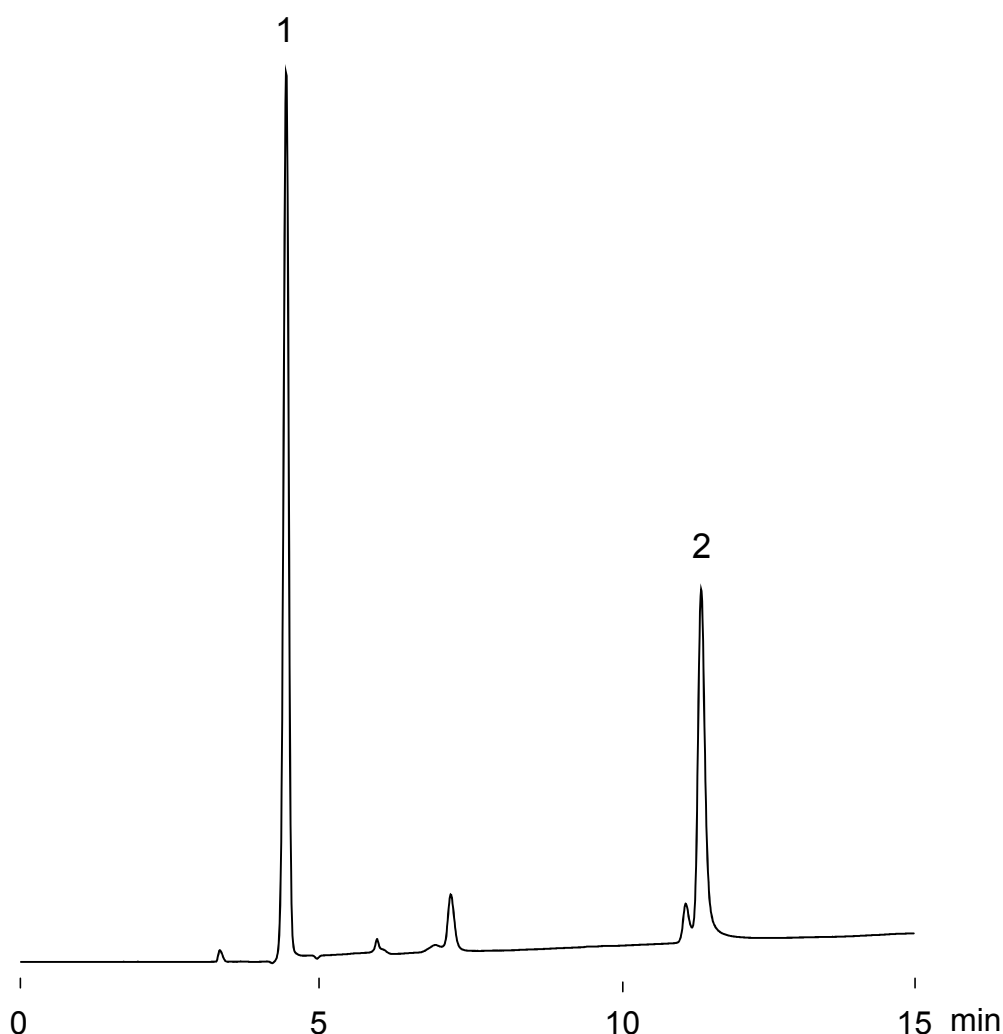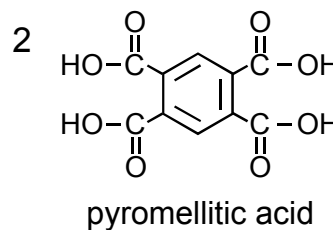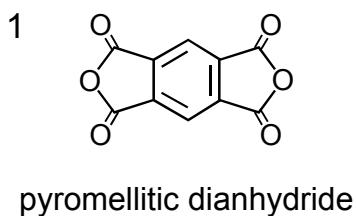

Unison UK-Silica, 250 x 4.6 mm

A: hexane / tetrahydrofuran / trifluoroacetic acid = 90 / 10 / 0.1

B: hexane / tetrahydrofuran / trifluoroacetic acid = 50 / 50 / 0.1

50-100%B (0-15min)

1mL/min (5MPa), 37°C, 260nm

5 µL (1mg/mL, tetrahydrofuran / trifluoroacetic acid = 100 / 0.1)