





Unison UK-C18
Cadenza CD-C18

75 x 4.6 mm


Application

Parasympatholytics

副交感神経遮断薬

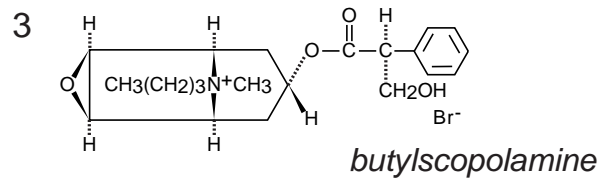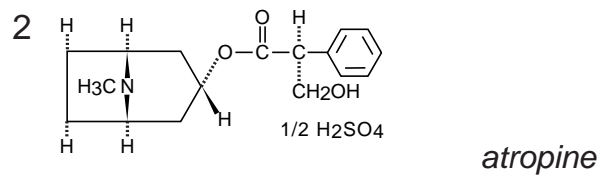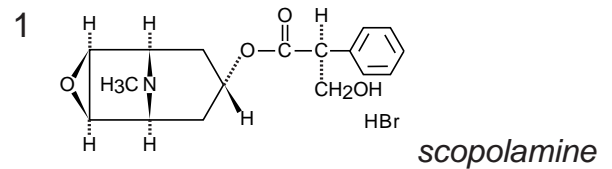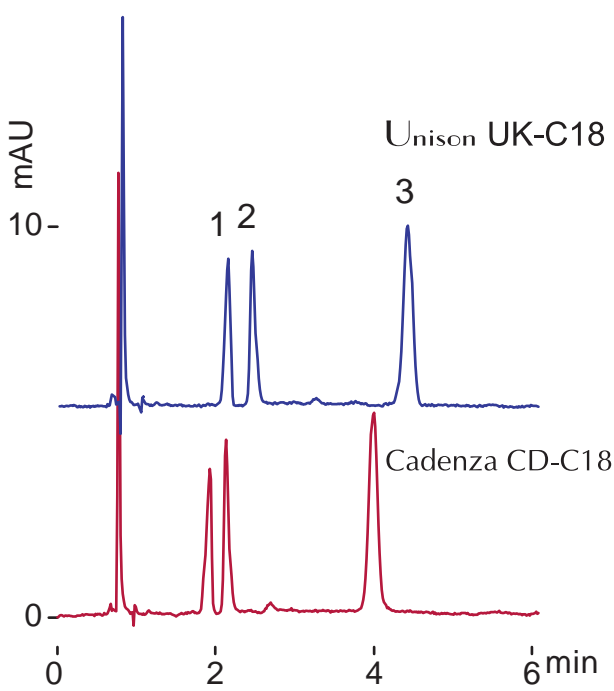

75 x 4.6 mm
100mM ammonium formate / acetonitrile
= 80 / 20
1 mL/min (6MPa), 37 deg.C, 260 nm
4uL (1-2ug)

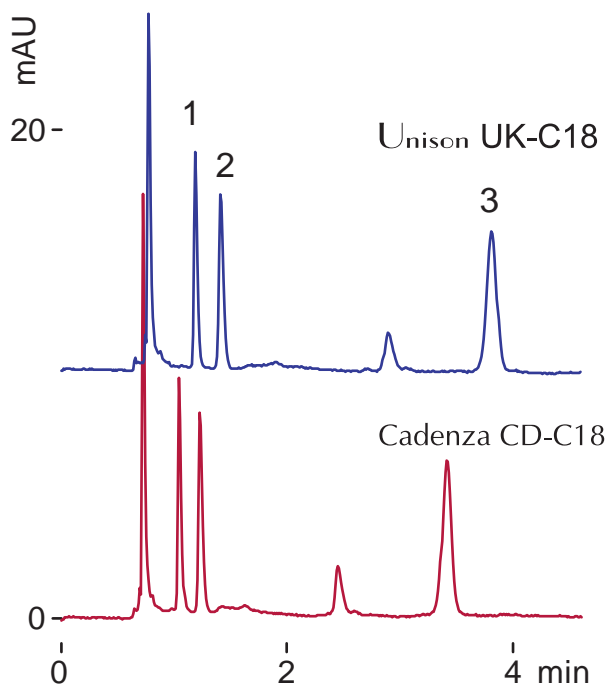

75 x 4.6 mm
water / acetonitrile / trifluoroacetic acid
= 70 / 30 / 0.1
1 mL/min (5MPa), 37 deg.C, 260 nm
4uL (1-2ug)