



Cadenza CD-C18

150 x 3 mm


Application

Aromatic Boronic Acids


芳香族ボロン酸類

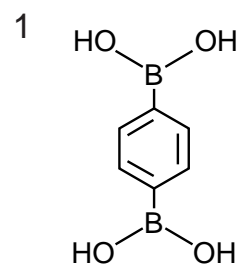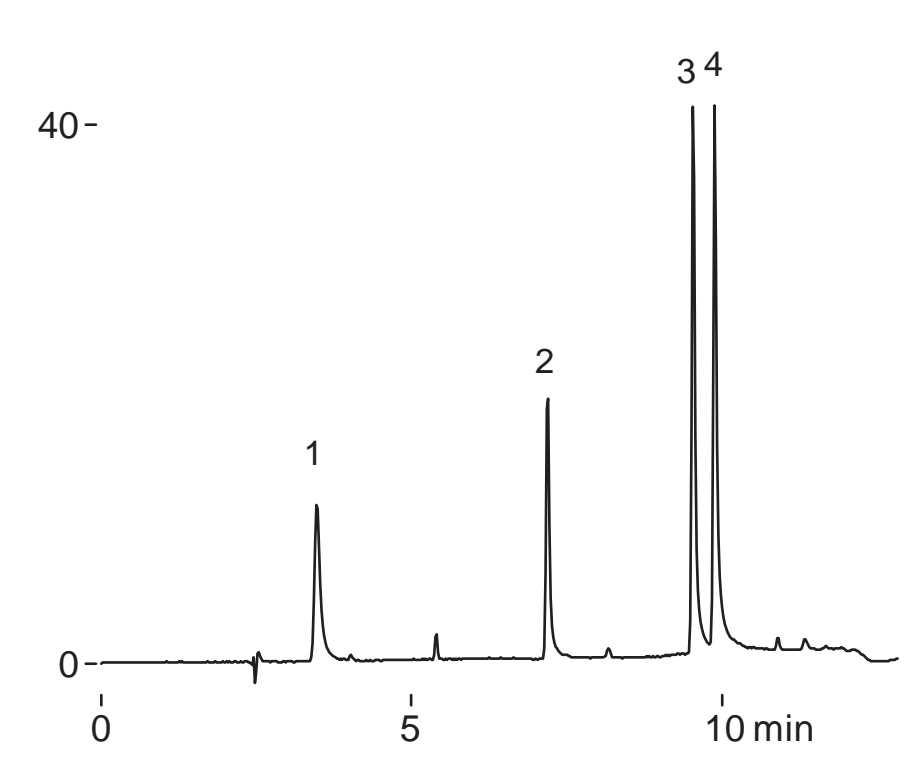


1,4-phenylenediboronic acid

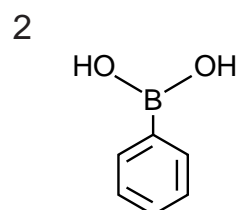

benzeneboronic acid

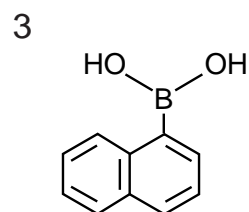

1-naphthaleneboronic acid

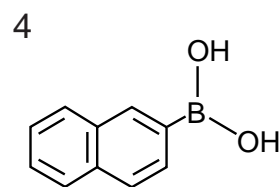

2-naphthaleneboronic acid

Cadenza CD-C18, 150 x 3 mm
 A: water /TFA = 100 /0.1
 B: acetonitrile /TFA = 100 /0.1
 35-75%B (0-10 min)
 1 mL/min, 37 deg.C, 260 nm
 1 uL (0.2-0.6ug)