


Unison UK-C18

250 x 2 mm

Application

## LC-MS/MS Application of Water Soluble Vitamins

LC-MS/MSアプリケーション(水溶性ビタミン)

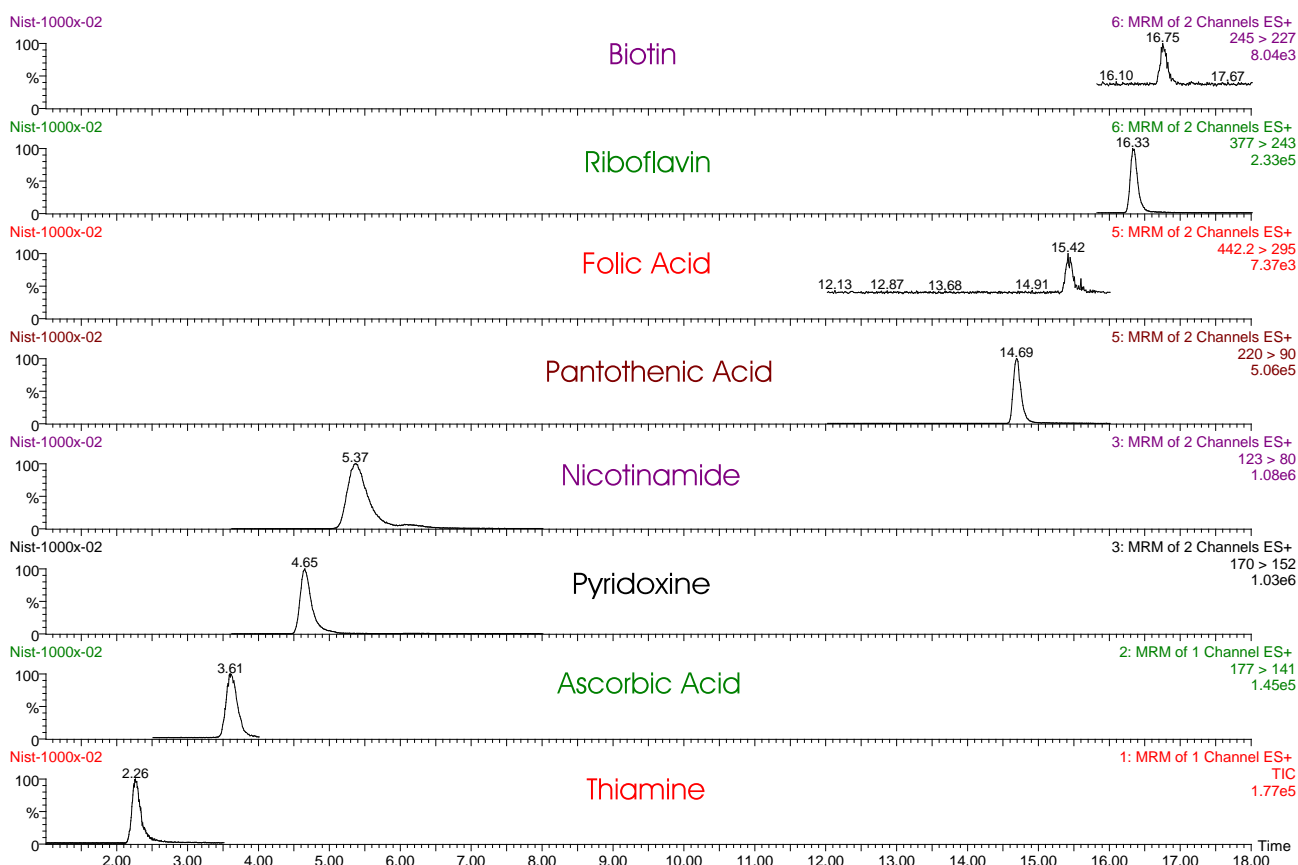

Unison UK-C18, 250 x 2 mm

A: 0.1% HCOOH in water

B: 0.1% HCOOH in acetonitrile

0%B (0-5min), 0-50%B (5-15min), 50-95%B (15-17min), 0%B (17-25min)

0.25mL/min, ambient, 5uL

ESI-MS, MRM Positive

Courtesy of Dr. Pei Chen, USDA, ARS, BHNRC, Food Composition Lab - Beltsville, MD, USA