


Unison UK-Amino  
Unison UK-C18

250 x 4.6 mm



Application

**Steroids**

ステロイド

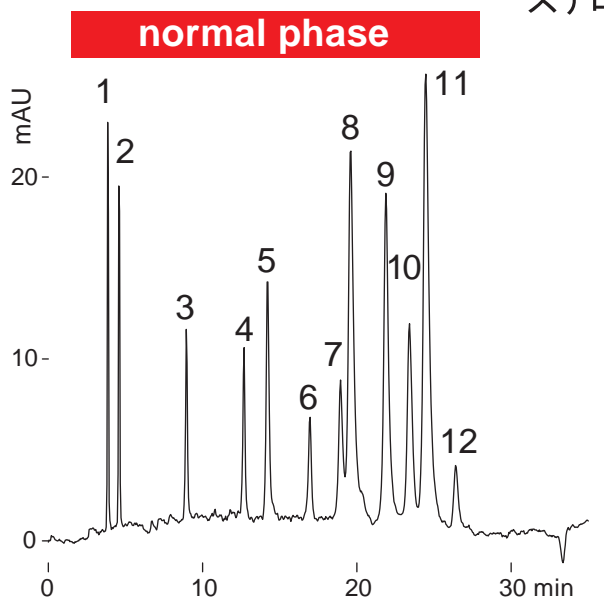

**Unison UK-Amino**, 250 x 4.6 mm  
A: hexane / acetic acid = 100 / 0.1  
B: ethanol / acetic acid = 100 / 1  
10-30%B (0-30 min), 1 mL/min  
37 deg.C, 260 nm, 2 uL (1-2ug)

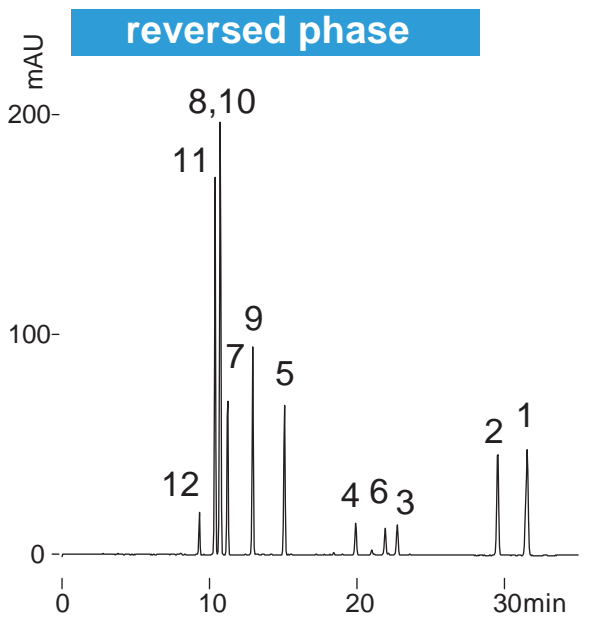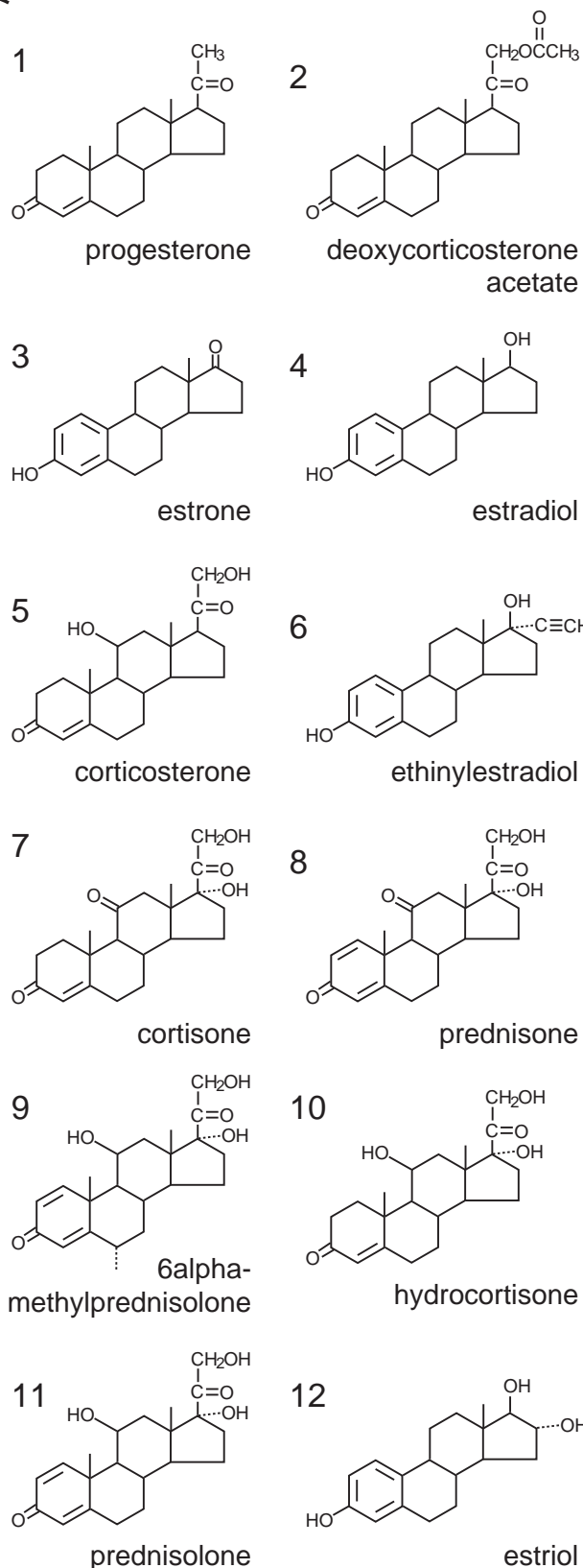

**Unison UK-C18**, 250 x 4.6 mm  
A: water, B: acetonitrile  
30-70%B (0-30 min), 0.8 mL/min  
37 deg.C, 260 nm, 2 uL (1-2ug)