


Cadenza CL-C18

250 x 4.6 mm


Application

Steroid Hormones


ステロイドホルモン

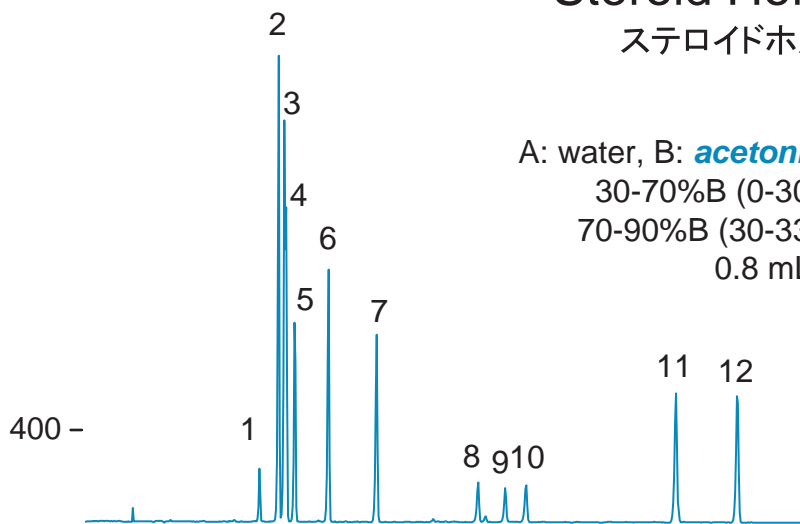


A: water, B: **acetonitrile**
30-70%B (0-30min)
70-90%B (30-33min)
0.8 mL/min

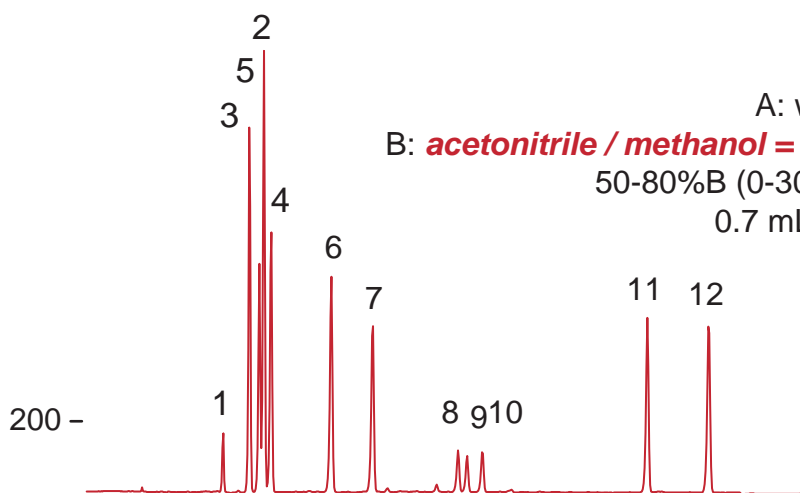

A: water
B: **acetonitrile / methanol = 1 / 1**
50-80%B (0-30min)
0.7 mL/min

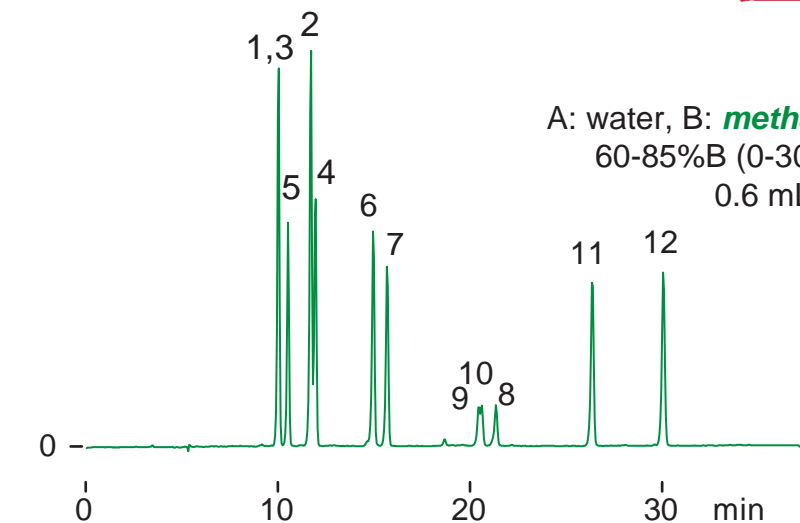

A: water, B: **methanol**
60-85%B (0-30min)
0.6 mL/min

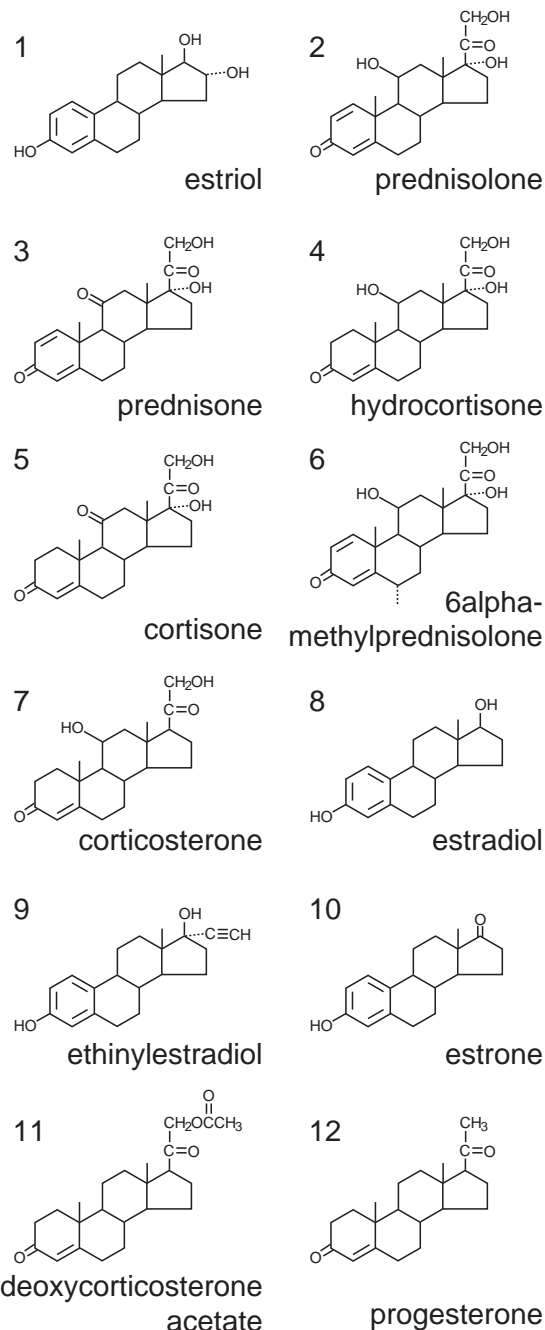

Cadenza CL-C18, 250 x 4.6 mm
37 deg.C, 260 nm