


Unison UK-Amino
Cadenza CL-C18

250 x 4.6 mm


Application

Surfactant (triton)

界面活性剤(トリトン)

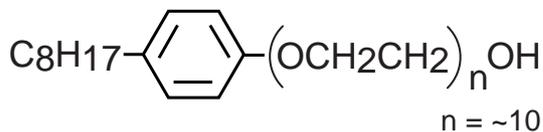


Triton X-100

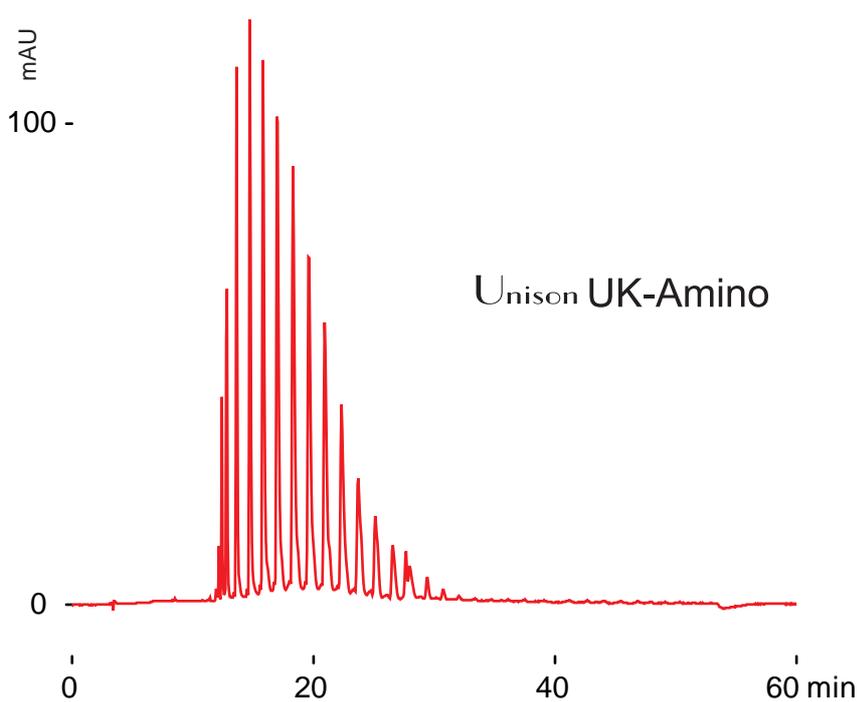

normal phase

Unison UK-Amino
250 x 4.6 mm,
A: hexane, B: ethanol
0.2-10%B(0-50min)
1 mL/min 37 deg.C
260 nm, 5 uL (10% solution)

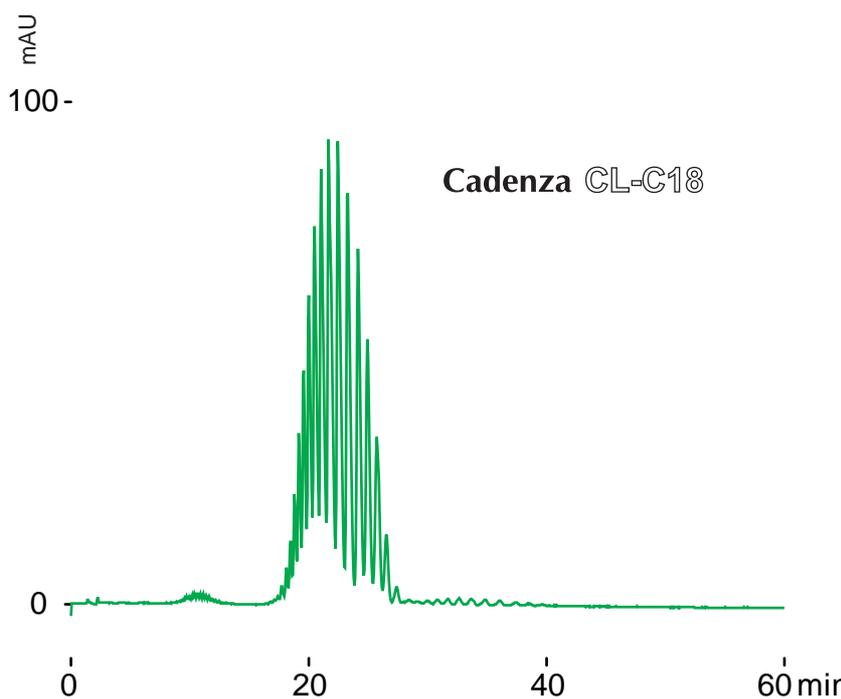

reversed phase

Cadenza CL-C18
250 x 4.6 mm
water / acetonitrile = 40 / 60
0.8 mL/min, 37 deg.C
260 nm, 5 uL (10% solution)