


Unison UK-Amino

250 x 4.6 mm


Application

Quercetin glycoside separation with an aminopropyl column


アミノプロピルカラムによるケルセチン配糖体の分離

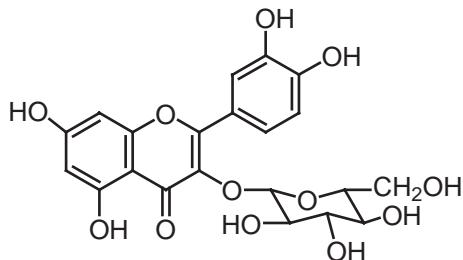

Quercetin-3-glucoside (Q3G)

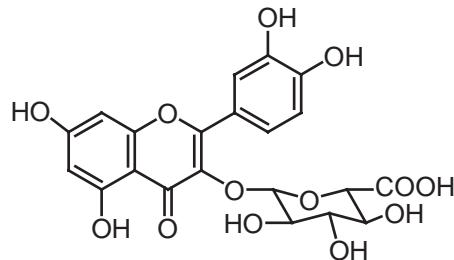

Quercetin-3-glucuronide (Q3GA)

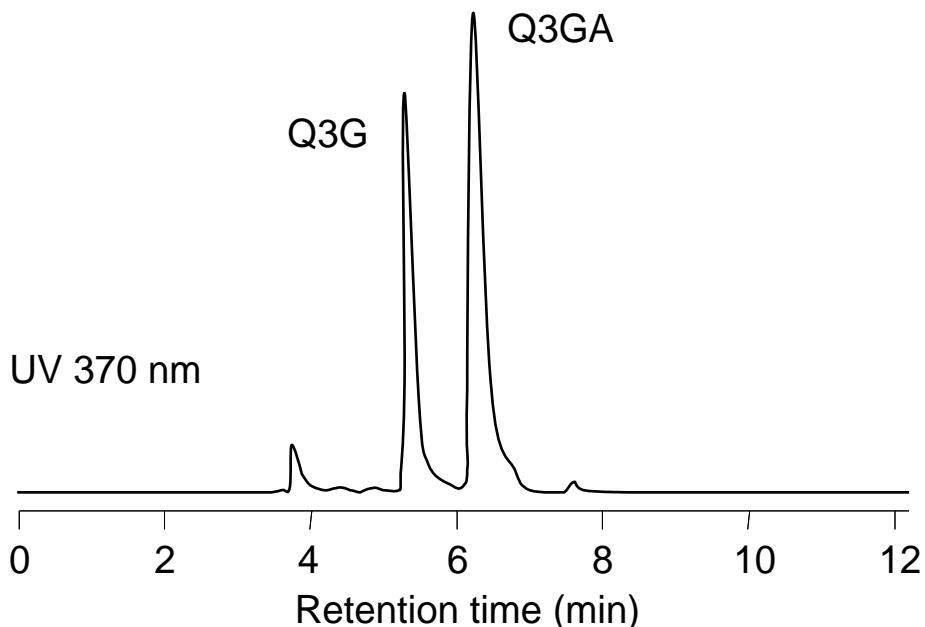

Column: UK-Amino, 250 x 4.6 mm
Flow: 0.8 ml/min
Solvent: 0.5% phosphoric acid/CH₃CN (15/85)

Courtesy of Dr. Y.KAWAI, The University of Tokushima, JAPAN