
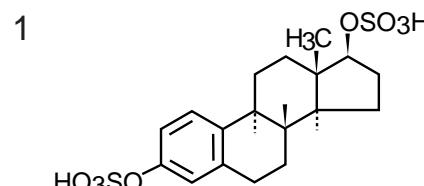


Cadenza CW-C18
 Unison UK-Phenyl
 Unison UK-C8
 Unison UK-C18
 Cadenza CL-C18
 Cadenza CD-C18


30 x 2 mm

Application


Estradiol and Metabolites エストラジオールと代謝物

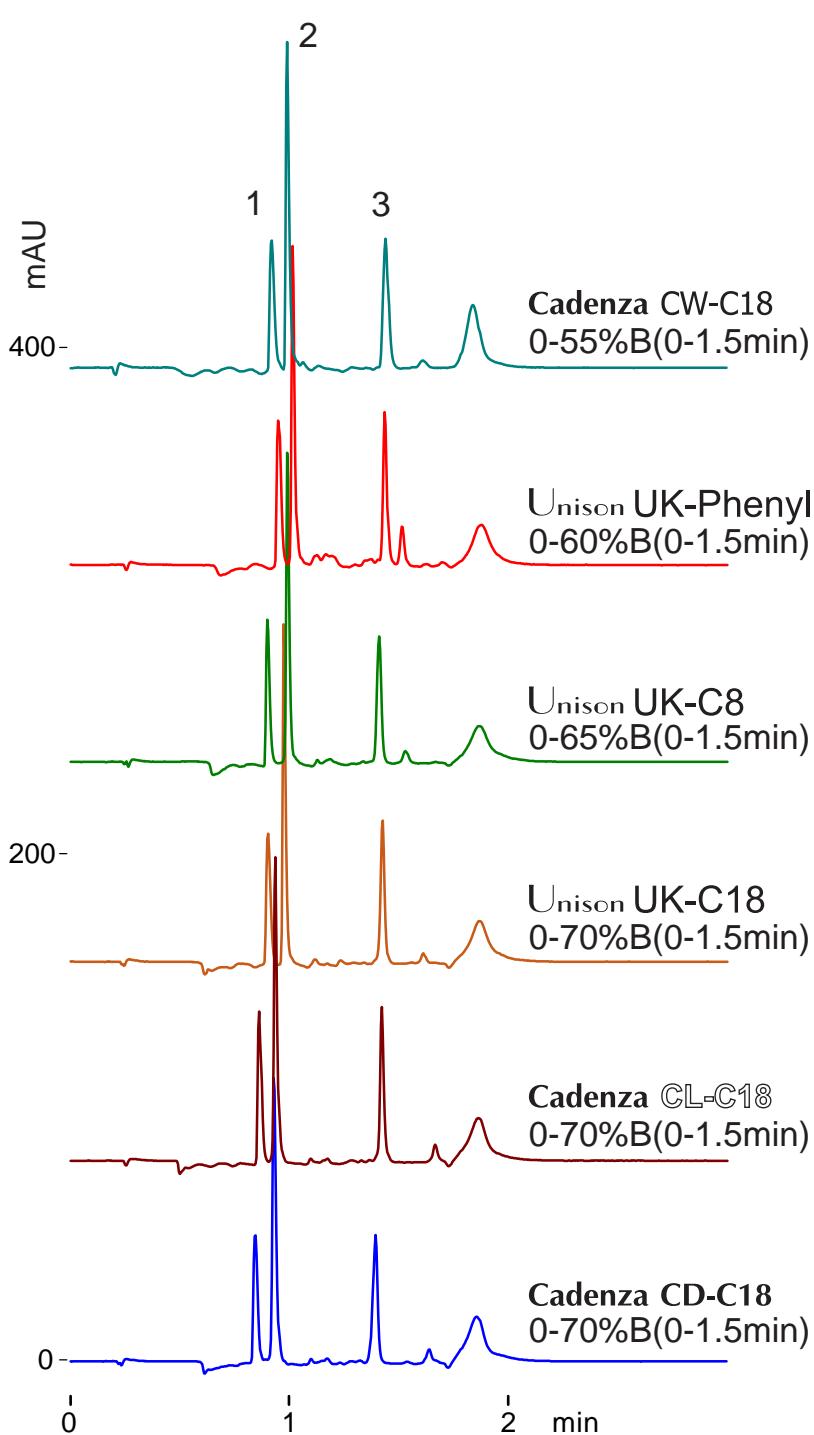

beta-Estradiol 3,17-disulfate

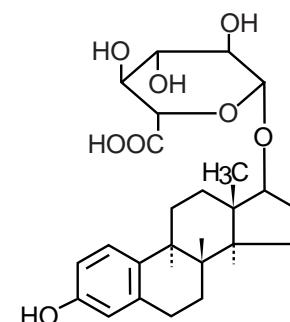

beta-Estradiol
17-(beta-D-glucuronide)

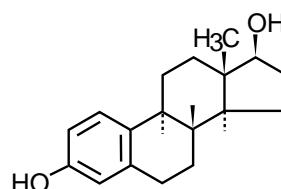

17beta-estradiol

30 x 2 mm,
 A: 50mM ammonium acetate
 B: acetonitrile
 0.5 mL/min, 37 deg.C
 280 nm, 1 uL (50-250ng)