


Cadenza CW-C18

150 x 4.6 mm


Application

Tetracycline antibiotics


テトラサイクリン系抗生物質

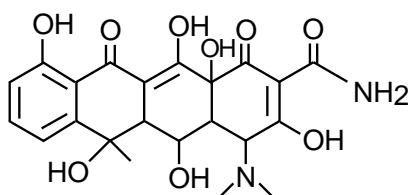


oxytetracycline (OTC)

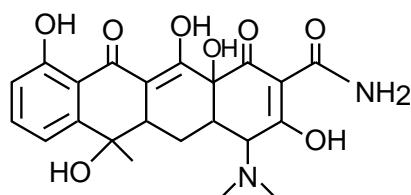

tetracycline (TC)

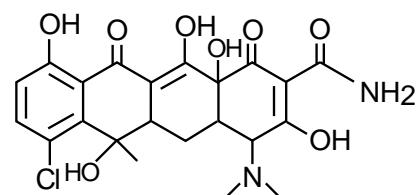

chlortetracycline (CTC)

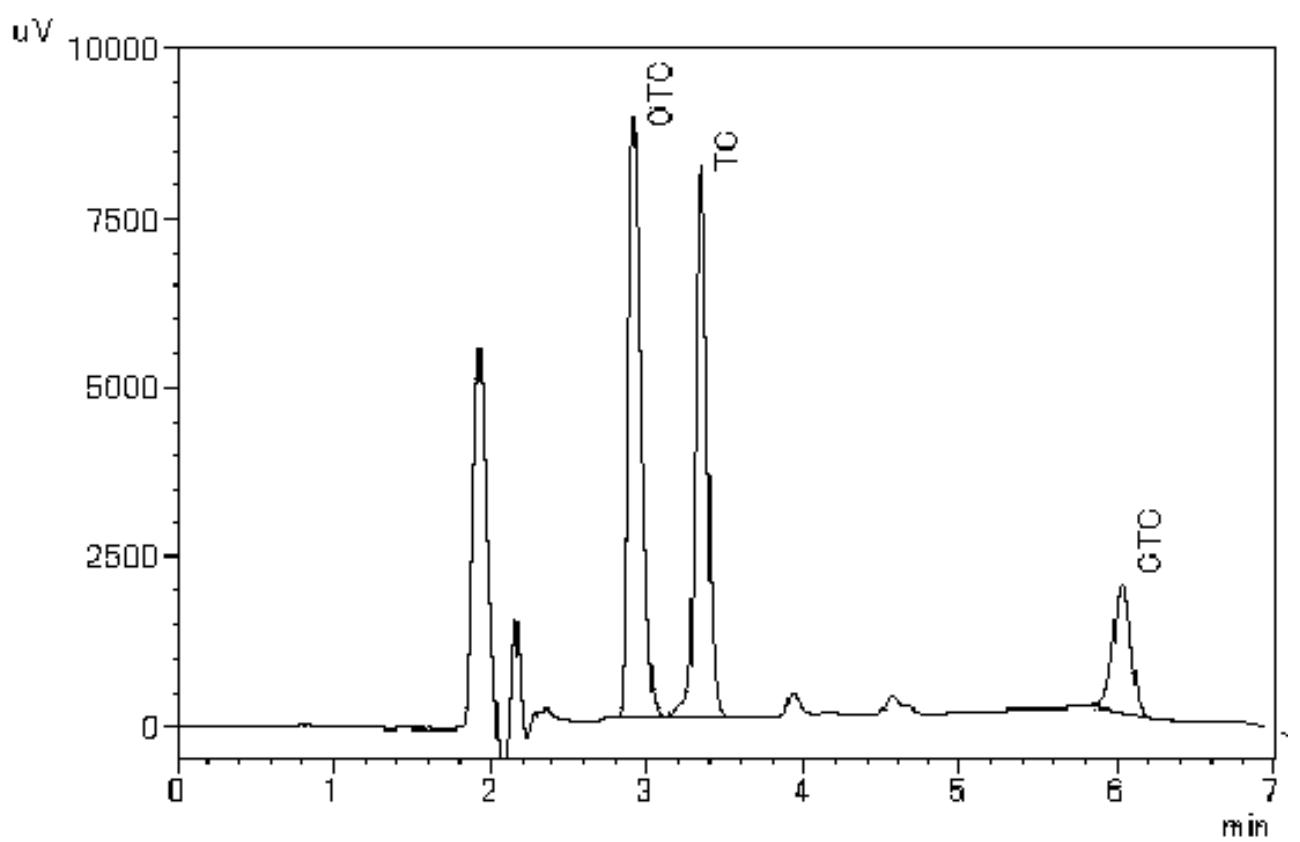

Cadenza CW-C18, 150 x 4.6 mm
water /acetonitrile /HCOOH = 80 /20 /0.1
1 mL/min, 40 deg.C, 260 nm, 20 uL (1ppm)

Courtesy of Mr. N. MAEDA, ASSOCIATION OF MEAT SCIENCE & TECHNOLOGY INSTITUTE

USA, Canada, and Mexico : Silvertone Sciences (USA) - info@silvertonesciences.com
Other Countries : Imtakt Corp.(JAPAN) - info@imtakt.com

TI410E-HE08