


Unison UK-Amino

250 x 3 mm


Application

Deoxynojirimycin, Glucose


デオキシノジリマイシン, グルコース

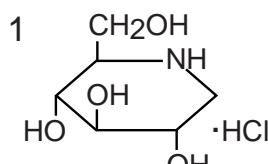

1-deoxynojirimycin
hydrochloride

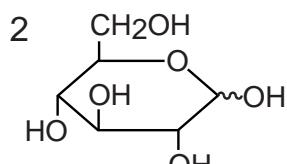

D-glucose

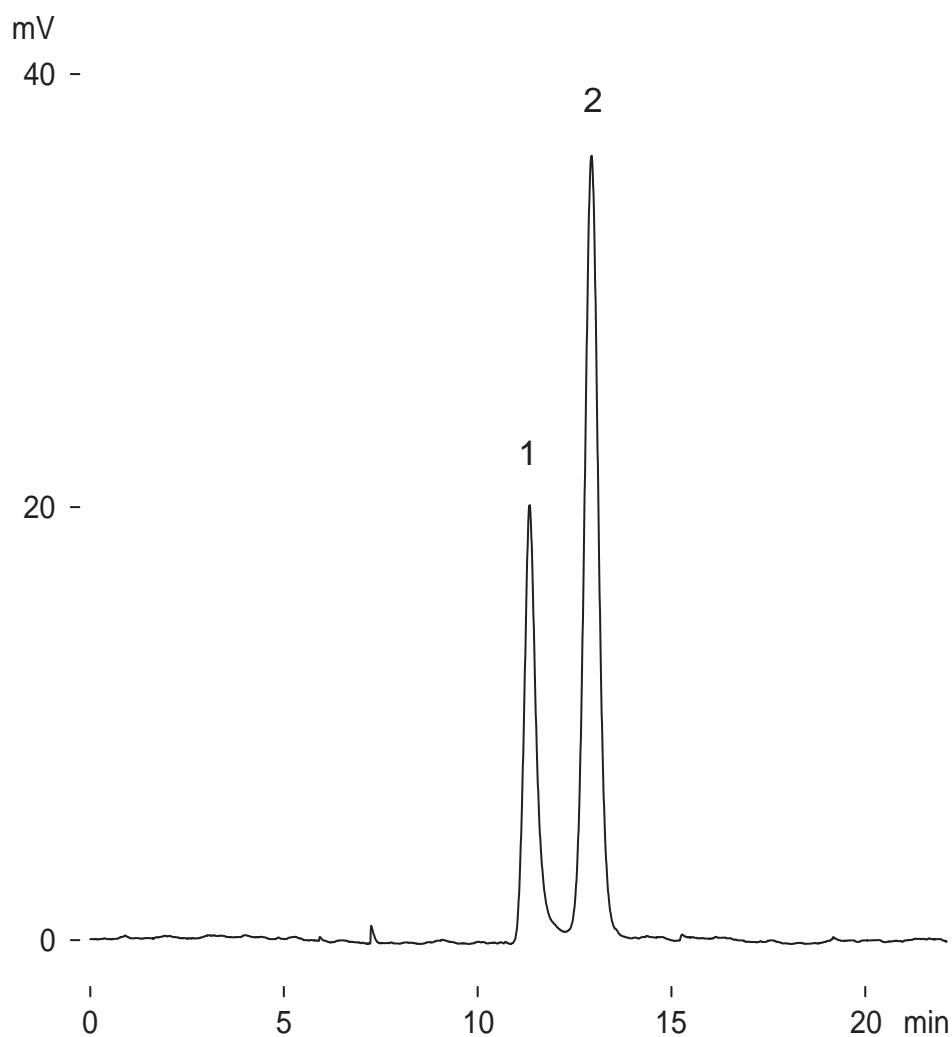

Unison UK-Amino, 250 x 3 mm
acetonitrile /10mM ammonium hydrogen carbonate = 90 /10
0.4 mL/min (6 MPa), 60 deg.C, ELSD, 3 uL (1.5 ug)

Courtesy of Prof. Takayuki Nakano, Kagoshima Immaculate Heart University, Japan