



Unison UK-Amino

100 x 3 mm


Application

### Sugar phosphates

糖リン酸

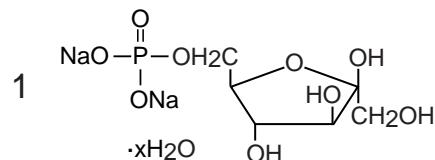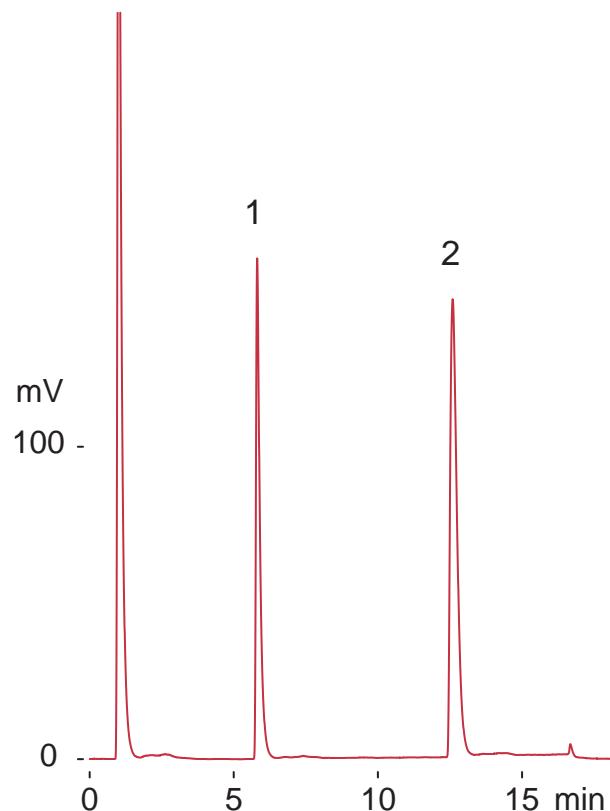


D-fructose 6-phosphate disodium hydrate  
(F6P)

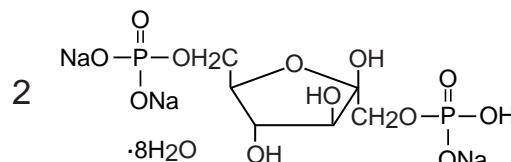

D-fructose 1,6-disphosphate trisodium octahydrate  
(F1,6P)

Unison UK-Amino, 100 x 3 mm  
A: water /acetic acid = 95 /5  
B: water /triethylamine = 94 /6  
0-40% B(0-15min), 0.4 mL/min (7MPa)  
37 deg.C, ELSD, 1.2 uL (4-8ug)

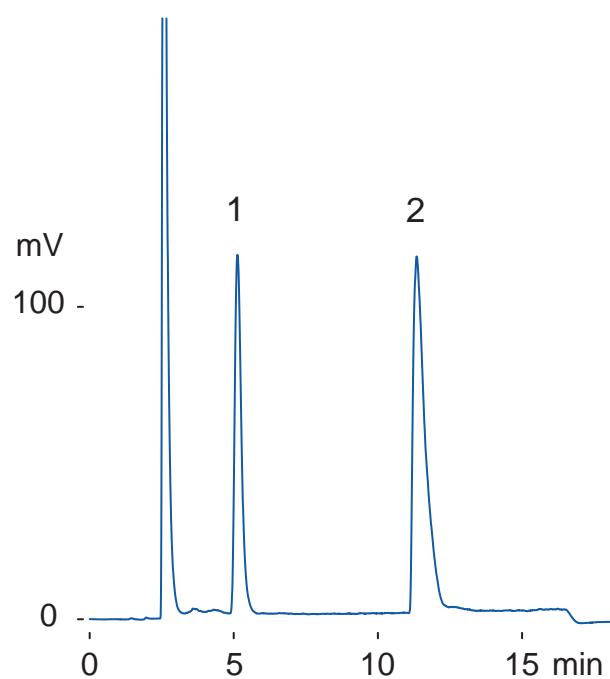

Unison UK-Amino, 100 x 3 mm  
A: acetonitrile  
B: water /triethylamine /acetic acid  
= 100 /4.1 /1.6  
35-80% B(0-15min), 0.4 mL/min (5MPa)  
37 deg.C, ELSD, 1.2 uL (4-8ug)