

Unison UK-Amino

150 x 3 mm




Cadenza CL-C18

250 x 3 mm

Application


Acrylamide and Methacrylamide

アクリルアミド, メタクリルアミド

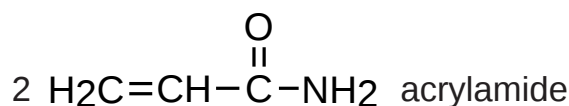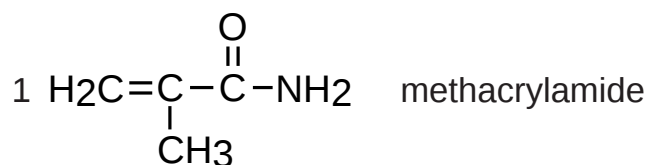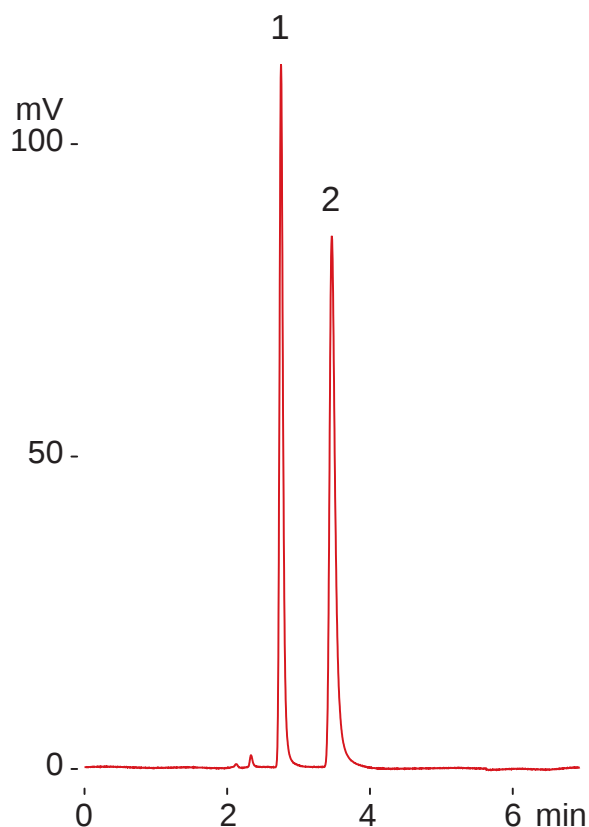

normal phase

Unison UK-Amino, 150 x 3 mm
THF /water /acetic acid = 100 /0.8 / 0.1
0.4 mL/min (5MPa), 37 deg.C, 230nm
1 uL (0.3-0.4ug, dissolved with eluent)

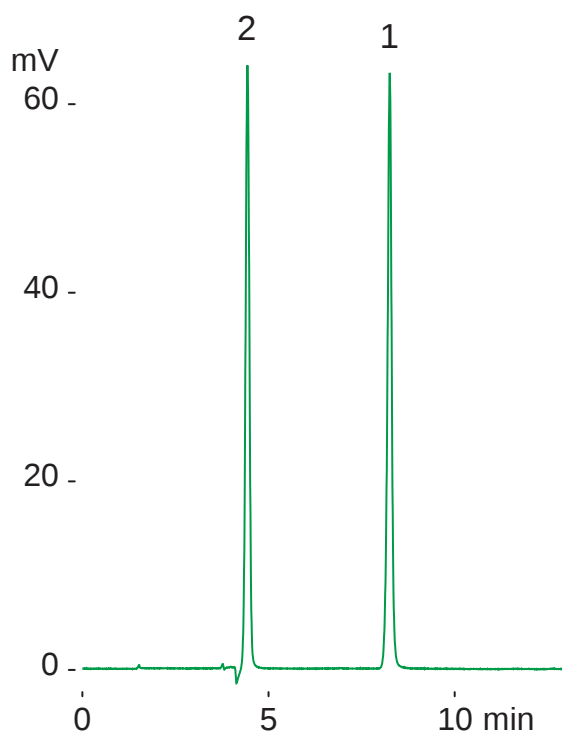

reversed phase

Cadenza CL-C18, 250 x 3 mm
water /acetonitrile /TFA = 97 /3 /0.1
0.4 mL/min (12 MPa), 45 deg.C, 230 nm
1 uL (0.3-0.4ug)