


Unison UK-Amino

250 x 3 mm

Application

Methimazole and potential metabolites

メチマゾールと代謝物

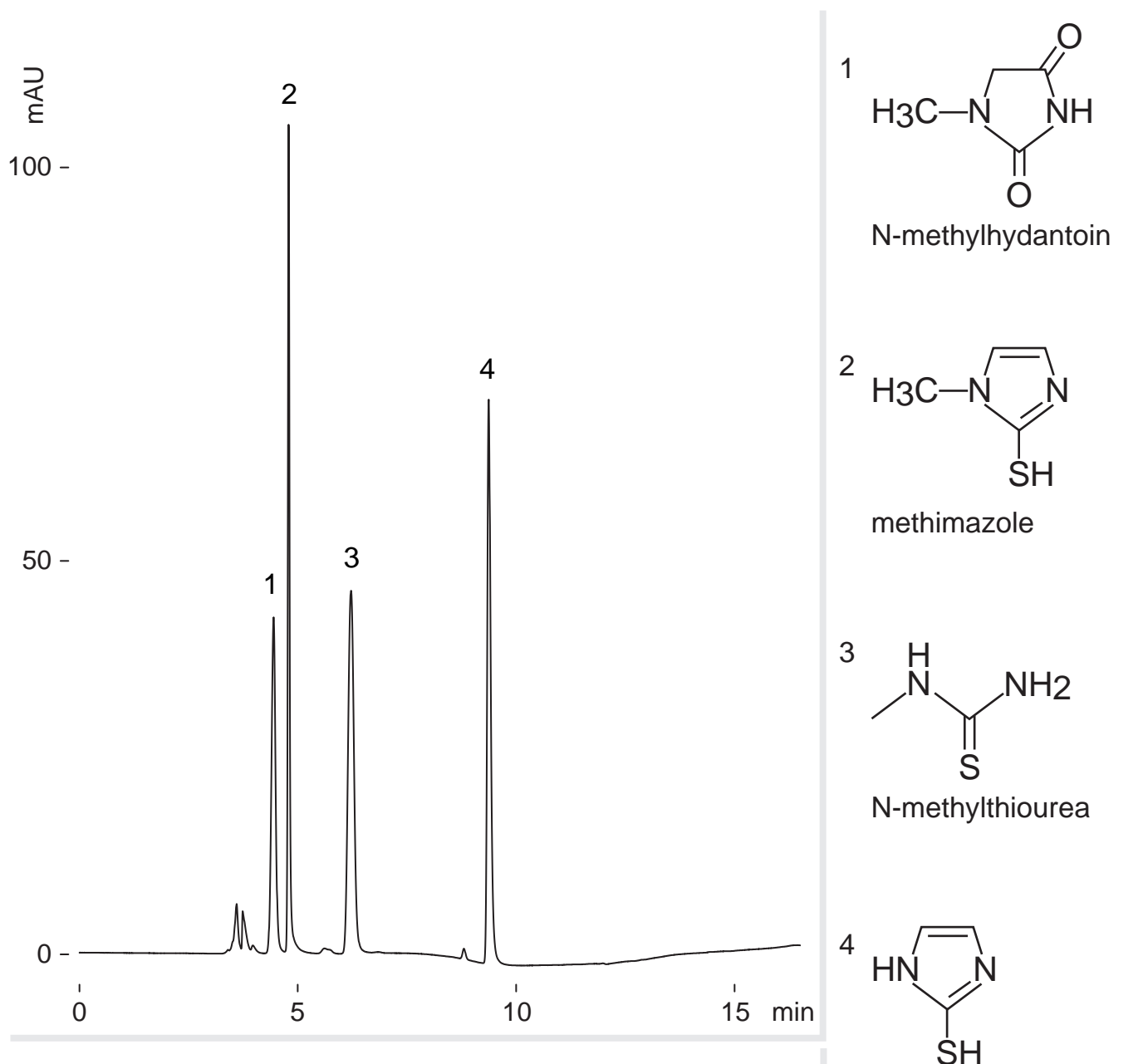

Unison UK-Amino, 250 x 3 mm  
 A: acetonitrile, B: 10mM ammonium formate  
 1-25 %B (0-15min)  
 0.4 mL/min (5MPa), 50 deg.C, 230 nm  
 2.4 uL (0.07-0.35 ug, diluted with ethanol)

Courtesy of Takeo Sakuma, Ph.D, MDS Analytical Technologies (SCIEX), Canada