


Scherzo SM-C18

250 x 3 mm


Application

Herbicide formulation (glyphosate)

除草剤製剤(グリホサート)

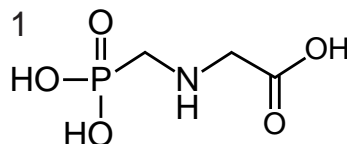

glyphosate

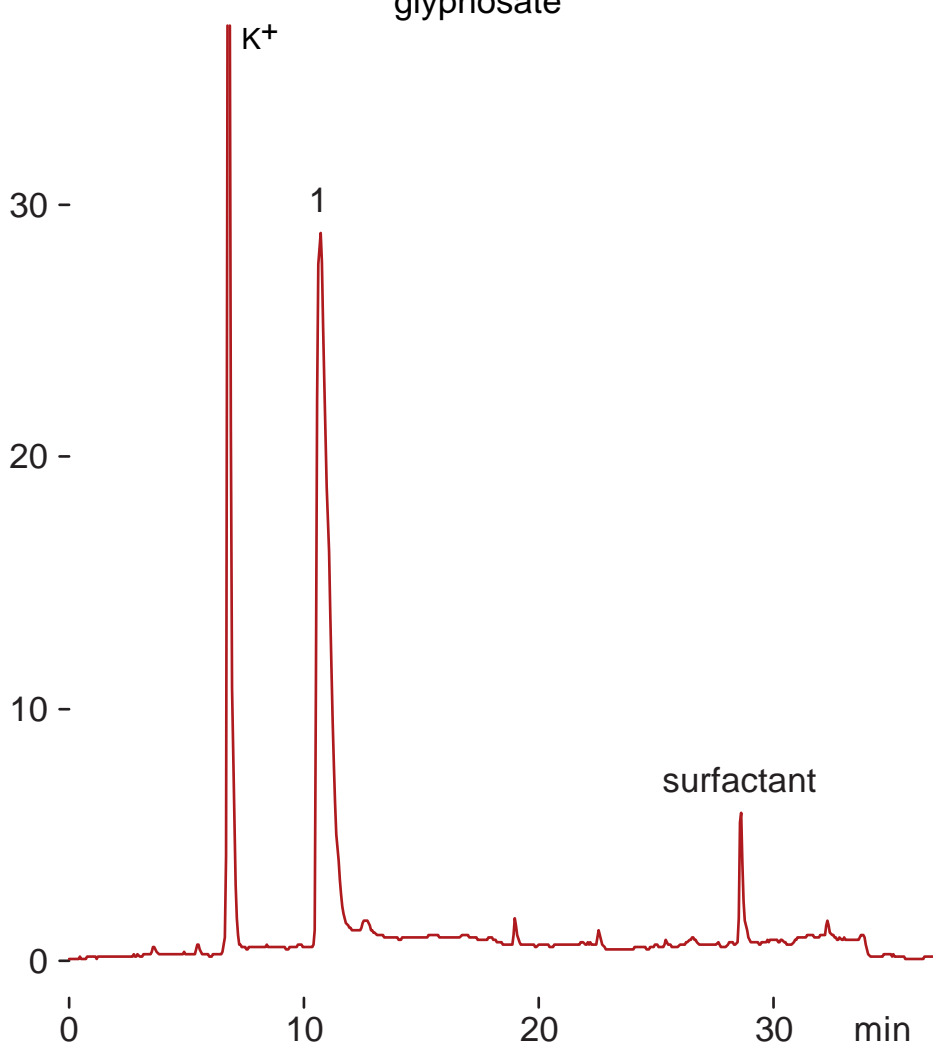

Scherzo SM-C18, 250 x 3 mm

A: 5 mM ammonium formate

B: 100 mM ammonium formate / acetonitrile = 25 / 75

0-100 %B (0-25min), 100 %B(25-30min)

0.4 mL/min (12MPa), 37 deg.C

ELSD (spray chamber 10 deg.C, drift tube 50 deg.C)

5 uL (0.1% with 10mM ammonium formate)