
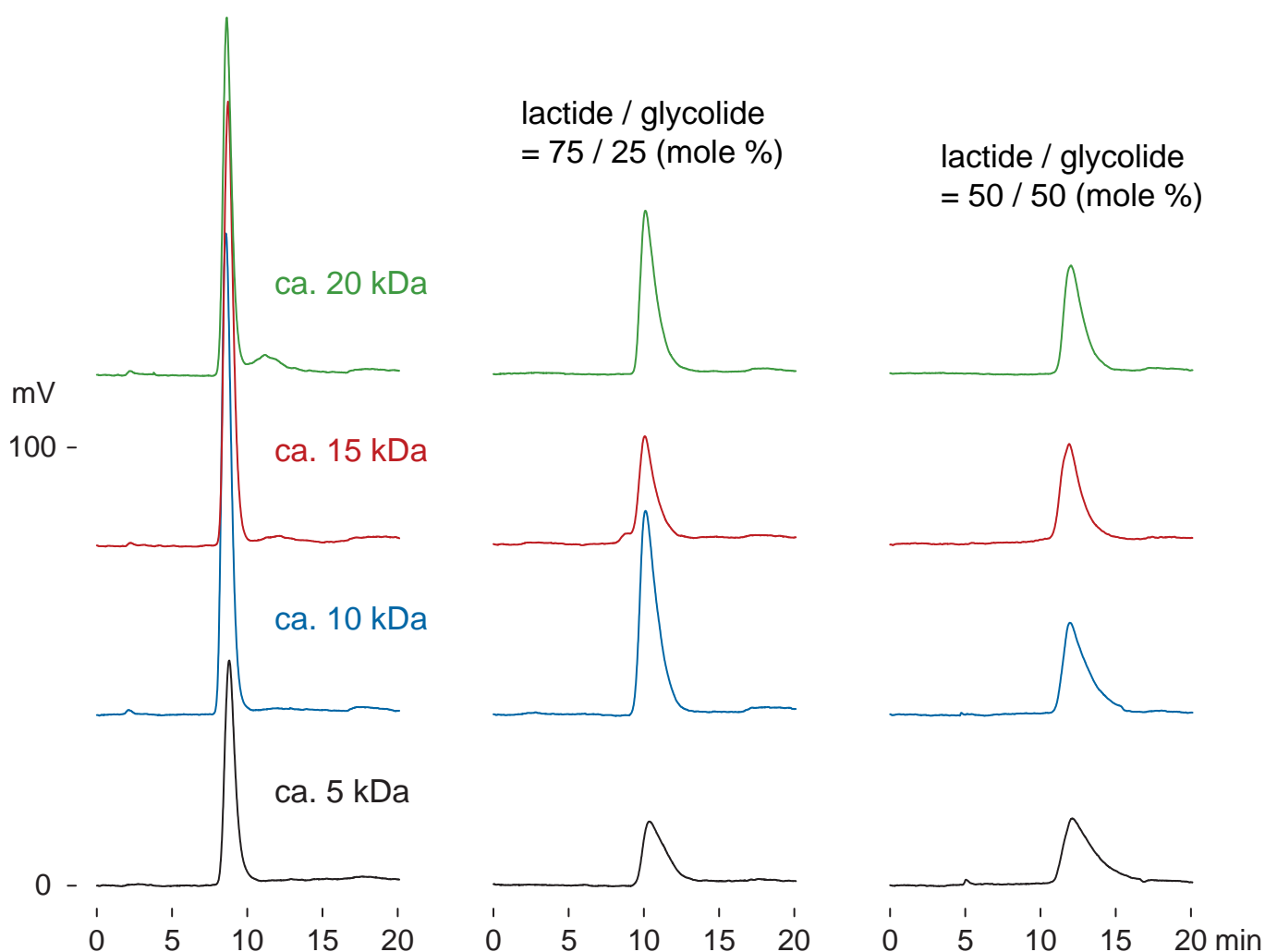
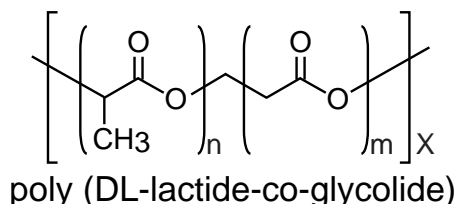


Unison UK-Amino

150 x 3 mm

Application

Biodegradable polymer
(poly (DL-lactide) and poly (DL-lactide-co-glycolide))
 生分解ポリマー (DL-乳酸ポリマー, DL-乳酸 / グリコール酸コポリマー)

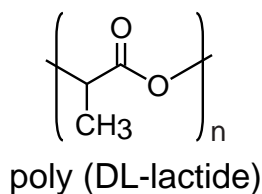

Unison UK-Amino , 150 x 3 mm

A: hexane / acetic acid = 100 / 0.1

B: acetone / acetic acid = 100 / 0.1

30-100 %B (0-10min), 100 %B (10-15min)

0.6 mL/min (4MPa), 60 deg.C

ELSD (spray chamber 50 deg.C, drift tube 100 deg.C), 0.6 uL (6 ug, THF)