


Presto FF-C18

150 x 4.6 mm

Application

**Starch, amylopectin, amylose (from potato)**

デンプン, アミロペクチン, アミロース (馬鈴薯由来)

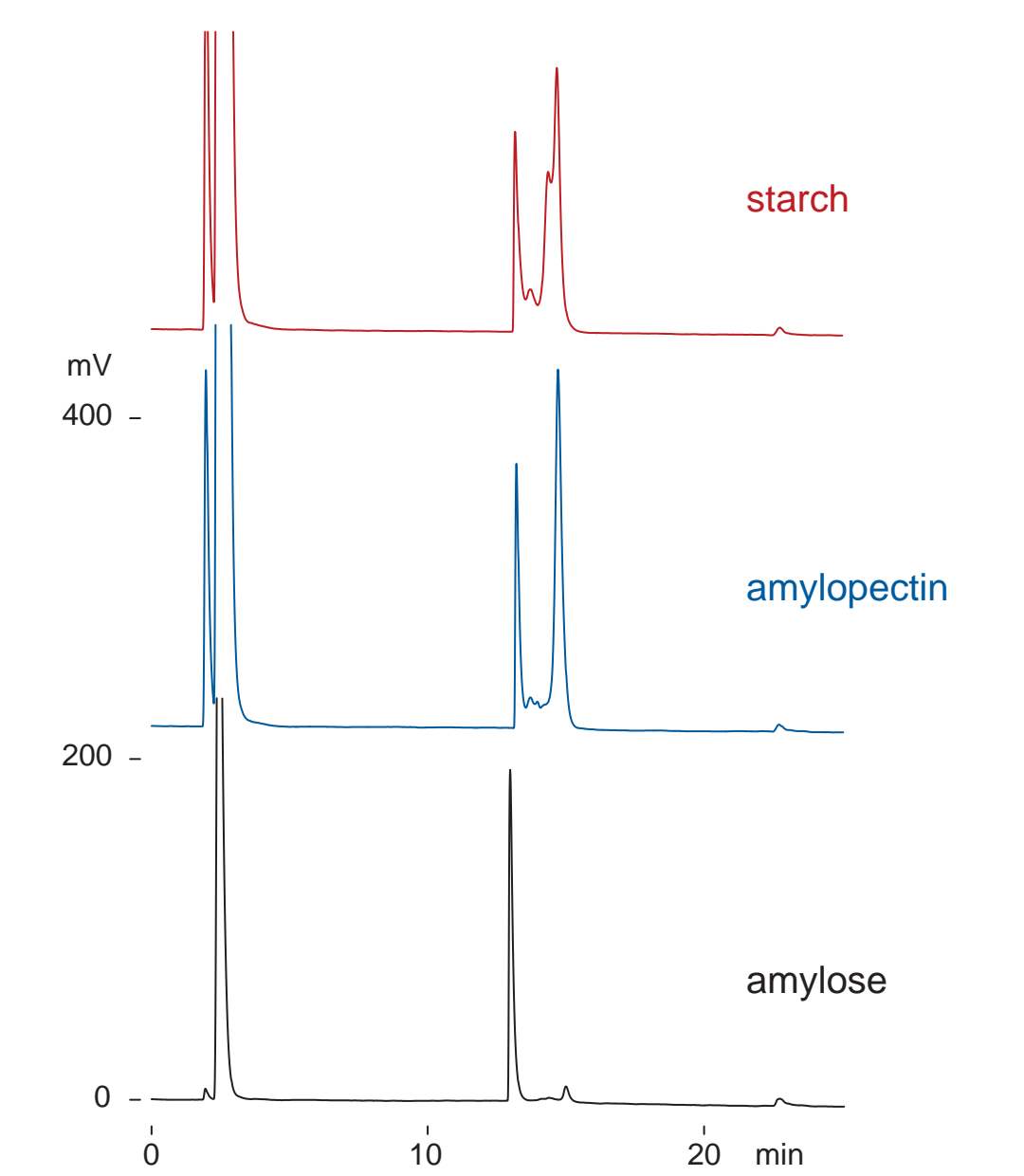

Presto FF-C18, 150 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-7 %B (0-20min), 90 %B (20-30min)

0.5 mL/min (15MPa), 60 deg.C

ELSD (spray chamber 50 deg.C, drift tube 100 deg.C)

1-5 uL (5-25 ug, 0.2 N-NaOH)