


Scherzo SM-C18

150 x 2 mm

Application

LC-MS/MS Application:  
Simultaneous determination of phytotoxins in urine

LC-MS/MSによる尿中植物性自然毒一斉分析

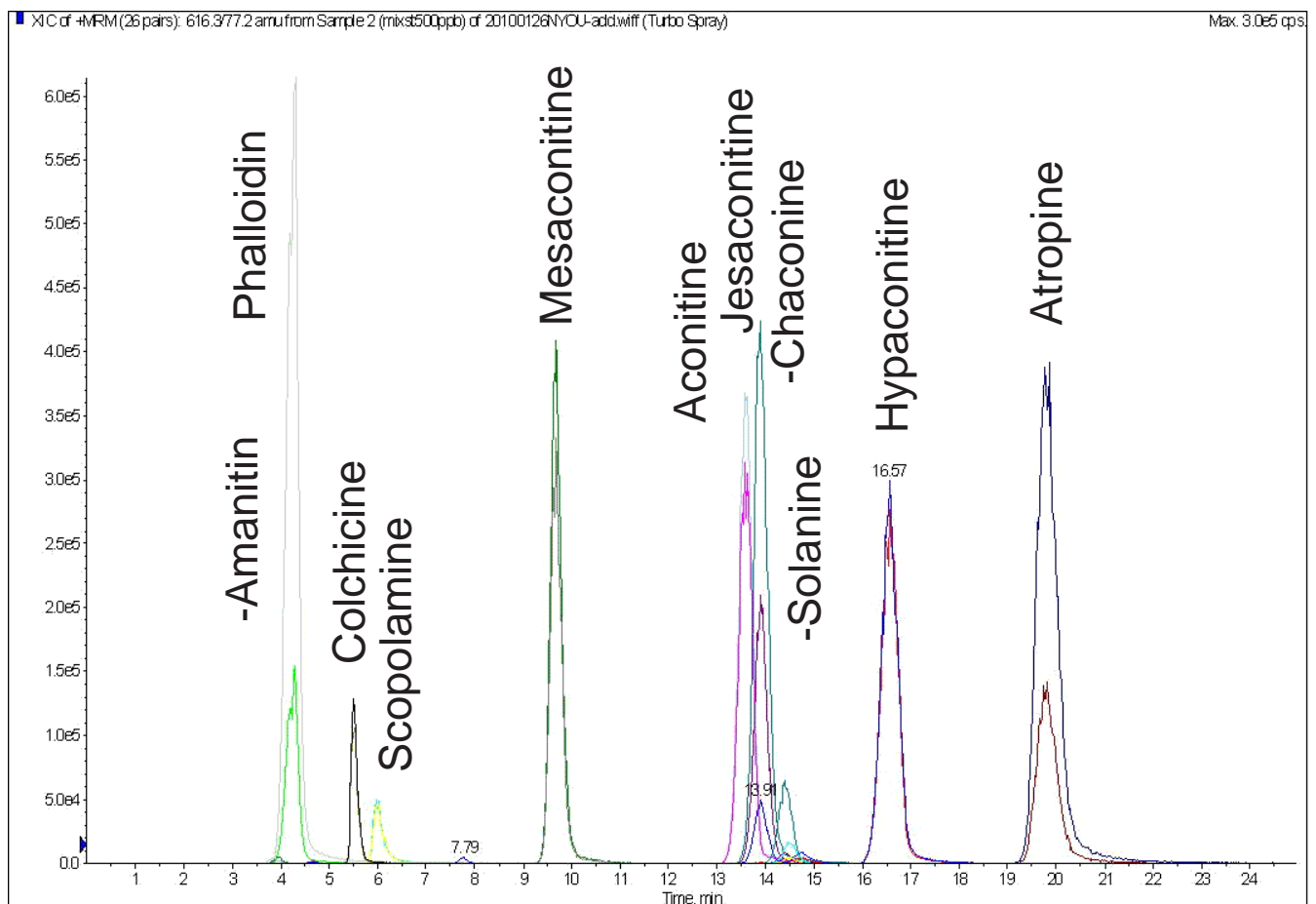

Scherzo SM-C18, 150 x 2 mm  
10mM ammonium formate / methanol = 10 / 90  
0.2 mL/min, 40 deg.C, 10uL  
ESI-MS (MRM, positive)

Ref.)  
Simultaneous Determination of phytotoxins in urine with LC/MS/MS  
Kouji Tachino, Michiko Fujiwara, Izumi Miura  
山口県環境保健センター所報 (ISSN 1882-4536), 第52号 (2010)

Courtesy of TACHINO Kouji, Yamaguchi Prefectural Institute of Public Health and Environment, JAPAN