



Cadenza CD-C18
 Unison UK-C18
 Cadenza CL-C18
 Scherzo SM-C18

150 x 3 mm


Application

Endosulfan isomers (organochlorine pesticide)

有機塩素系殺虫剤・エンドスルファンの異性体

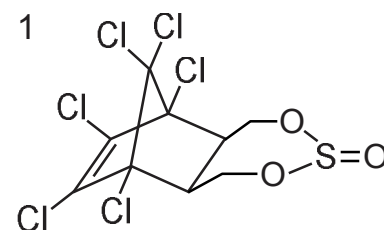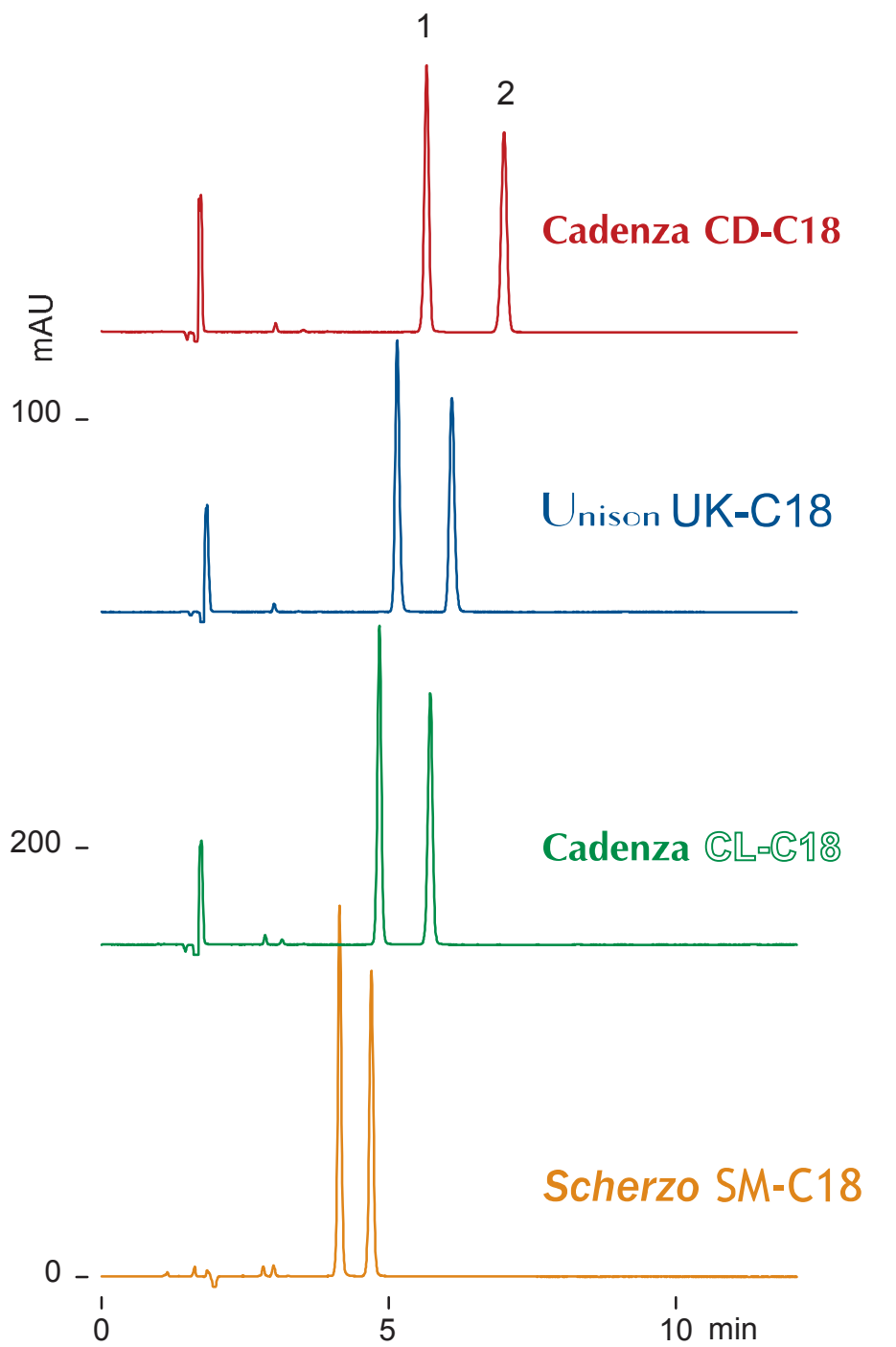

beta-endosulfan

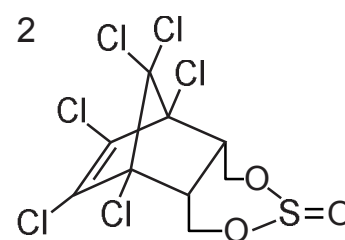

alpha-endosulfan

150 x 3 mm

water / acetonitrile / acetic acid = 20 / 80 / 0.1

0.4 mL/min (5-6 MPa), 37 deg.C, 220 nm, 1.5 uL (0.75 ug, methanol)