




Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18
 Unison UK-C18
 Cadenza CL-C18

150 x 3 mm

Application

Acetaminophen and metabolite

アセトアミノフェンと代謝物(中性条件)

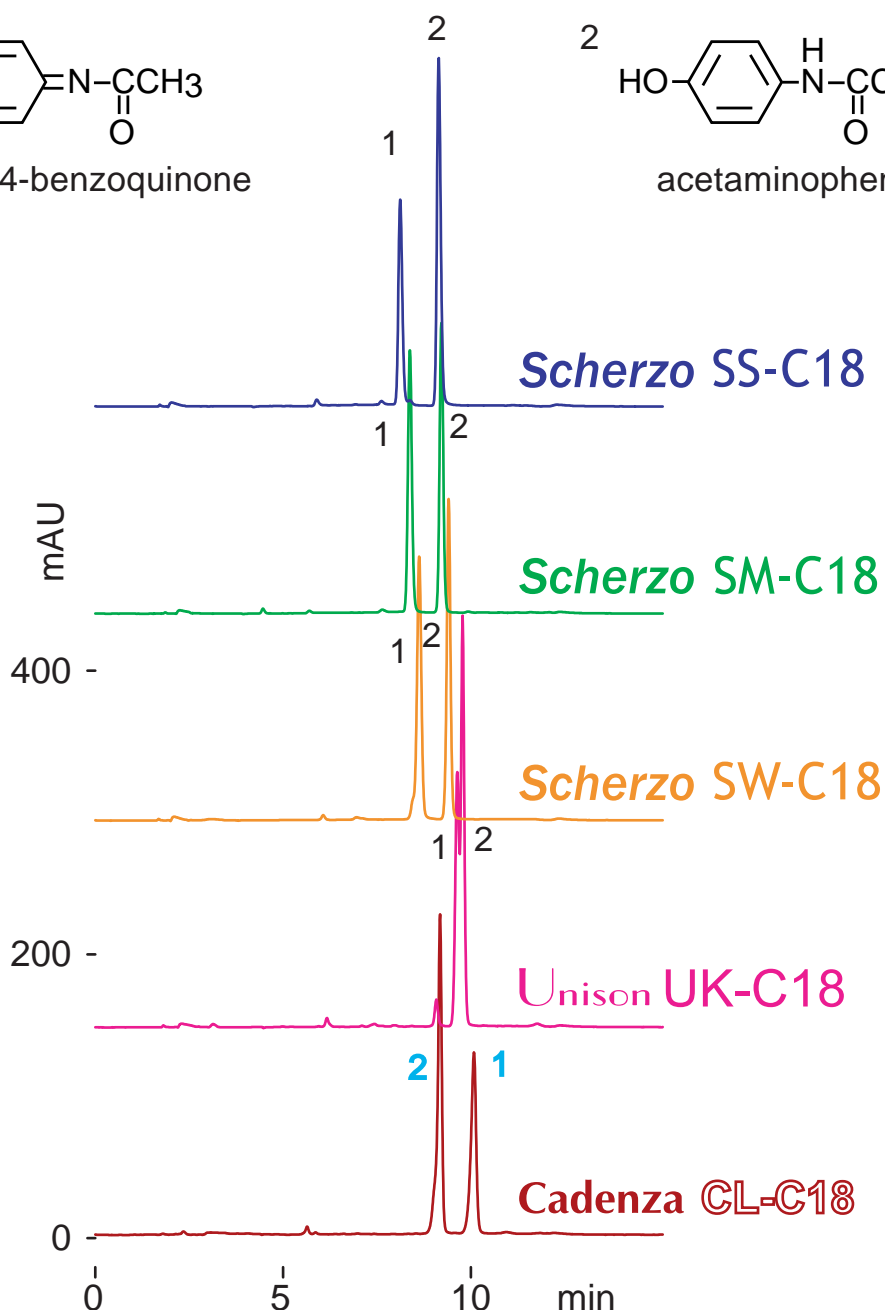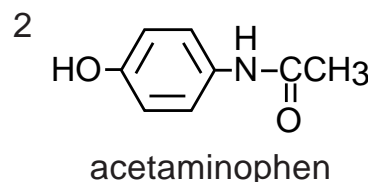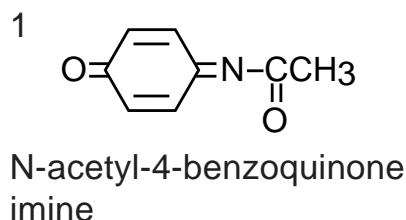

150 x 3 mm

A: 2mM ammonium acetate

B: 100mM ammonium acetate / acetonitrile = 50 / 50

0-20 %B (0-10min)

0.5 mL/min (9-11MPa), 37 deg.C, 260 nm, 1 uL (0.8 ug)