

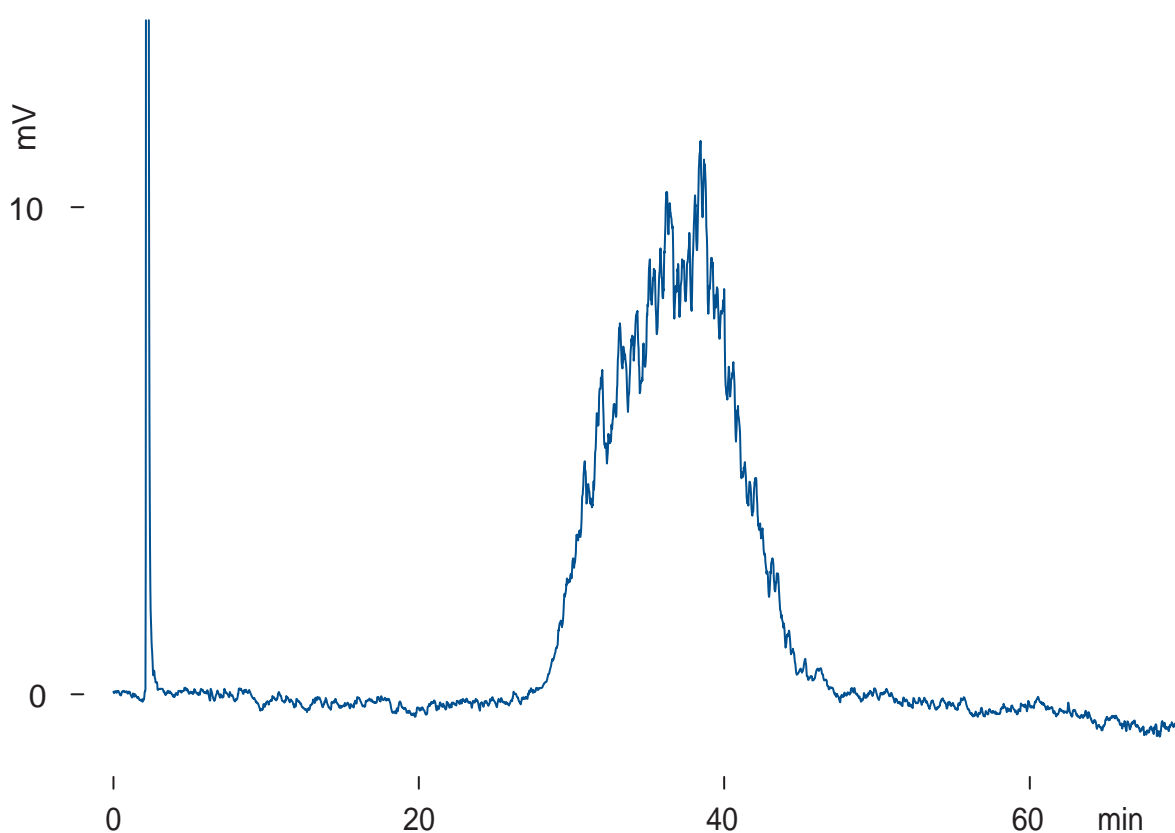
Presto FF-C18

150 x 3 mm

Application

Cytochrome c (50°C, formic acid)

ギ酸移動相によるシトクロム c

**cytochrome c
(horse heart, 12.4 kDa)**

Presto FF-C18, 150 x 3 mm

A: water / formic acid = 100 / 0.5

B: acetonitrile / formic acid = 100 / 0.5

20-30 %B (0-60 min)

0.2 mL/min (17 MPa), 50 deg.C, ELSD (spray chamber 50 deg.C, drift tube 100 deg.C)

5 uL (50 ug)