




Cadenza CD-C18 Unison UK-Phenyl
 Cadenza CL-C18 Unison UK-C8
 Cadenza CW-C18 Unison UK-C18

150 x 3 mm

Application

Prostaglandins : acidic condition

プロスタグランジン : 酸性条件

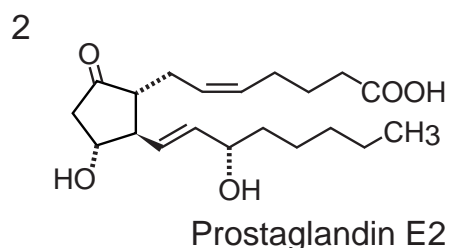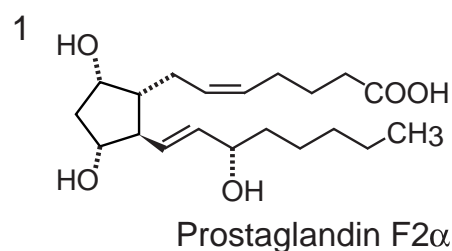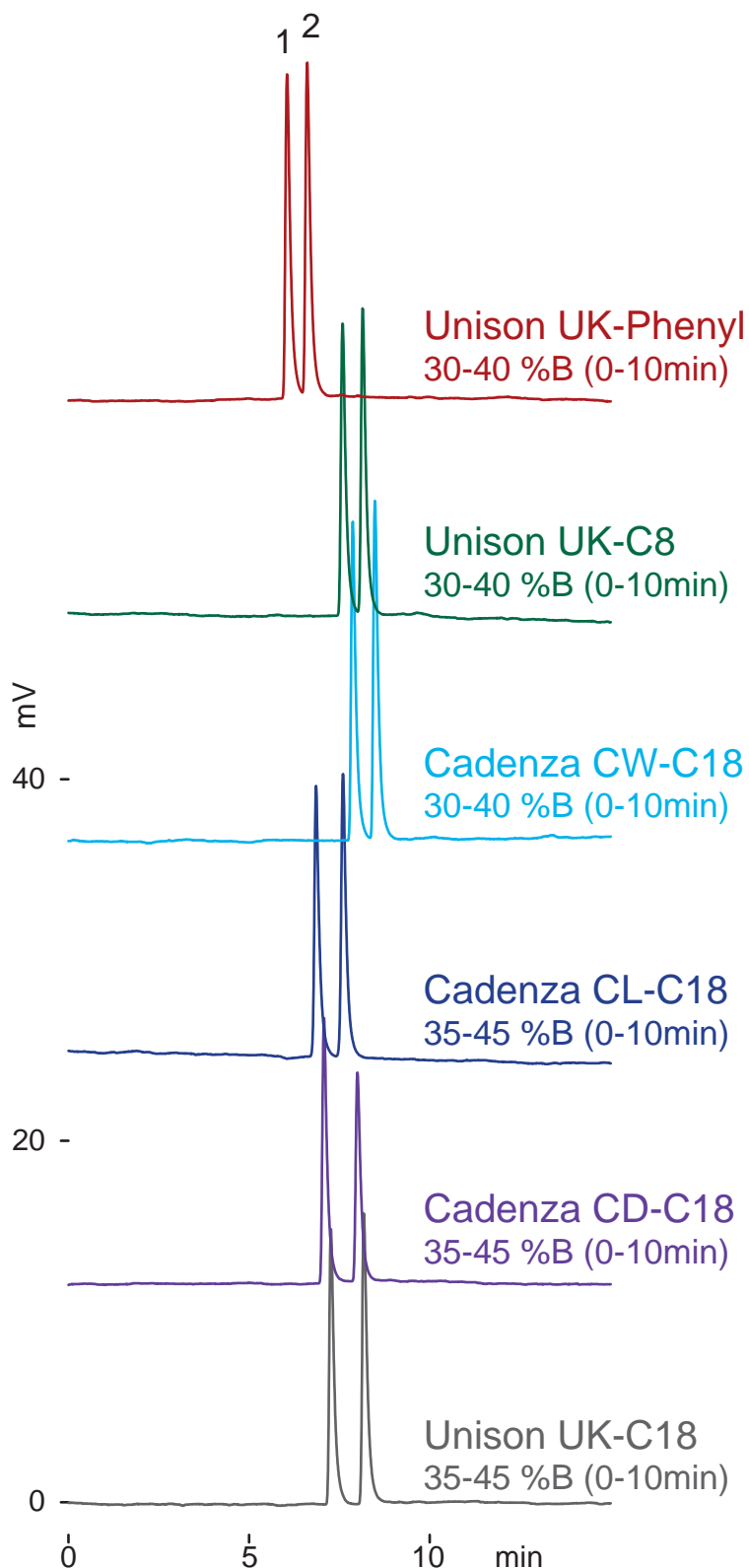

150 x 3 mm
 A: water / HCOOH = 100 / 0.5
 B: acetonitrile / HCOOH = 100 / 0.5
 0.4 mL/min , 37 deg.C, 7-9 MPa
 ELSD(spray chamber 15 deg.C,
 drift tube 50 deg.C)
 1.4 μ L (0.56 μ g, methanol)