




Cadenza CD-C18 Unison UK-Phenyl  
 Cadenza CL-C18 Unison UK-C8  
 Cadenza CW-C18 Unison UK-C18

150 x 3 mm

Application

Prostaglandins : neutral pH condition

プロスタグランジン : 中性条件

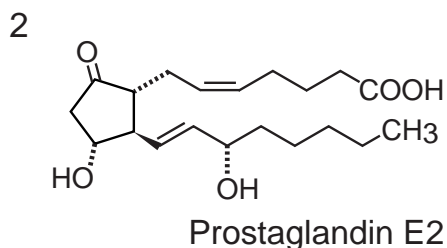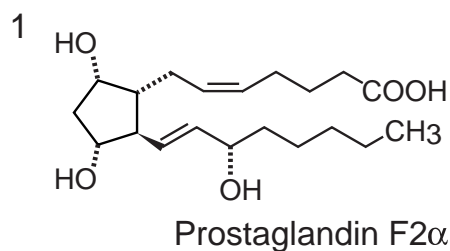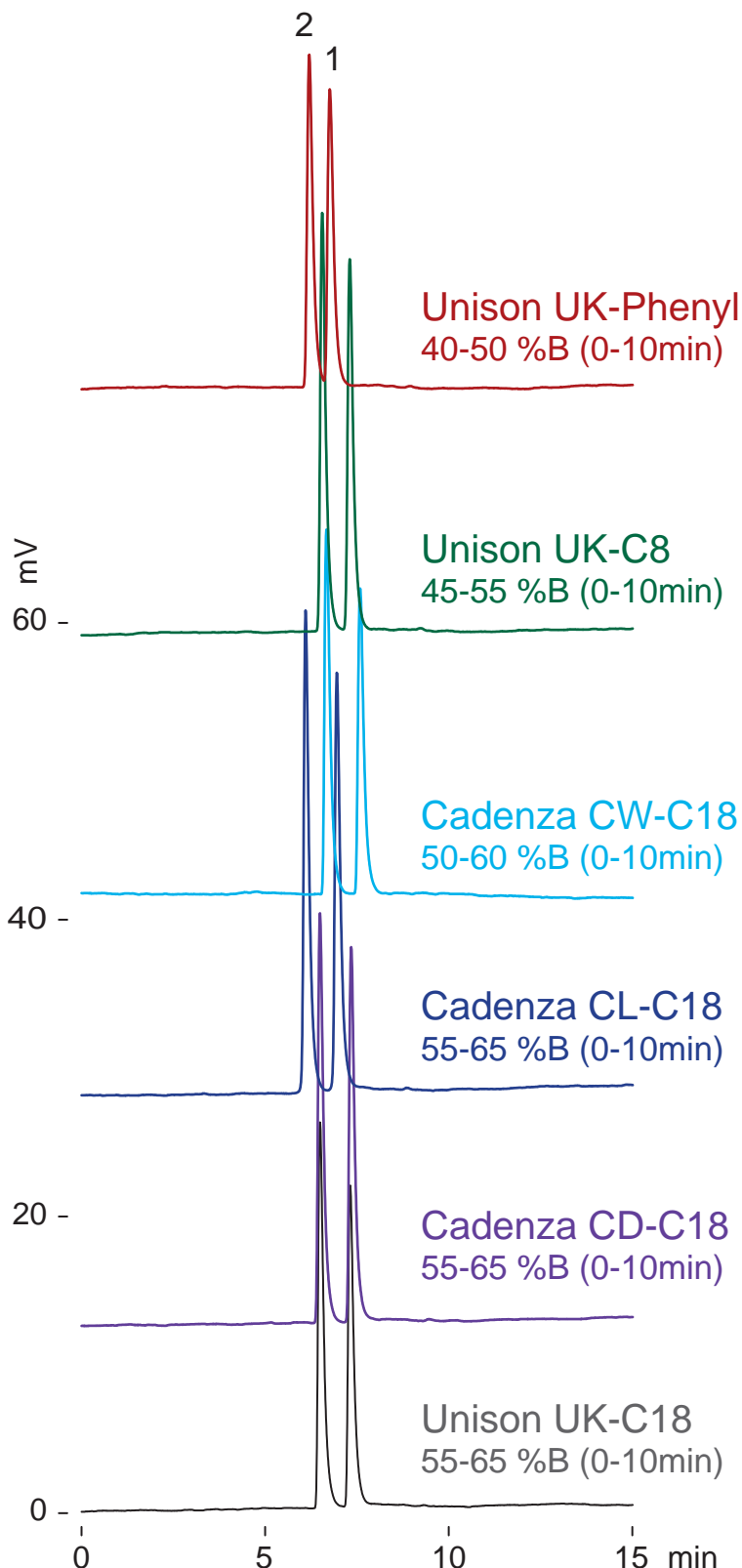

150 x 3 mm  
 A: 50 mM ammonium acetate  
 B: methanol  
 0.4 mL/min (11-14 MPa), 37 deg.C  
 ELSD(spray chamber 15 deg.C,  
 drift tube 50 deg.C)  
 1.4 uL (0.56 ug, methanol)