


**Scherzo SW-C18**  
**Scherzo SM-C18**  
**Unison UK-C18**

150 x 3 mm

Application

**Oligopeptides (angiotensins)**

オリゴペプチド (アンジオテンシン類)

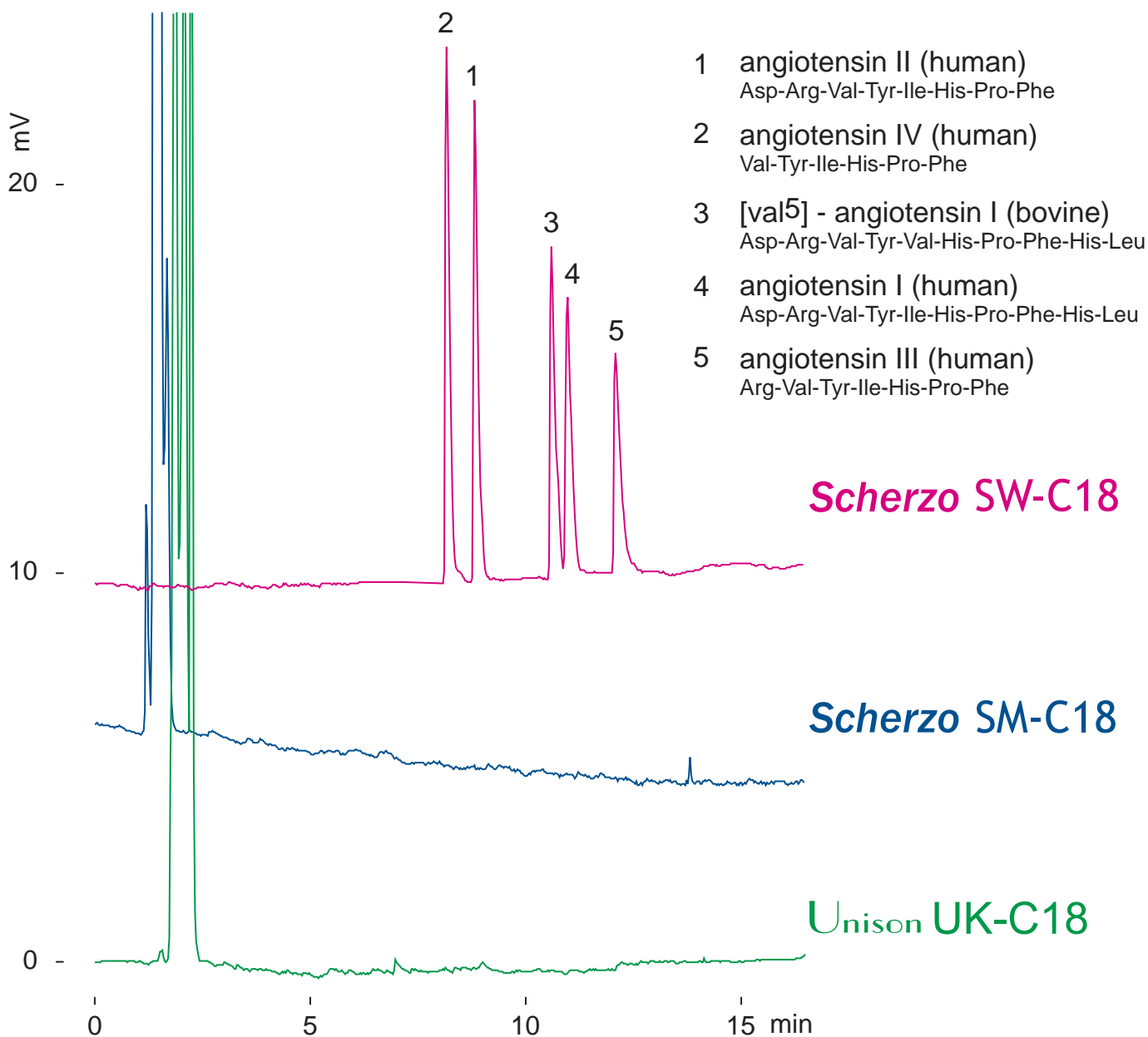

**Scherzo SW-C18**

**Scherzo SM-C18**

**Unison UK-C18**

150 x 3 mm

A: water / acetonitrile / formic acid = 70 / 30 / 0.2

B: 100 mM ammonium formate / acetonitrile = 40 / 60

0-100 %B (0-15 min)

0.4 mL/min (10-11 MPa), 37 deg.C, ELSD (spray chamber 50 deg.C drift tube 100 deg.C)

5 uL (1 ug)