

Unison UK-Phenyl


150 x 2 mm

Application

### LC-MS: Anxiolytics and metabolites


LC-MS: 抗不安薬(エチゾラム, トリアゾラム)と代謝物

**8OH-etizolam**

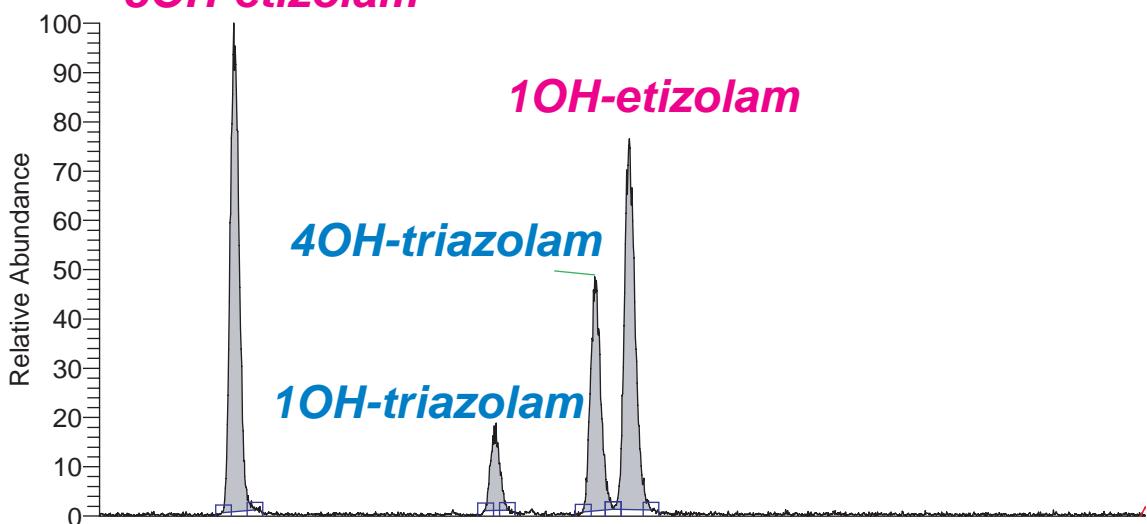

**triazolam**

**etizolam**

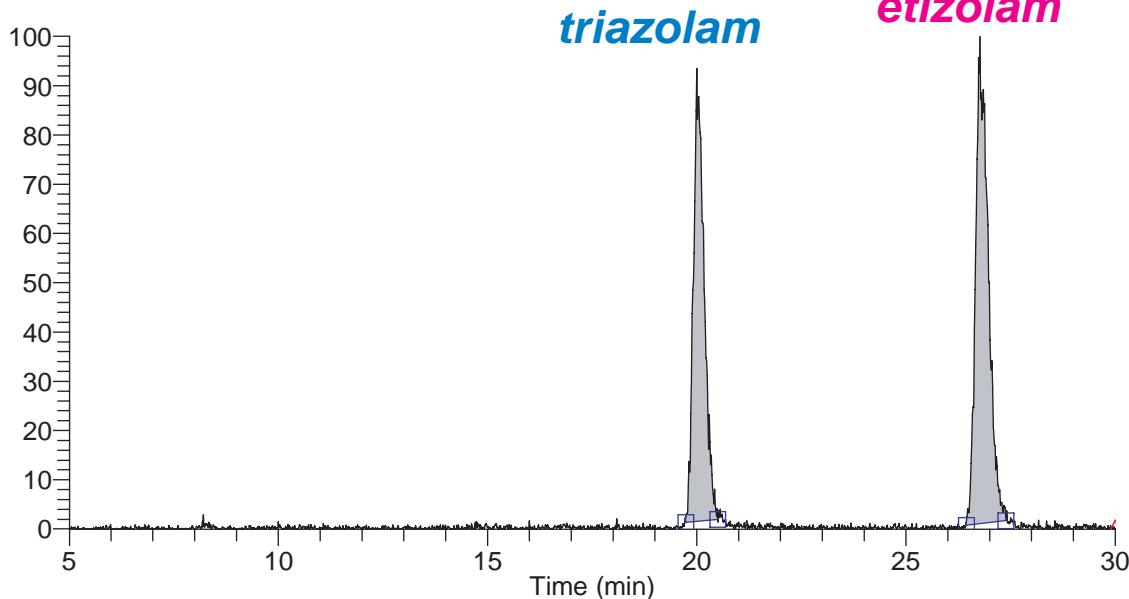

Unison UK-Phenyl, 150 x 2 mm

A: 0.1% HCOOH in water

B: 0.1% HCOOH in methanol

40-50% B (0-20 min), 50% B (20-30min)

0.25 mL/min (15-16 MPa), 40 deg.C, 10uL (100ng/mL, water/ACN=1/1)

MS: ESI, positive

Courtesy of Dr. Koji Yamaguchi, Nippon Medical School, JAPAN