


Intrada Amino Acid

150 x 3 mm

Application

Batch-to-Batch Reproducibility (Intrada Amino Acid)

バッチ間再現性 (Intrada Amino Acid)

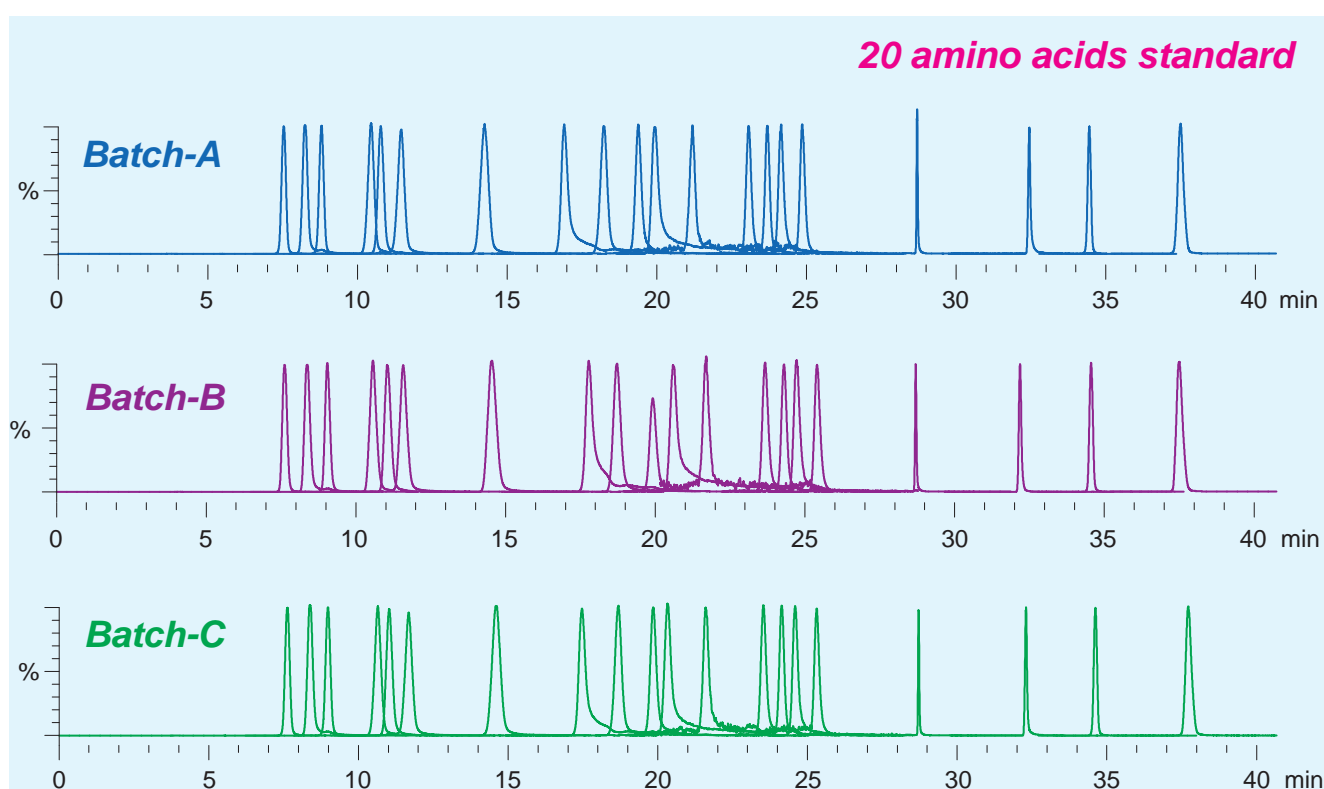

Intrada Amino Acid, 150 x 3 mm

A: ACN /THF /25mM HCOONH₄ /HCOOH = 9 / 75 / 16 / 0.3 (v/v/v/v)

B: ACN /100mM HCOONH₄ = 20 / 80 (v/v)

0 %B (0-12 min), 0-17 %B (12-26 min), 100 %B (26-45 min)

0.4 mL/min (7MPa), 35 deg.C, 1.2 uL (1 umol/mL in 0.1N HCl)

ESI, positive