



Intrada Amino Acid

100 x 3 mm

Application

LC-MS: 20 Amino acids and 5 isotopes (standard)

20種アミノ酸(標準試料)と安定同位体5種のLC-MS/MS分析

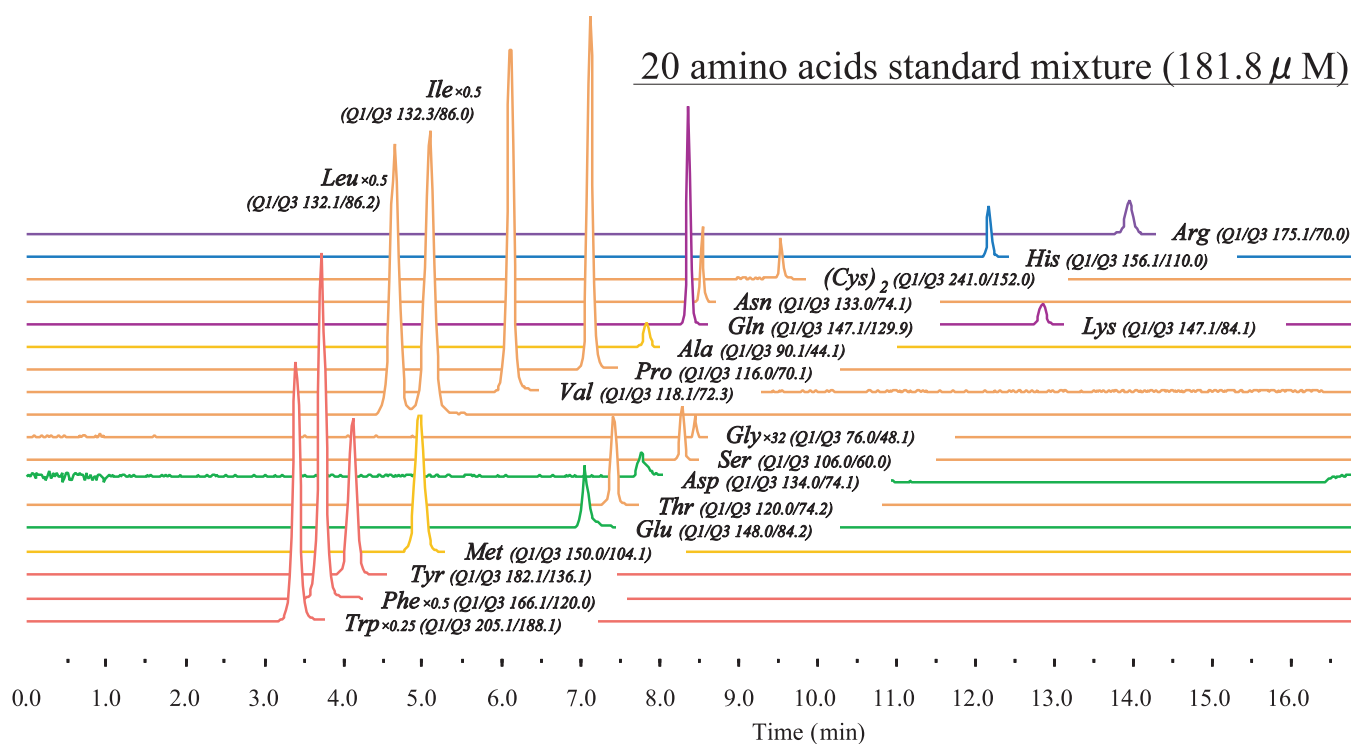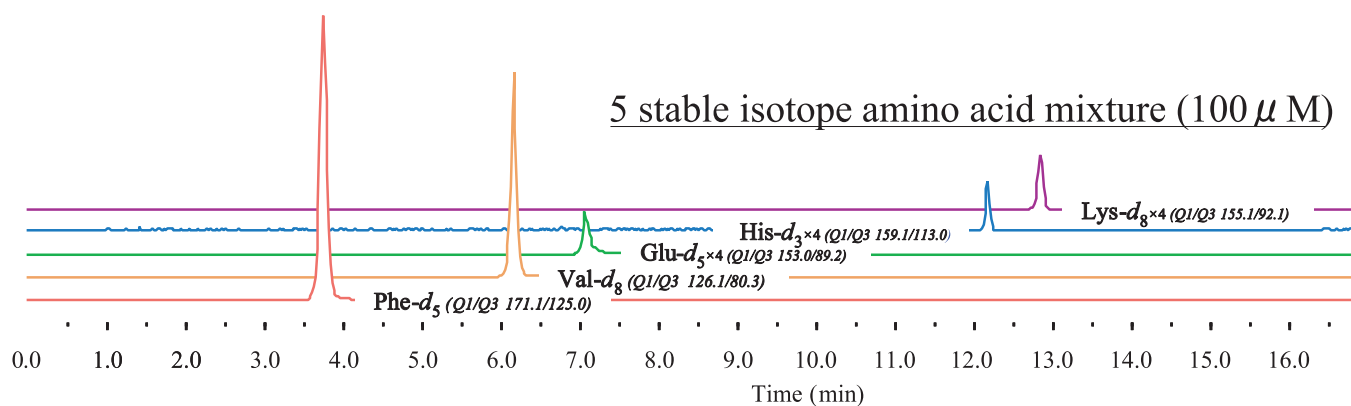

Intrada Amino Acid, 100 x 3 mm

A: acetonitrile /tetrahydrofuran /25mM ammonium formate /formic acid = 9 /75 /16 /0.3

B: acetonitrile /100mM ammonium formate = 20 / 80

0 %B (0-3min), 0-17 %B (3-9 min), 100 %B (9-16 min)

0.6 mL/min, 40 °C, 1μL

INSTRUMENT:

Shimadzu Prominence UFLC System - AB Sciex API 4000 Qtrap MS (ESI, MRM, positive)

Courtesy of Dr.TODOROKI Kenichiro, University of Shizuoka, JAPAN

Miss TOMITA Ryoko, Fukuoka University, JAPAN