

Unison UK-Phenyl






Unison UK-C18

100 x 3 mm

Application

LC-MS: antimigraines (selective serotonin agonists)

トリプタン系頭痛薬のLC-MS分析

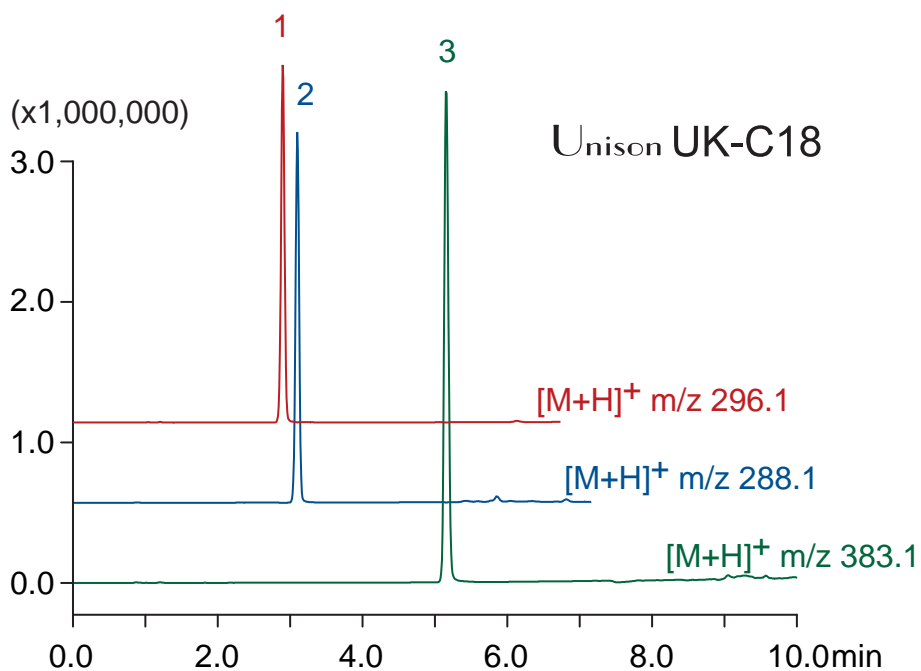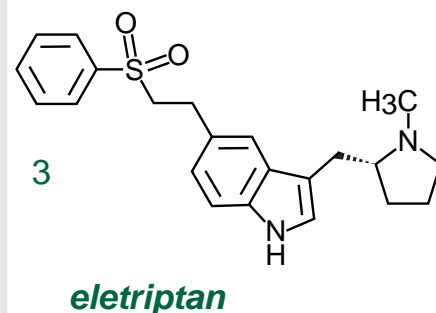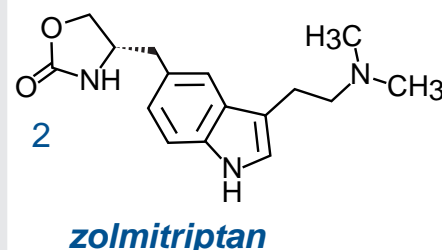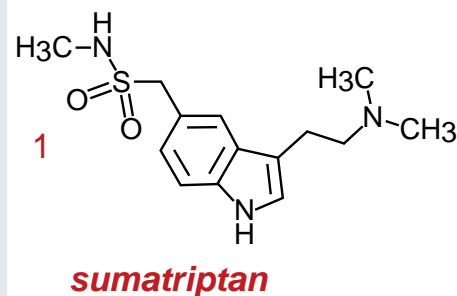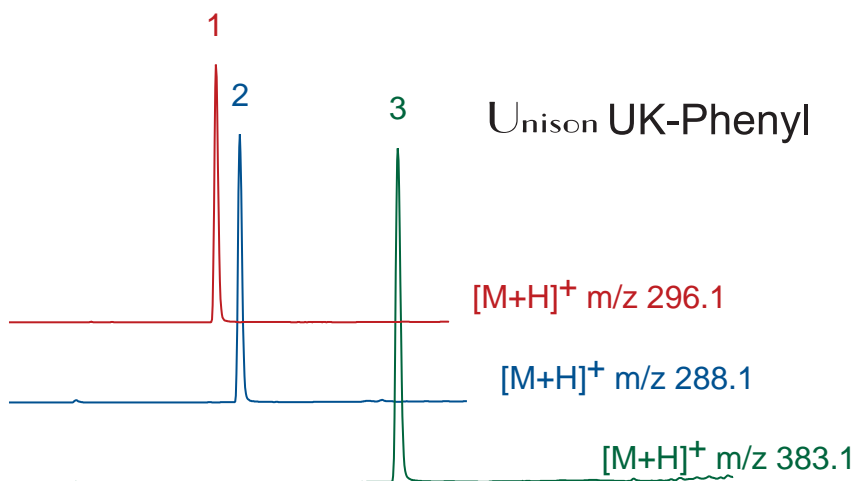

100 x 3 mm

A: water / formic acid = 100 / 0.1

B: methanol / formic acid = 100 / 0.1

10-90 %B (0-8 min), 90 %B (8-10 min)

0.4 mL/min (7 MPa), 37°C, 4 µL (1µg/mL, 0.1% HCOOH)

Single Quad. MS (ESI)