

Cadenza CW-C18

Unison UK-C8

Unison UK-Phenyl

Scherzo SM-C18








Unison UK-C18

50 x 3 mm

Application

LC-MS: Thyroid hormones

甲状腺ホルモン類

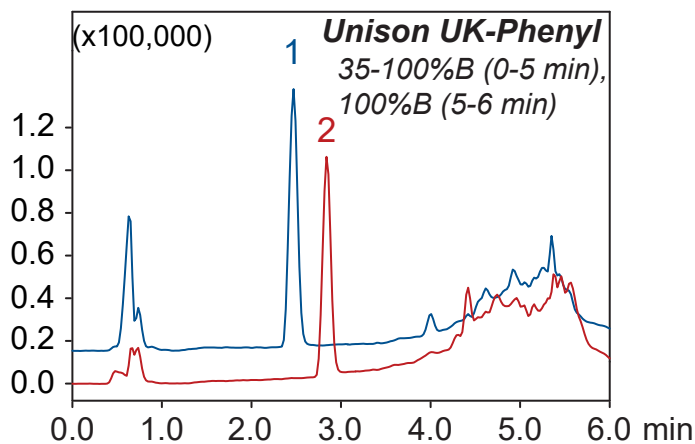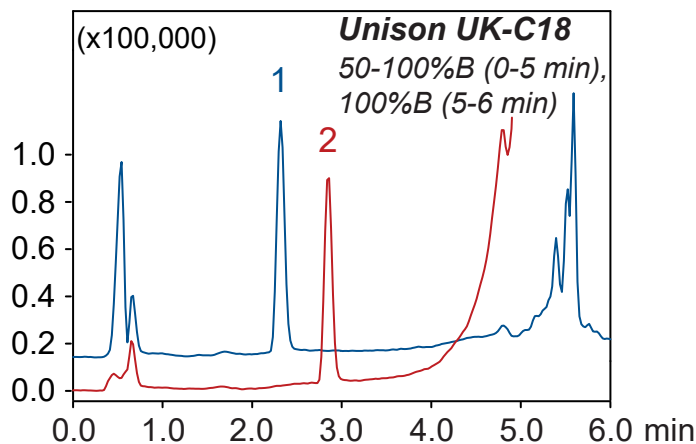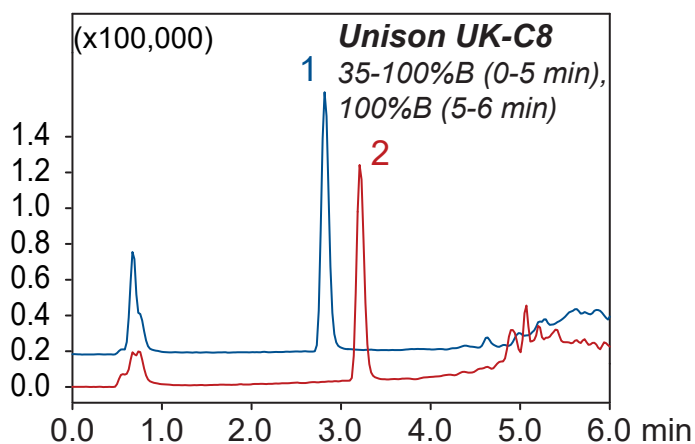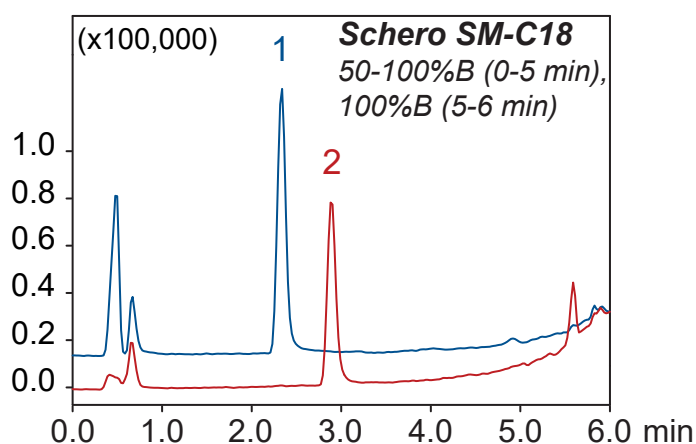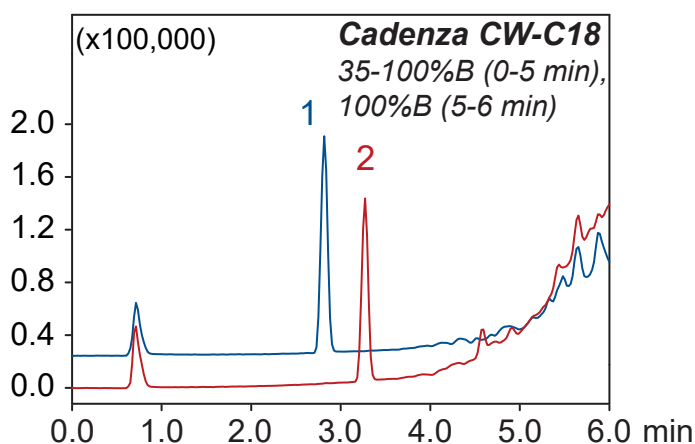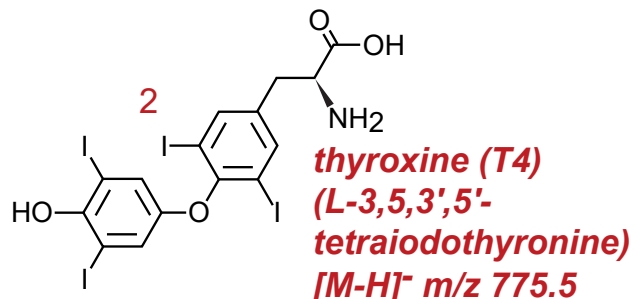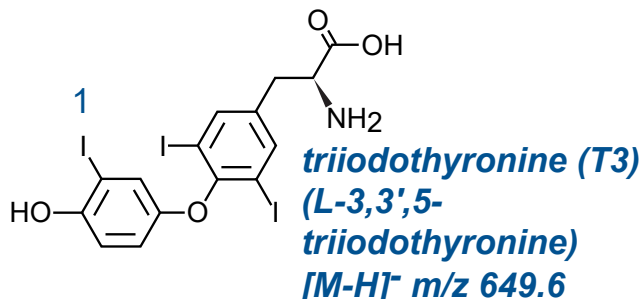

50 x 3 mm  
A: water / AcOH = 100 / 0.1  
B: MeOH / AcOH = 100 / 0.1  
0.4 mL/min (6-7 MPa), 37°C  
20 µL (50ng/mL)  
MeOH / water / AcOH = 50 / 50 / 0.1)  
Single Quad. MS (ESI)