


Cadenza CD-C18

150 x 3 mm

Application

**Ginsenoside**

高麗人參の生理活性物質 ジンセノサイド

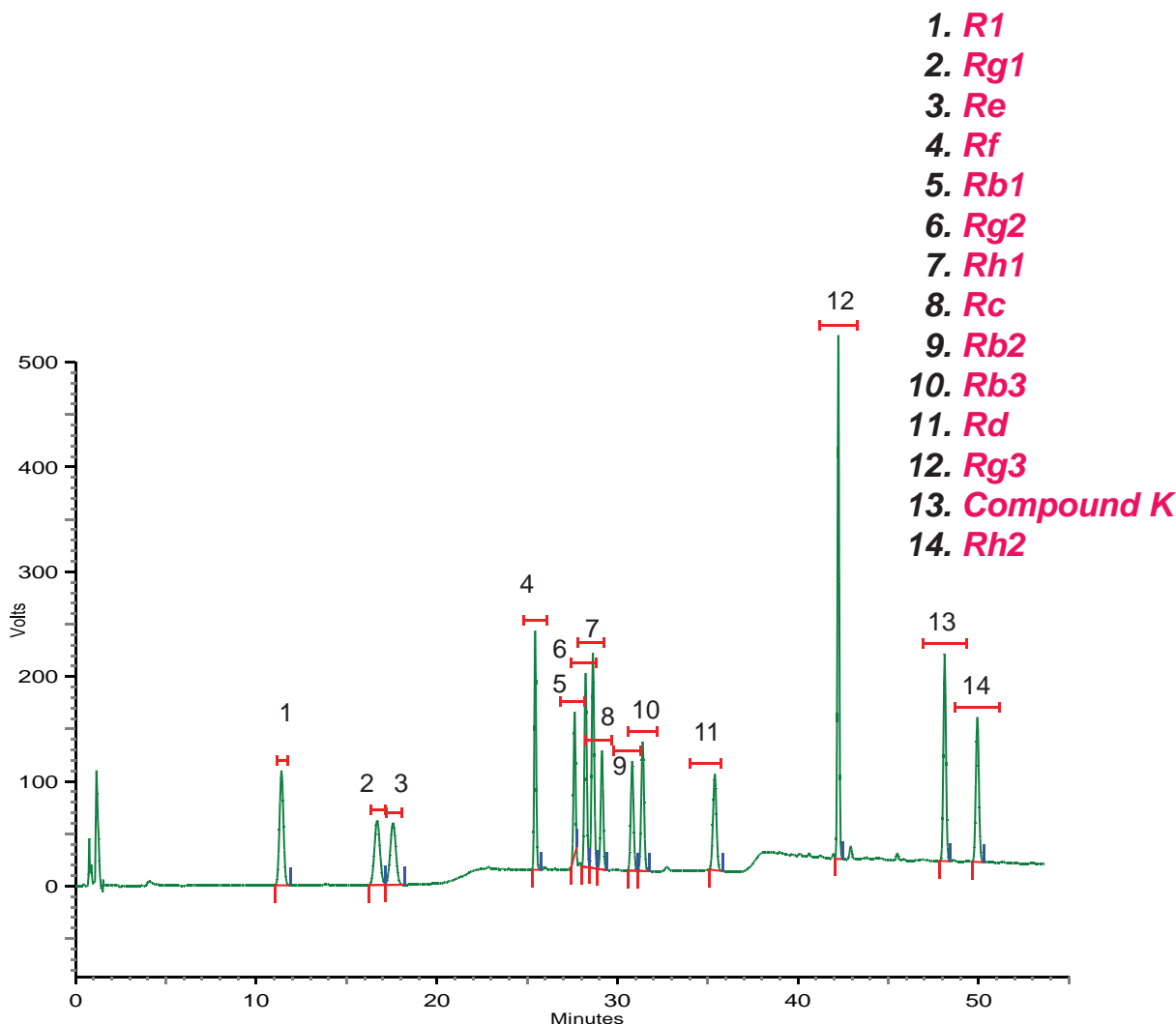

Cadenza CD-C18, 150 x 3 mm

A: 10% acetonitrile

B: 90% acetonitrile

13%B (0-15min), 13-15%B (15-18min), 15-27%B (18-20min),

27-31%B (20-35min), 31-45%B (35-36min), 45-57%B (36-45min),

57%B (45-52min), 57-100%B (52-55min), 100%B (55-65min),

100-13%B (65-66min), 13%B (66-80min)

0.6 mL/min, 40°C

5 µL (50 µg/mL, methanol)

203nm

Courtesy of PEAKMAN SP, KOREA