


Intrada SEC

250 x 4.6 mm

Application


SEC: Proteoglycan

プロテオグリカンのSEC分離

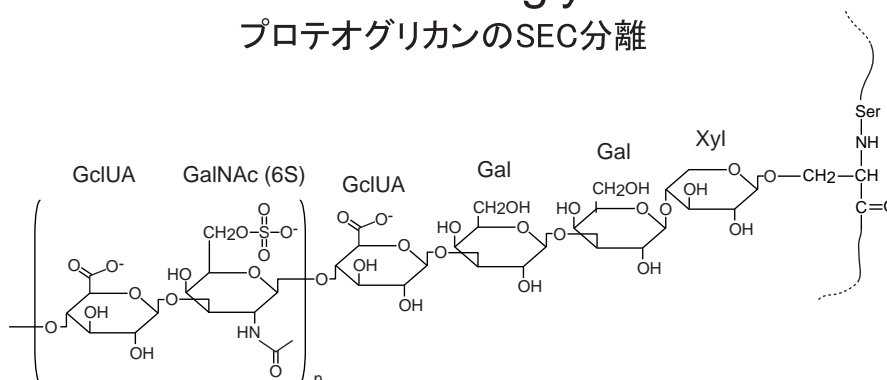

proteoglycan

(from nasal cartilage of salmon head, ca.450 kDa)

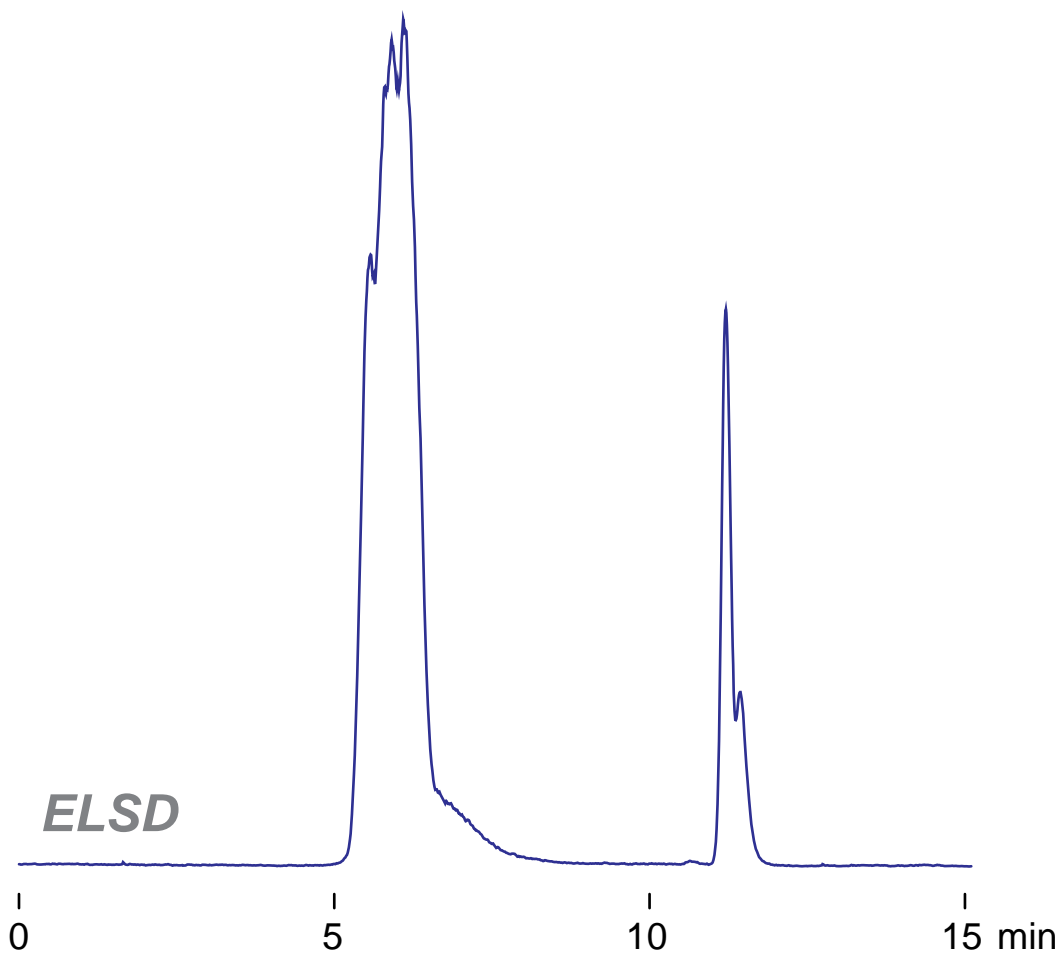

*Intrada SEC, 250 x 4.6 mm
 100mM ammonium formate
 0.3 mL/min (5MPa), 37°C
 10 µL (5mg/mL, eluent)*

Courtesy of Mr. M. Tsuboi, ICHIMARU PHARCOS Co., Ltd. & KAKUHIRO Co., Ltd., Japan