



Scherzo SM-C18 MF

100 x 2 mm


Application

Metal-Free Column: phosphoserine

メタルフリーカラム: ホスホセリン

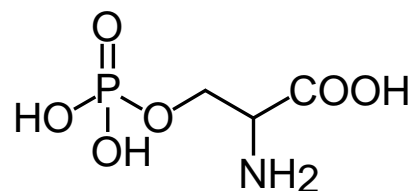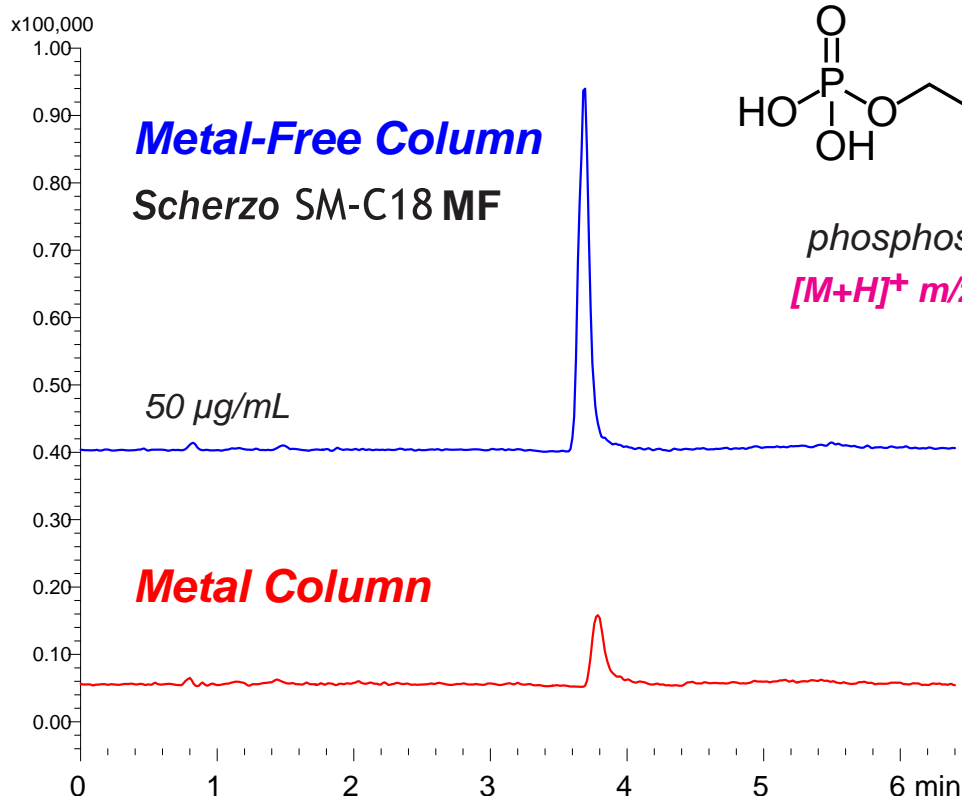

phosphoserine
[M+H]⁺ m/z 186.0

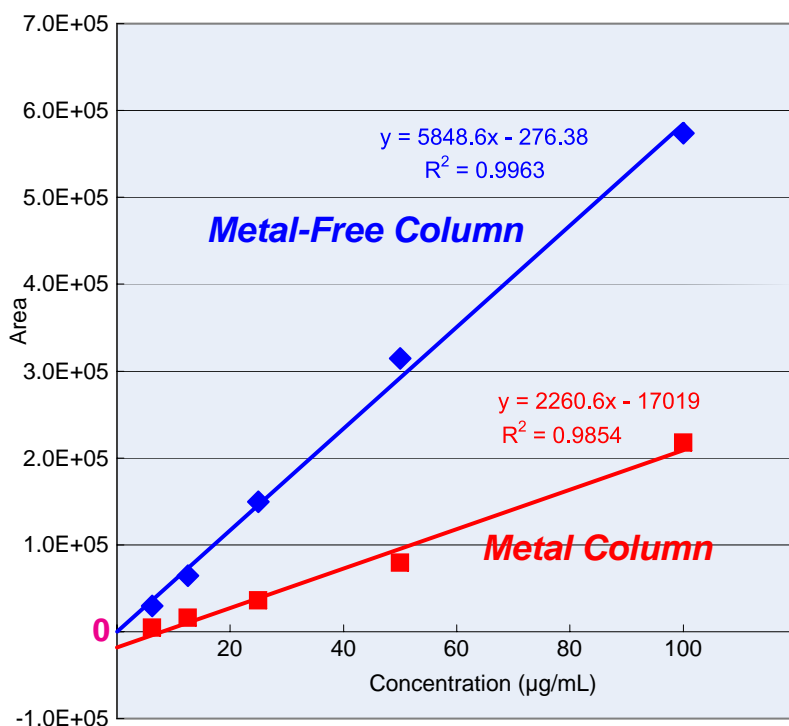

Scherzo SM-C18 MF, 100 x 2 mm

A: water / HCOOH = 100 / 0.1
B: ACN / 100mM HCOONH4 = 20 / 80

0%B (0-2min)
0-100%B (2-4 min)
100%B (4-7min)

0.2 mL/min (6MPa), 37°C
5 µL (1% HCOOH)

Single Quad. MS (ESI, positive)