


Unison UK-Phenyl

100 x 3 mm


Application

LC-MS: Carbonic anhydrase inhibitor (antiglaucoma drug)


炭酸脱水酵素阻害剤 (緑内障治療薬)

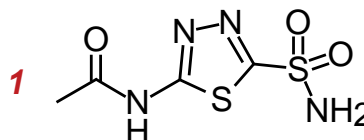

**acetazolamide**

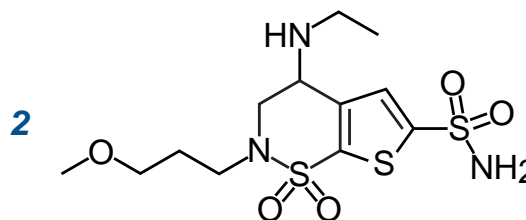

**brinzolamide**

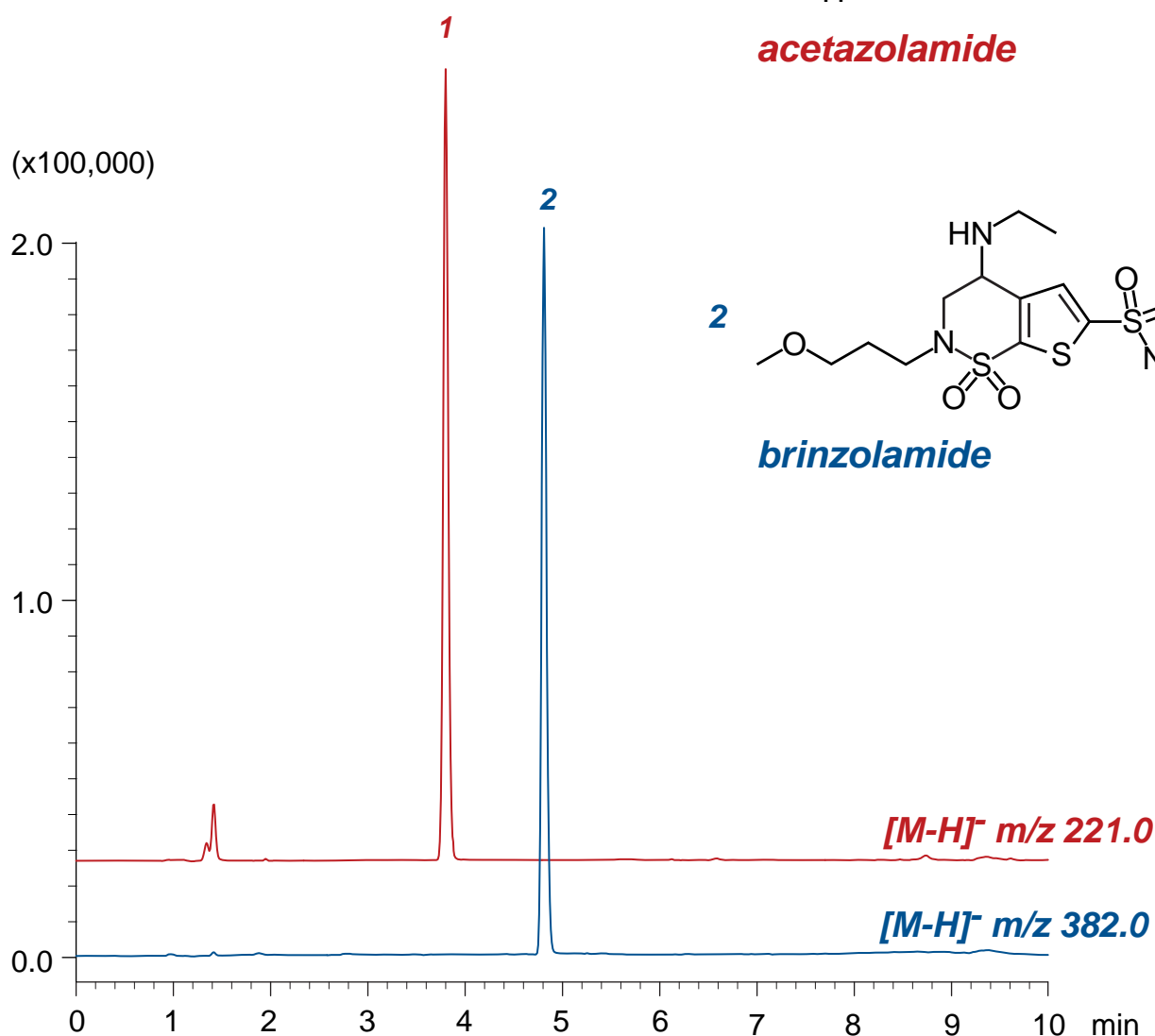

Unison UK-Phenyl , 100 x 3 mm  
 A: water/ formic acid = 100 / 0.1  
 B: acetonitrile / formic acid = 100 / 0.1  
 0-50%B (0-8 min), 90%B (8-10 min)  
 0.4mL/min (6 MPa), 37°C  
 10 µL (500ng/mL, 0.1% formic acid )  
 Single Quad. MS (ESI)