

Cadenza CW-C18


100 x 3 mm

Application

LC-MS: Carbonic anhydrase inhibitor (antiglaucoma drug)

炭酸脱水酵素阻害剤 (緑内障治療薬)

(x100,000)

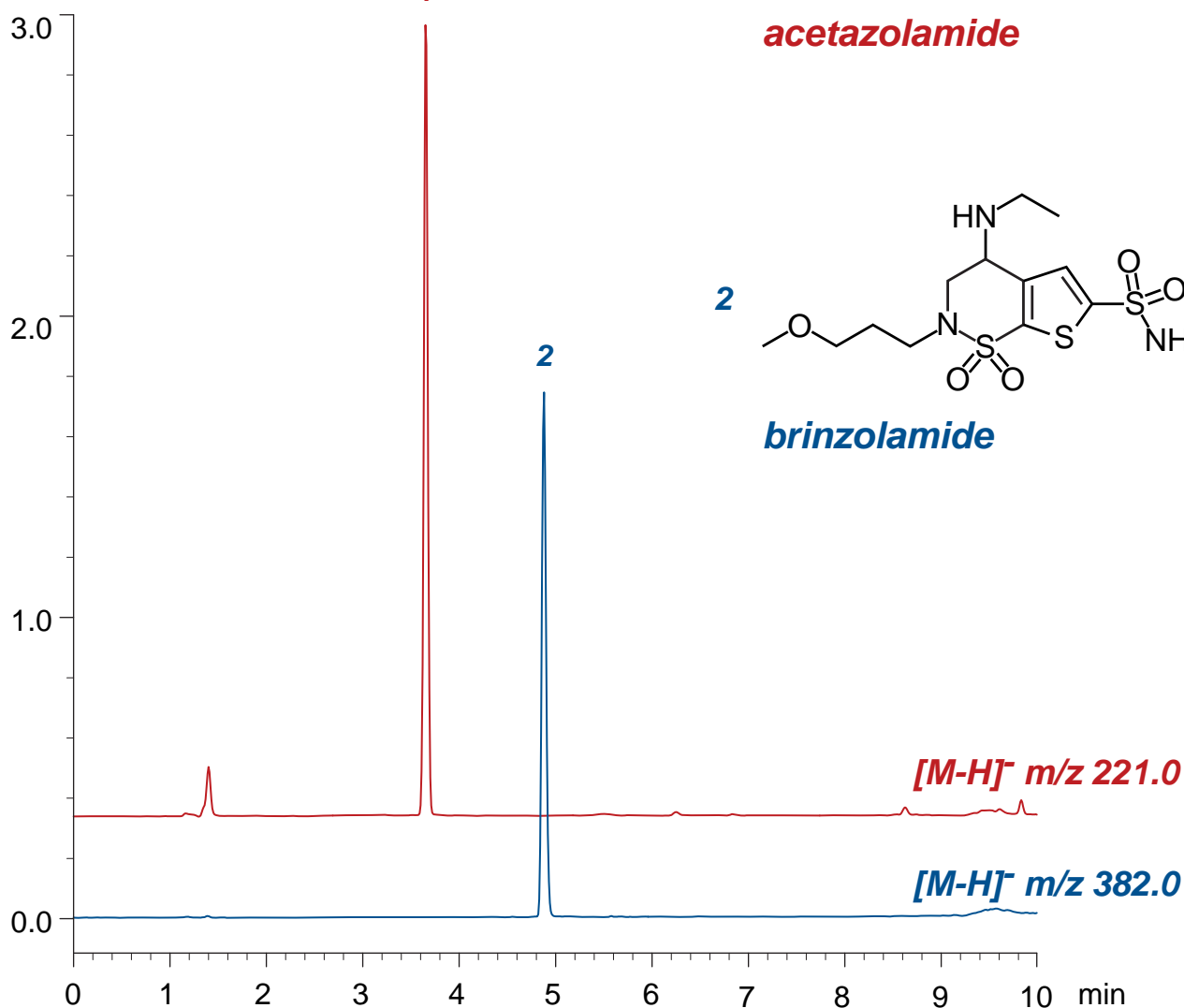

Cadenza CW-C18, 100 x 3 mm  
 A: water/ formic acid = 100 / 0.1  
 B: acetonitrile / formic acid = 100 / 0.1  
 0-50%B (0-8 min), 90%B (8-10 min)  
 0.4mL/min (7 MPa), 37°C  
 10 μL (500ng/mL, 0.1% formic acid)  
 Single Quad. MS (ESI)