








Scherzo SM-C18

150 x 3 mm

Application

LC-MS: biogenic polyamines

生体ポリアミンのLC-MS分析

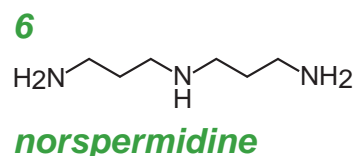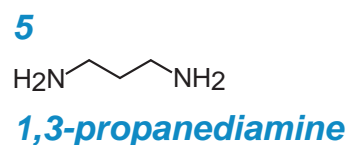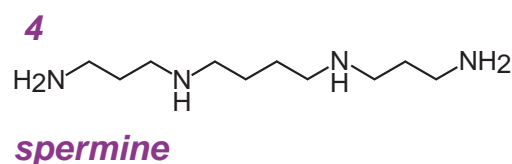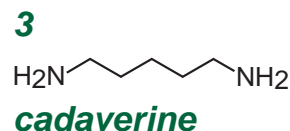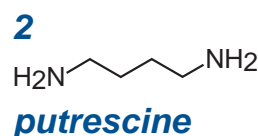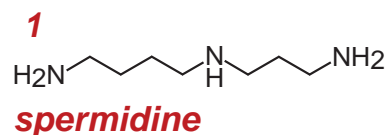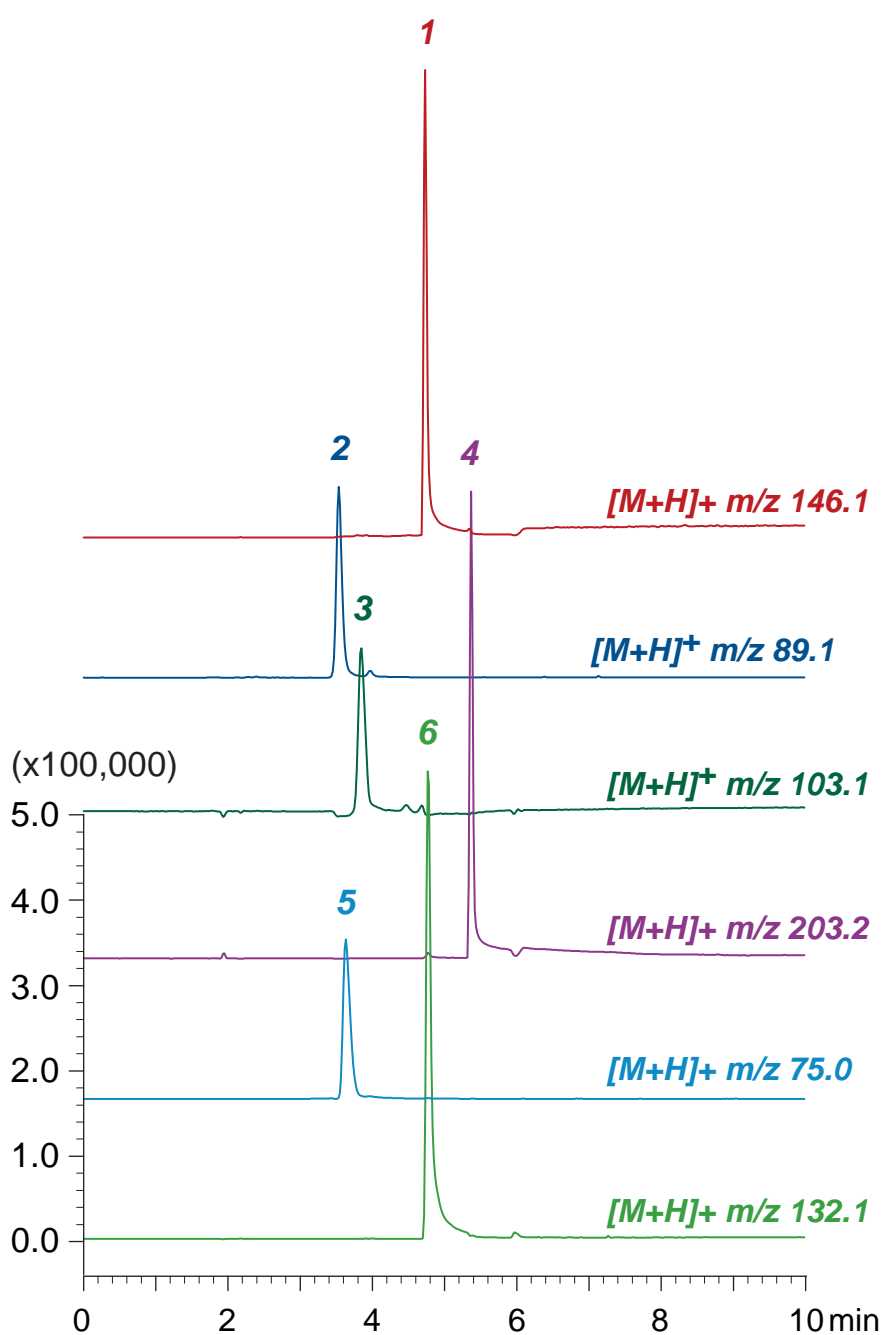

Scherzo SM-C18, 150 x 3 mm

A: 30mM ammonium acetate / acetic acid = 100 / 0.02

B: acetonitrile / water / acetic acid = 20 / 80 / 0.8

0-20%B (0-3.5min), 100%B (3.5-10min)

0.4 mL/min (8 MPa), 37°C, 10 $\mu$ L (5 $\mu$ g/mL, water)

Single Quad. MS (ESI)