


Scherzo SM-C18

150 x 3 mm

Application

LC-MS: aminoglycoside antibiotics

アミノグリコシド系抗生物質のLC-MS分析

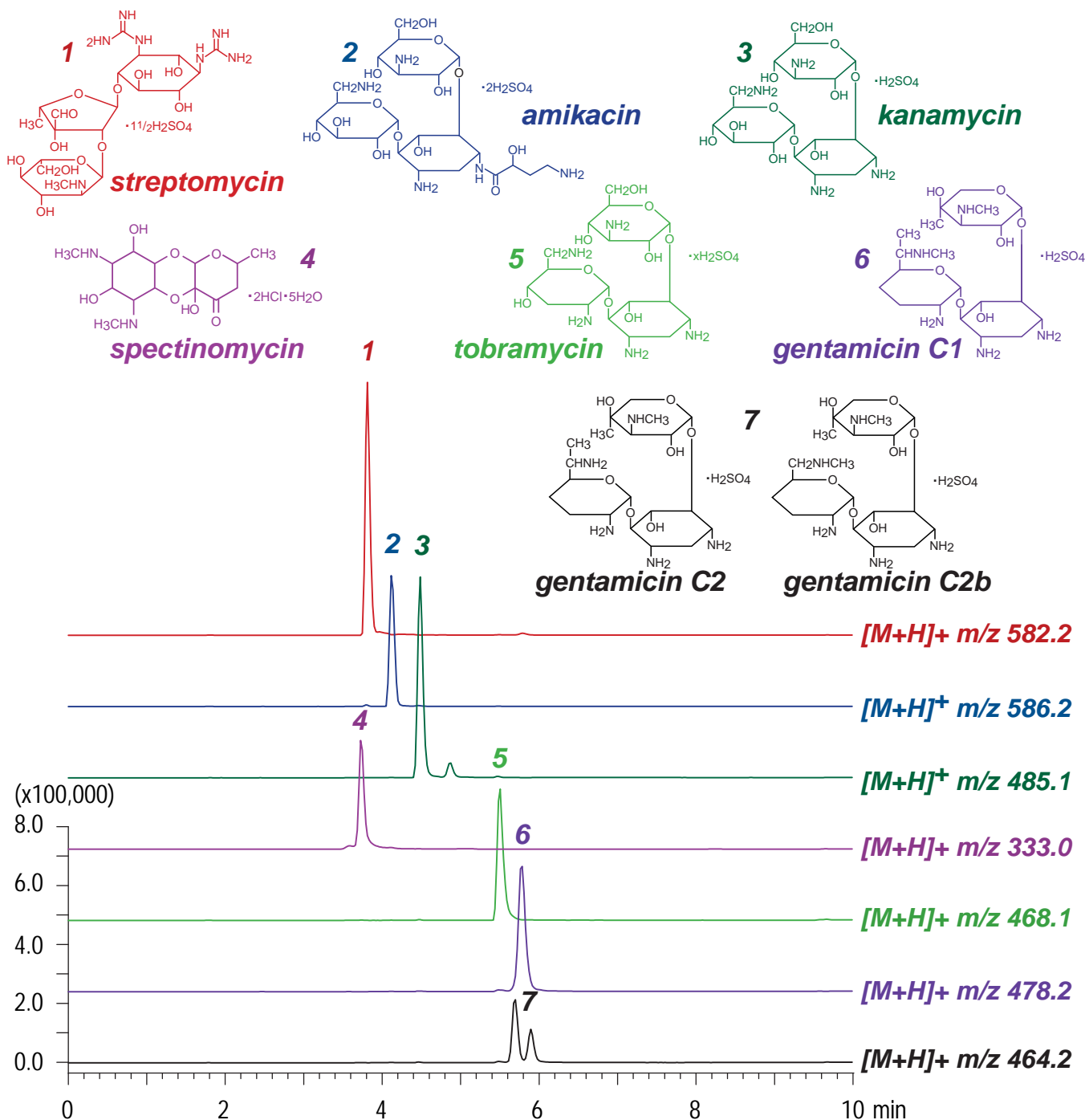

Scherzo SM-C18, 150 x 3 mm

A: 30mM ammonium acetate

B: acetonitrile / water / acetic acid = 20 / 80 / 0.8

5-70%B (0-8min), 100%B (8-10min)

0.4 mL/min (8 MPa), 37°C, 10µL (50µg/mL, water)

Single Quad. MS (ESI)