



Cadenza CD-C18

150 x 4.6 mm

Application

**Rotigotine (Antiparkinsonian drug)**

ロチゴチン (抗パーキンソン病薬)

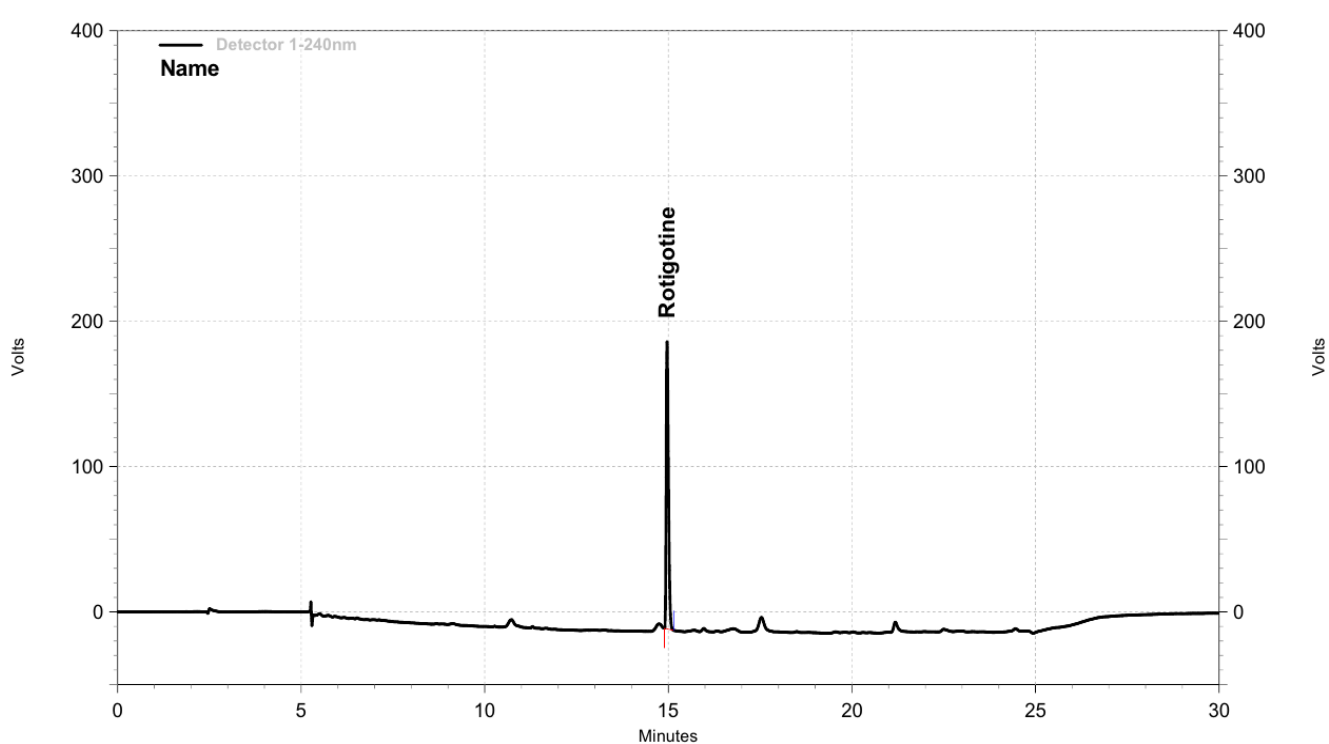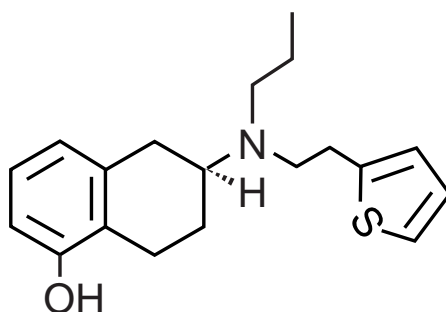

Cadenza CD-C18, 150 x 4.6 mm  
 A: 0.1% phosphoric acid  
 B: acetonitrile  
 0-70%B (0-23,im), 70-0%B (23-23.1min), 0%B (23.1-30min)  
 1 mL/min, 35°C  
 5 µL  
 240nm

Courtesy of PEAKMAN SP, KOREA