



**Cadenza HS-C18**

250 x 4.6 mm

Application

**Nitrovin (animal medicine)**

ニトロビン (動物用医薬品)

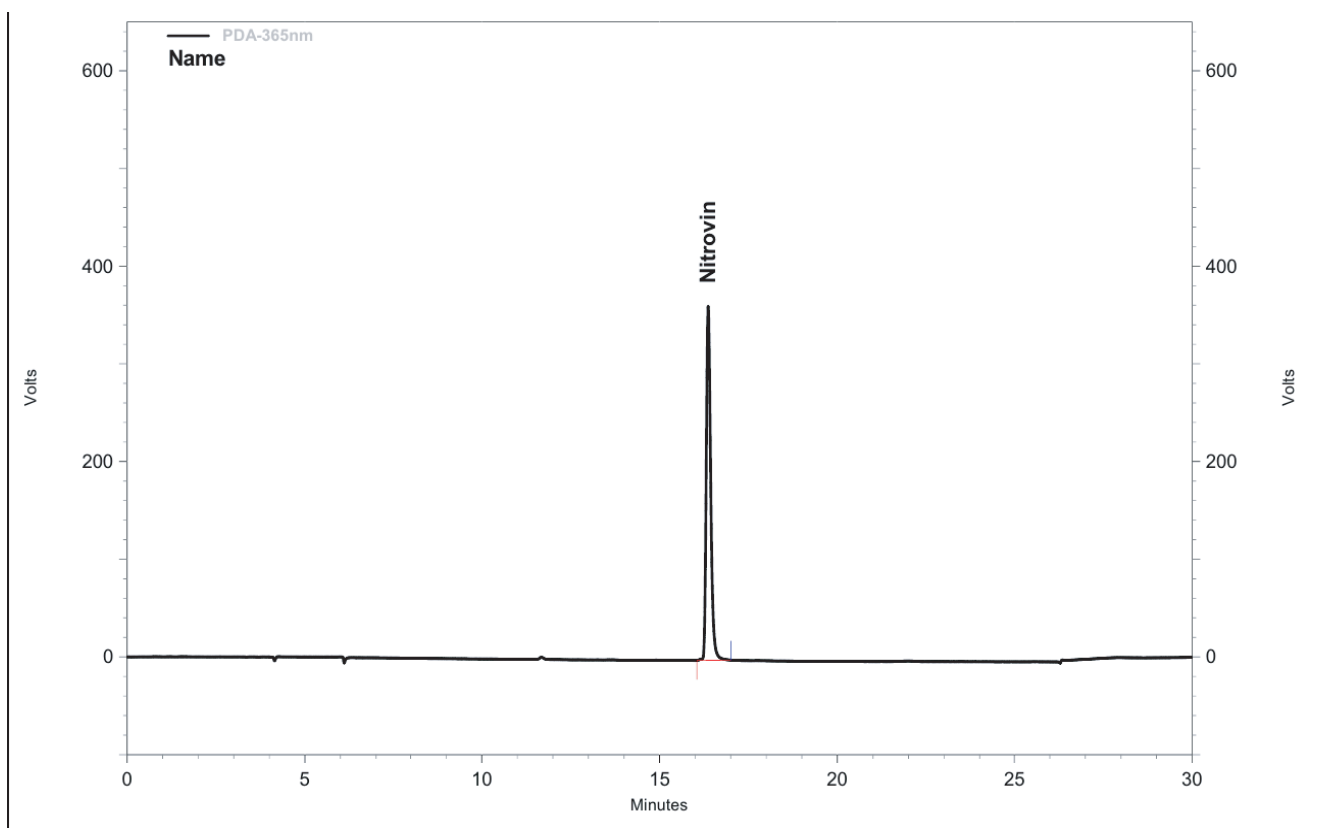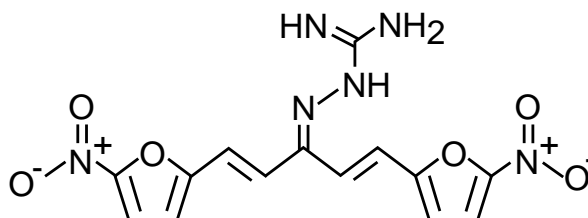

*Cadenza HS-C18, 250 x 4.6 mm*

*A: 0.1% formic acid in water*

*B: 0.1% formic acid in acetonitrile*

*0-80 % B(0~23 min), 80-0% B(23~23.1 min), 0% B(23.1~30 min)*

*0.8 mL/min, 35°C*

*5 µL*

*365nm*

Courtesy of PEAKMAN SP, KOREA