


Intrada Organic Acid

150 x 2 mm

Application

LC-MS: TCA (tricarboxylic acid) cycle / Citric acid cycle

TCA回路 (トリカルボン酸回路) / クエン酸回路のLC-MS分析

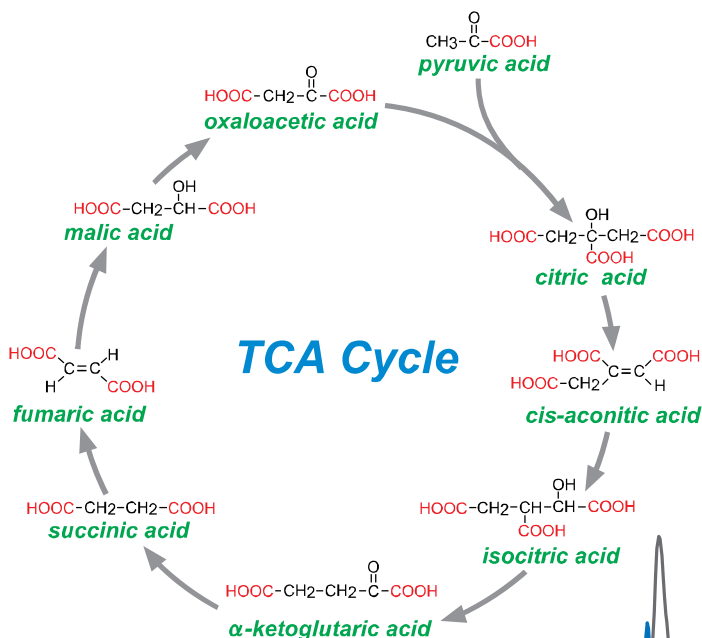

Intrada Organic Acid, 150x2mm

A: ACN / water / HCOOH = 10 / 90 / 0.1

B: ACN / 100mM HCOONH4 = 10 / 90

0%B (0-1min)

0-100%B (1-7min)

100%B (7-10min)

0.2mL/min, 37°C

MS (ESI, SIM, Negative)

5 μ L (100 μ g/mL with 0.1% HCOOH)

