






Scherzo SM-C18 MF

100 x 2 mm

Application

LC-MS: Organophosphorus and quaternary ammonium pesticides/herbicides

有機リン系および四級アンモニウム系農薬のLC-MS分析

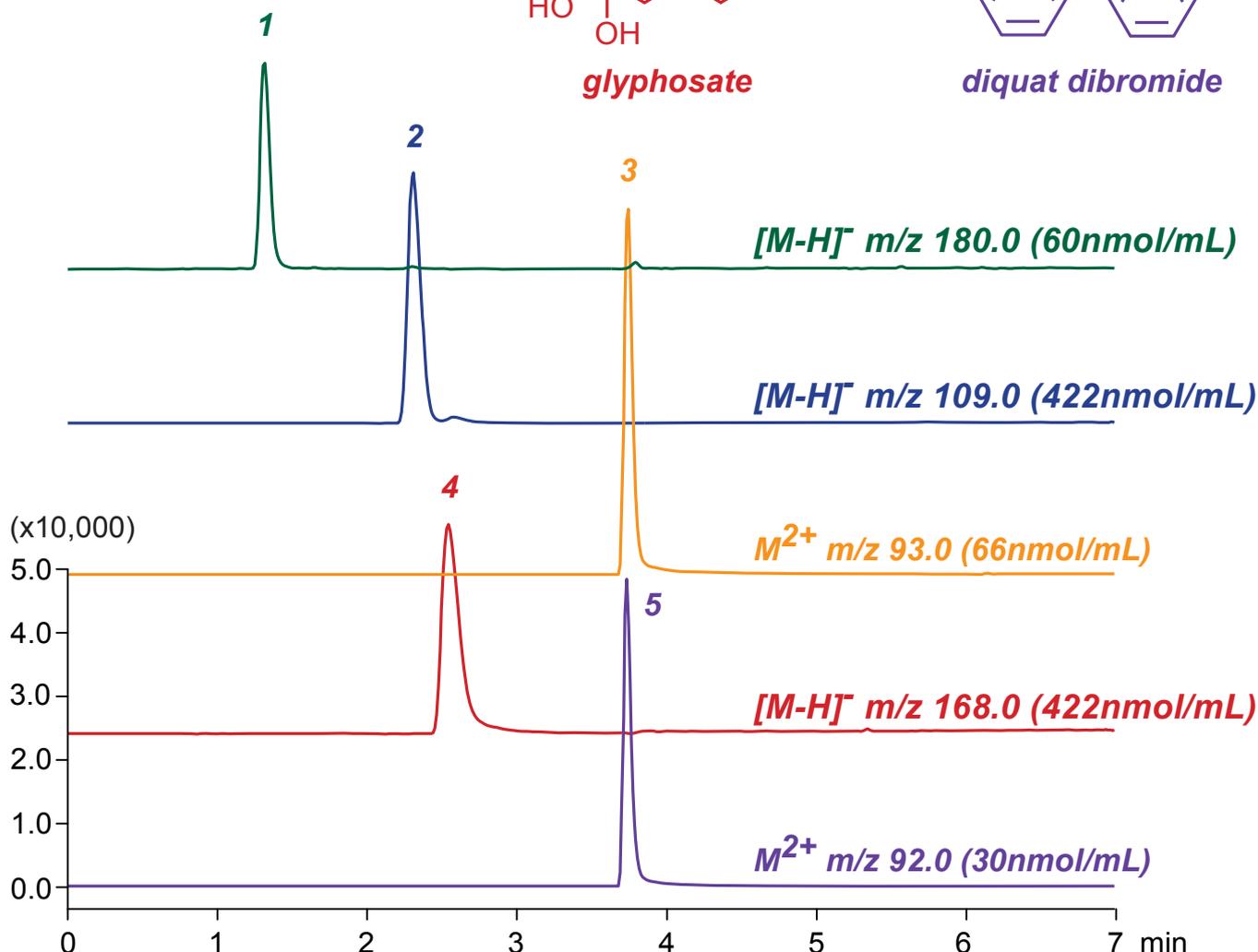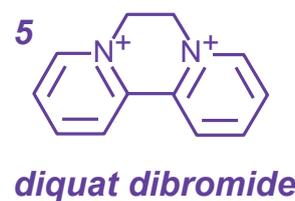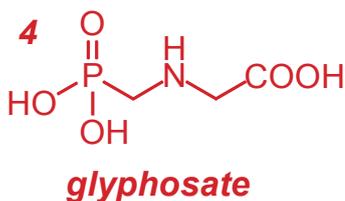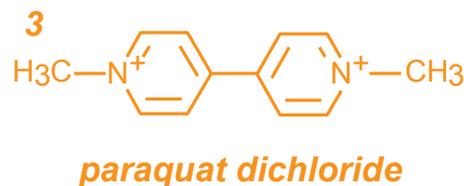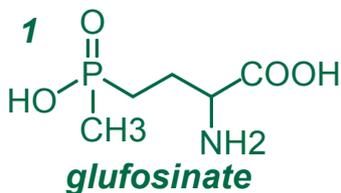

Scherzo SM-C18 MF, 100 x 2 mm

A: 15mM ammonium formate

B: acetonitrile / water / formic acid = 20 / 80 / 0.5

0%B (0-1min), 0-100%B (1-5min), 100%B (5-7min)

0.2 mL/min (6MPa), 37°C, 5 μL (30-422nmol/mL, 15mM HCOONH₄)

Single Quad. MS (ESI)