



Cadenza CX-C18
 Cadenza CD-C18
 Unison UK-C18
 Unison UK-C8
 Unison UK-Phenyl

75 x 4.6 mm

Application

Tadalafil (phosphodiesterase type 5 inhibitor, PDE5I)

タダラフィル(ホスホジエステラーゼ5阻害剤)

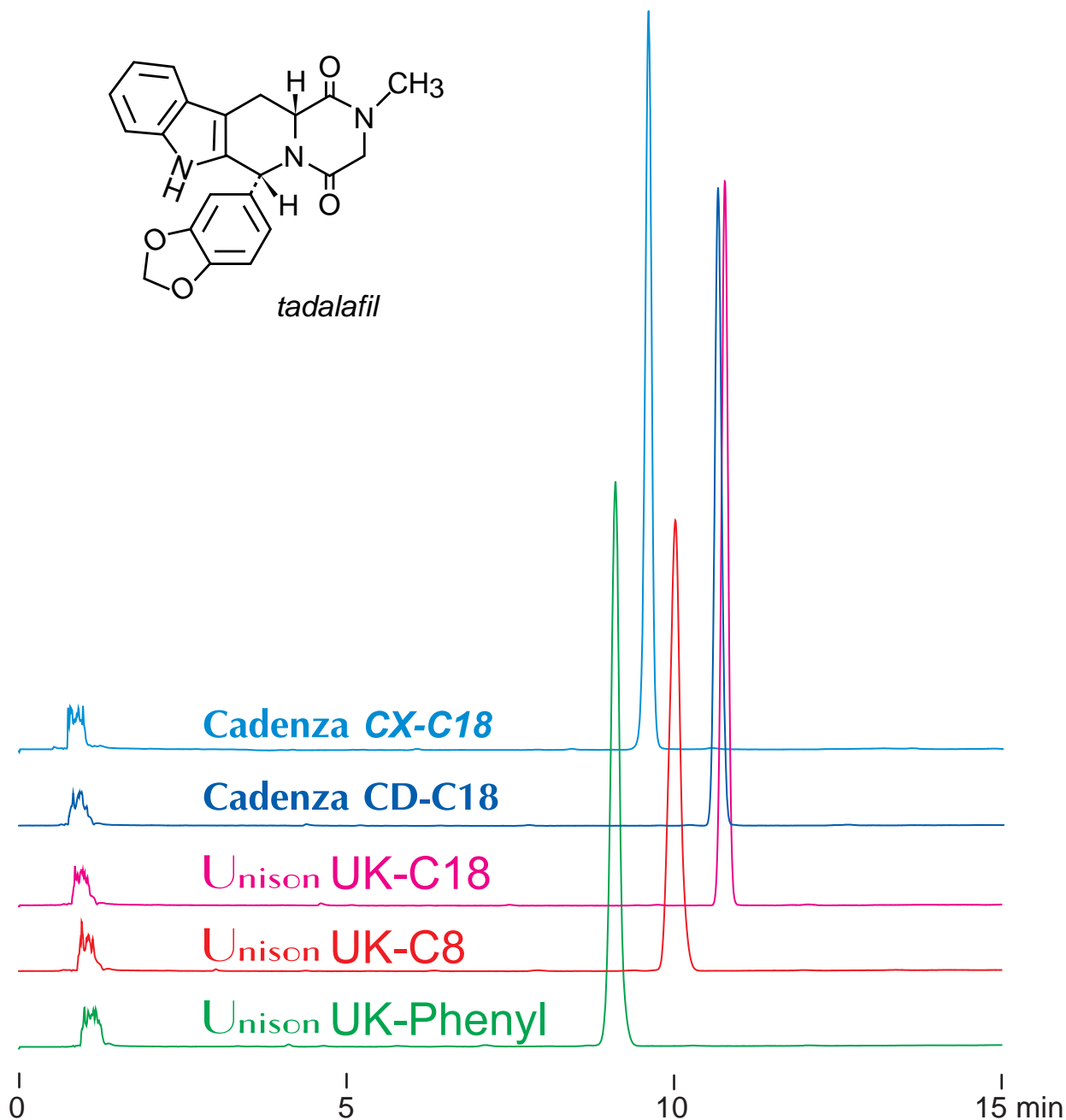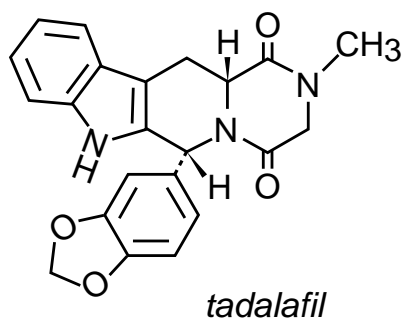

75 x 4.6 mm
 A: water / formic acid= 100 / 0.1
 B: acetonitrile / formic acid = 100 / 0.1
 20-50%B (0-15min)
 1 mL/min (7MPa), 37°C
 10 µL (0.1mg/mL, DMSO), 280nm