

**Cadenza CD-C18**



**Unison UK-C18**

50 x 3 mm


Application

LC-MS: Varenicline (Smoking cessation adjuvant)

バレニクリン(禁煙補助薬)のLC-MS分析

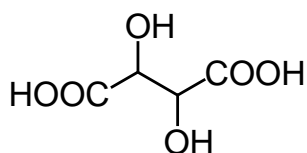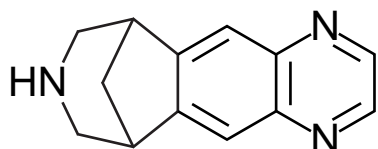

varenicline tartrate

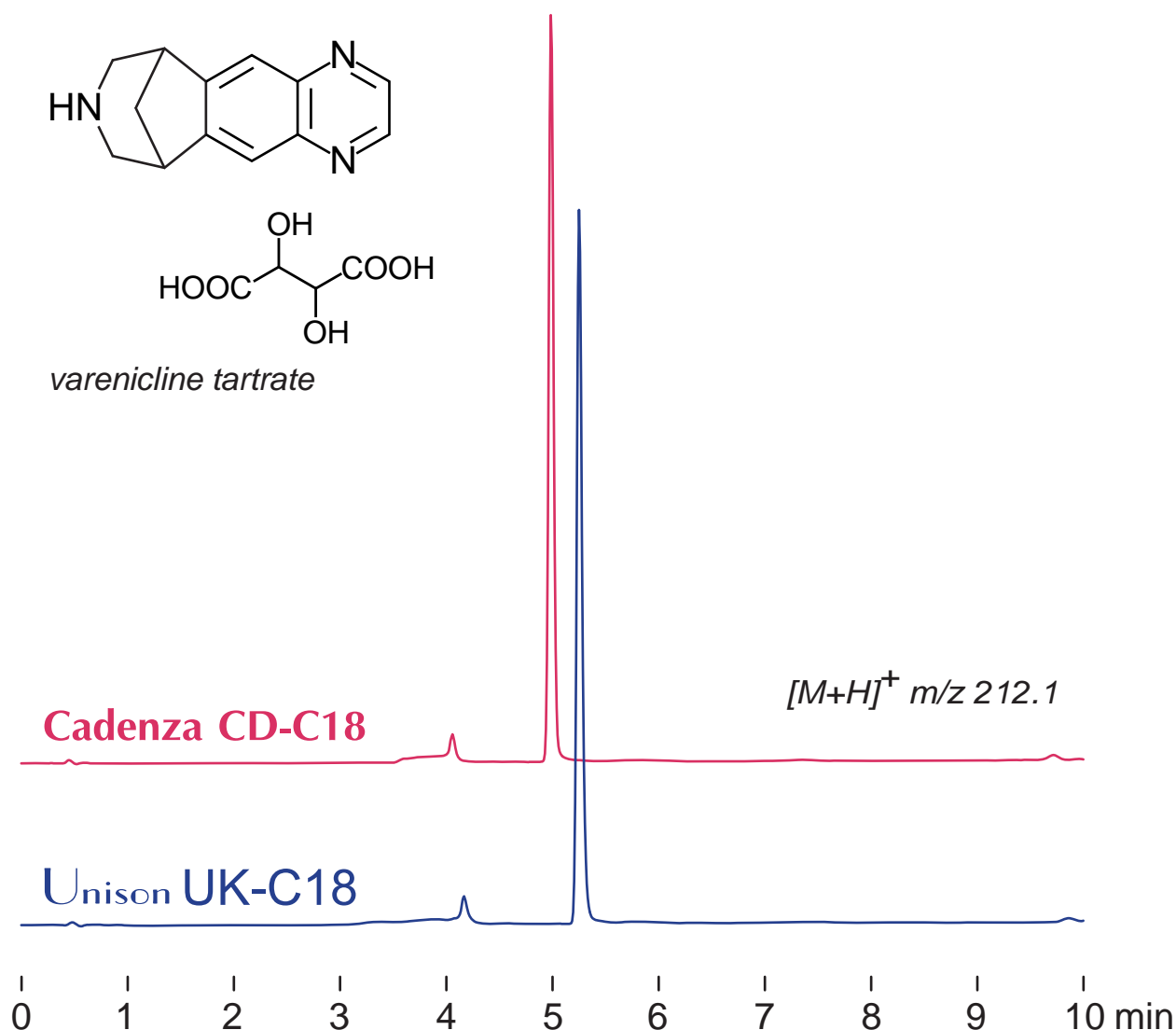

50 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0%B (0-1min), 0-30%B (1-10min)

0.4mL/min (4MPa), 37°C, MS (ESI, SIM, positive)

5µL (1µg/mL, 0.1% HCOOH)